

# COMPLEX APPROACH *with*




Deta  
Elis  
Holding  
Health  
Innovation  
Human



# HEALTH COMPONENTS

BY DEHOLDING COMPANY

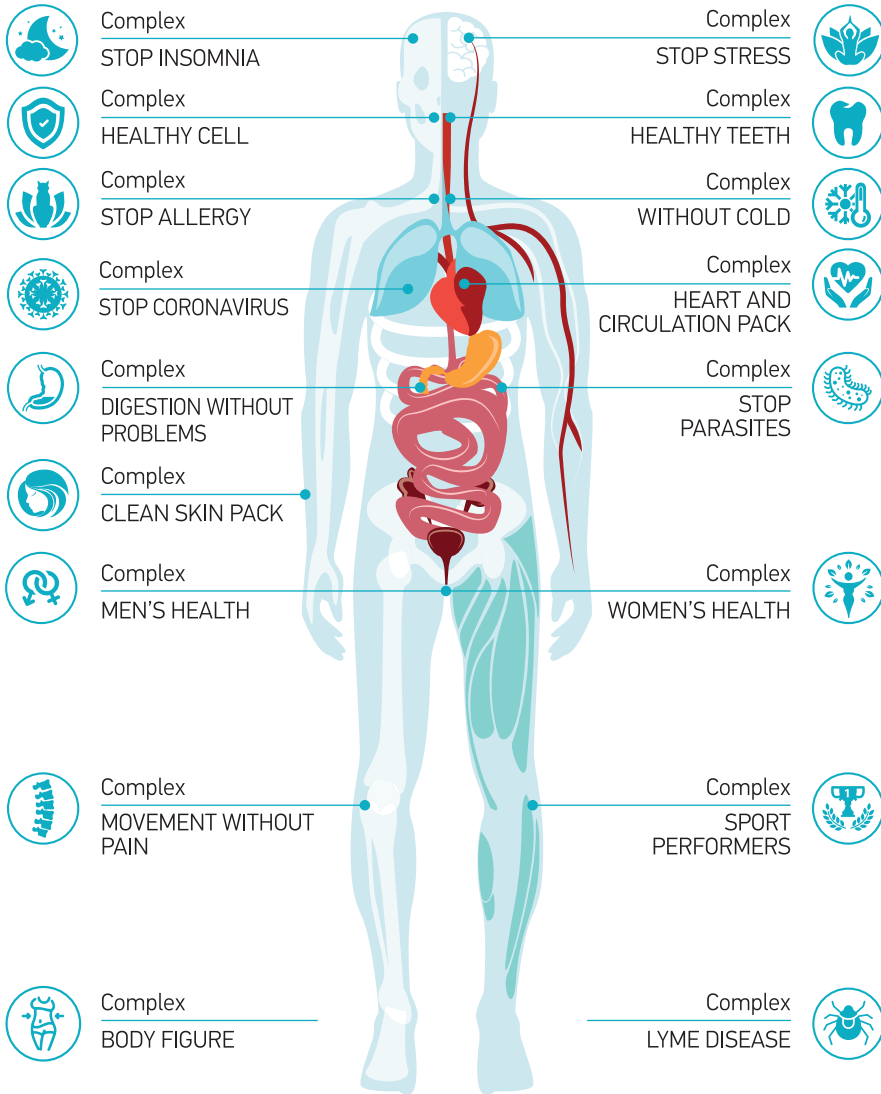

FOR BEST RESULTS, IT IS RECOMMENDED TO USE THE DEVITA DEVICES WITH THE DELIXIR SMART CELL NUTRITION SUPPLEMENTS

# COMPLEX HEALTH FOR YOUR ORGANISM

## COMPLEX "WITHOUT COLD"

The action of the complex is aimed at enhancing the action of the immune system, increasing resistance and reducing recovery time.

## COMPLEX "STOP ALLERGY"

The action of the complex is aimed at eliminating biological allergens, increasing immunity, cleansing of toxins.

## COMPLEX "STOP STRESS"

The complex is aimed at harmonizing the nervous system, relieving chronic fatigue and improving overall well-being.

## COMPLEX "HEART AND CIRCULATION"

The action of the complex is aimed at improving blood microcirculation and harmonizing heart rate.

## COMPLEX "STOP INSOMNIA"

The action of the complex is aimed at improving resistance to mental and psychological stress, reducing nervous excitement.

## COMPLEX "DIGESTION WITHOUT PROBLEMS"

The action of the complex is aimed at improving the secretory and motor functions of the digestive system and restoring the functions of the excretory system.

## COMPLEX "WOMEN'S HEALTH"

The action of the complex is aimed at restoring the functions of the excretory system and reducing psychoemotional stress.

## COMPLEX "MAN'S HEALTH"

The action of the complex is aimed at protecting against infections, regulating sexual activity and restoring physiological parameters.

## COMPLEX "MOVEMENT WITHOUT PAIN"

The action of the complex is aimed at relieving joint inflammation, recovering from injuries and improving the functioning of the musculoskeletal system.

## COMPLEX "CLEAN SKIN"

The action of the complex is aimed at getting rid of peeling and inflammation, increasing the immunity of the skin and improving the external condition of the skin, hair and nails.

## COMPLEX "SPORT PERFORMERS"

The action of the complex is aimed at increasing strength indicators, normalizing muscle tone and relieving pain after training.

## COMPLEX "STOP PARASITES"

The action of the complex is aimed at detoxifying the body, increasing immunity, harmonizing the work of internal organs and recovering from cleansing from parasites.

## COMPLEX "BODY FIGURE"

The action of the complex is aimed at normalizing metabolic processes, reducing body weight, improving the psychoemotional state and quality of life.

## COMPLEX "HEALTHY CELL"

The action of the complex is aimed at reducing the negative effects of the environment and toxic load, preventing the development of undesirable consequences for the body.

## COMPLEX "HEALTHY TEETH"

The action of the complex is aimed at reducing bleeding of the gums, improving the health of the mucous membranes, strengthening local immunity, and early recovery from dental interventions.

## COMPLEX "LYME DISEASE"

The action of the complex is aimed at reducing parasitic and toxic loads, maintaining harmonious functioning of the nervous system, increasing the energy state and improving overall well-being.

## COMPLEX "STOP CORONAVIRUS"

The action of the complex is aimed at enhancing the defences of the organism, increasing resistance, reducing the likelihood of complications and shortening the recovery time..

## INTRODUCTION

The complex approach to health by the DEHolding Company is based on the homotoxicology theory by Dr. H.H. Reckeweg. Dr. Reckeweg introduced the concept of "homotoxin", meaning any substance toxic to humans, regardless of their nature. They can be endotoxins (metabolic products or unsplit/not removed in sufficient quantities) or exotoxins (harmful environmental factors, pathogens of infectious diseases, parasites, etc.). Due to the accumulation of such substances in the human organism, certain diseases can appear. Homotoxicology considers health as "A state of freedom of an organism from homotoxins and the absence of functional or organic lesions caused by them".

The complex approach makes it possible to gently, without stress for the organism, normalize the internal environment of the body, increase energy, and improve the quality of life. To achieve this goal, before using any of the special complexes, it is necessary to take a course of a gastrointestinal tract restoration using the basic complexes "Digestion without problems" or "Stop parasites". They are fundamental for releasing the organism from toxic load and normalizing the activity of all excretory organs.

# CONTENTS

## BASIC COMPLEXES

Digestion without problems.....	7
Stop parasites.....	11

## THEMED COMPLEXES

Heart and circulation.....	17
Stop stress.....	21
Stop insomnia.....	25
Women's health.....	29
Men's health.....	33
Movement without pain.....	37
Clean skin.....	41
Sport performers.....	45
Stop allergy.....	49

## PARTICULAR COMPLEXES

Active immune system.....	55
Stop Coronavirus.....	59

## SPECIAL COMPLEXES

(Used only after one of the basic complexes)

Body figure.....	65
Healthy Cell.....	69
Lyme Disease.....	73
Healthy teeth.....	77


## PROGRAMS AND DELIXIR

The list of manual programs DeVita AP.....	81
The list of manual programs DeVita Ritm.....	91
Elixir of healthy life.....	98
DeLixir Multi-Energy, DeLixir Detox.....	99
DeLixir Collagen +, DeLixir pH Balance.....	100
DeLixir Ultra Slim +, Delixir Top Form.....	101

## REFERENCE INFORMATION

Explanations for the use of programs DeVita Ritm.....	102
Explanations for the use of programs DeVita AP.....	108
The biological role of DeLixir.....	118
Recommended DeVita programs for use with various functional problems.....	120

# BASIC COMPLEXES

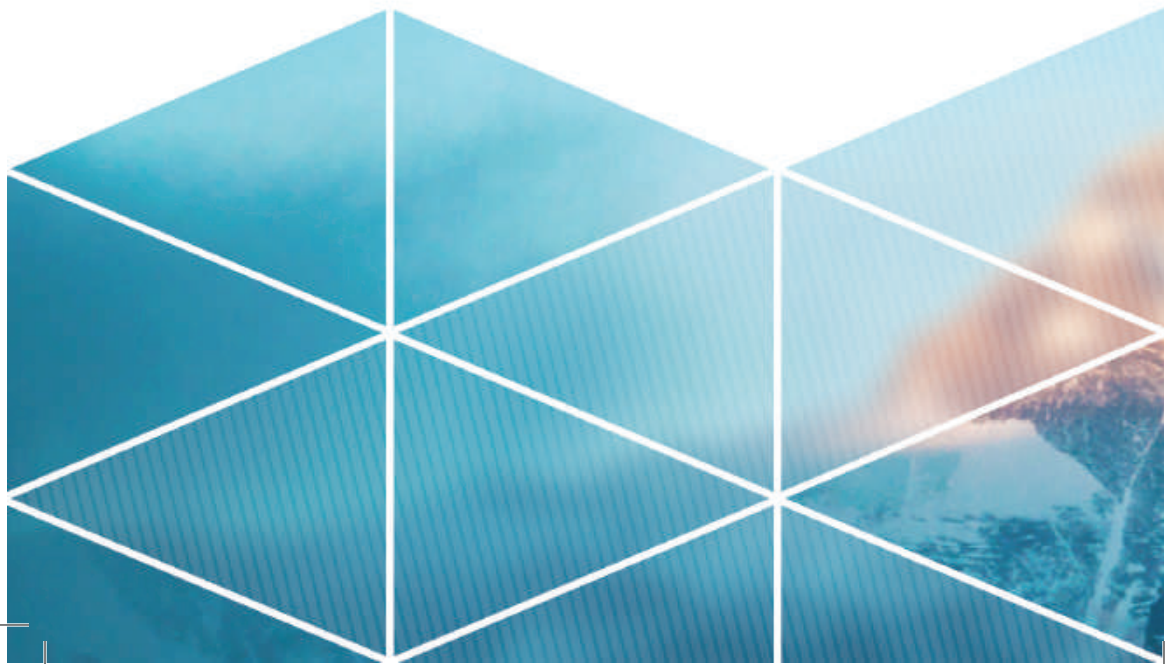



## COMPLEX DIGESTION WITHOUT PROBLEMS

WHY IS IT IMPORTANT MAINTAIN HEALTH OF THE GASTROINTESTINAL TRACT?



The gastrointestinal tract is a complex system on which the vital activity of the whole organism depends. In this system the food is processed into substances necessary for the functioning of organs. Also, the processes of obtaining energy and building material for new cells begin in the gastrointestinal tract.



More than half of the world's population suffers from digestive disorders. Their main causes are an unbalanced diet, poor ecology, poor eating habits (overeating, diets, nutrition on the go, convenience foods, etc.). One in four people has a history of a disease of the gastrointestinal tract.



This leads to serious consequences - from problems with body weight to diseases that require surgery. In addition, problems with the gastrointestinal tract also affect overall well-being, performance and appearance. Improper metabolism causes skin problems - acne, dryness or greasiness, dull complexion, wrinkles.



Common digestive problems are pain after eating and on an empty stomach, heartburn, stool disorders, nausea and flatulence. A person is forced to follow a special diet. Discomfort associated with the stomach and intestines affect the lifestyle, making it less active, and puts a limit on travelling.

## MECHANISM OF ACTION




**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Parasites disrupt the secretory and motor functions of the stomach, pancreas, liver and intestines. Parasites can block the ducts of the liver and pancreas with their bodies, causing their inflammation. Helminths are also conductors for other microbes. DeVita AP basic antiparasitic programs affect most known pathogens, quickly, effectively and without side effects, removing the parasitic load from the organism.

The special program "No schistosomes (blood-flukes) & trematods" relieves from the parasitic burden of the liver and blood, restoring their functions and reducing intoxication.

Special programs "Healthy Gastrointestinal", "Help against poisoning", "Drainage Antitox" improve the condition of intoxication or intestinal flu, neutralize and remove toxins.

"Drainage -General" program is used to cleanse the body after applying basic antiparasitic programs, reduces parasitic intoxication and improves overall well-being.

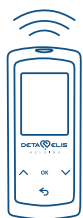

**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

"Healthy Day" automatic program normalizes the work of the stomach and intestines, which is especially important for people with a violation of the diet, helps to increase the body's defenses and harmonize the work of all systems.

Special programs "Digestion normalization" and "Intestine without problems" are aimed at improving the functioning of the pancreas, normalizing the functions of the intestines, which is important not only to restore overall digestion, but also to improve absorption, reduce severity and bloating, restore immunity, metabolism normalization.



**DeVita Energy** restores physiological processes in the digestive system, starts the mechanism of self-regulation.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating "junk food" (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- pain, bloating;
- heartburn, a feeling of heaviness;
- irregular bowel movements.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir pH BALANCE and DeLixir DETOX smart cell nutrition supplements.



## APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Digestion normalization and / or Intestine without problems	Daily in the afternoon, if needed.
	4	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	5			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	6			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
	7			Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.









SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Digestion normalization and / or Intestine without problems	Daily in the afternoon, if needed.
	4	DeVita AP	Special programs	Cleansing-automatic	Daily at night.
	5			No schistosomes (blood-flukes) & trematods	Daily in the evening after DeVita Ritm programs for 1 <sup>st</sup> -2 <sup>nd</sup> weeks.
6	Drainage Antitox			Daily in the evening after DeVita Ritm programs for 3 <sup>rd</sup> -4 <sup>th</sup> weeks.	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed, at night.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Digestion normalization and / or Intestine without problems	Daily in the afternoon, if needed.
4	DeVita AP	Special programs	No nematodes and ascarides, No parasites - General, No parasites - Digestive track	Daily in the evening after the DeVita Ritm programs, alternating 1 program every 1-2 weeks.	

IN CASE OF INTOXICATION	DEVICE	PROGRAM		USAGE
	DeVita Ritm	Special programs	Detoxification, Intestine without problems	
DeVita AP	Special programs	Healthy Gastrointestinal, Help against poisoning, Drainage Antitox		Daily for 10 days.

# PREVENTIVE MEASURES FOR DISEASES OF THE GASTROINTESTINAL TRACT

-  Drink fluids 15-30 minutes before or 2-4 hours after eating, but not with meals. Fluids dilute the gastric acid and interfere with digestion.
-  Eat fruits separately from other foods: either 30 minutes before or 3-4 hours after eating.
-  Obtain moderation in eating. Get up from the table with a feeling of slight hunger. The signal of the stomach about saturation does not immediately go to the brain. By eating up to satiety, you risk overeating.
-  Eat freshly cooked foods. It is advisable to cook food immediately before eating it.
-  Eat small portions 4-5 times a day. You should not experience severe hunger in between meals, but you should not eat if you don't feel hungry.
-  When eating, do not get distracted by reading, TV, long conversations. Also, do not eat in a bad mood. All this interferes with the proper digestion of food.
-  Eat slowly, chew for a long time and carefully so that the food is saturated with saliva as much as possible.
-  In the evening, it is advisable not to eat 2 hours before bedtime. A feeling of hunger can be quenched with something light - for example, low-fat cottage cheese or kefir.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before use



Pain, discomfort in the abdomen, heartburn, heaviness, bloating.

After the first use of DeVita



Reduced pain, relieve from heartburn and heaviness, improved mood.

After a month of using DeVita



Normalization of digestion, elimination of exacerbation of chronic diseases, excellent health.



Use the scheme of the "DIGESTION WITHOUT PROBLEMS" complex together with the DeLixir pH BALANCE and DeLixir DETOX smart cell nutrition supplements, and your results will be noticeable much faster.



## COMPLEX STOP PARASITES

### WHY PARASITE INFESTATIONS ARE DANGEROUS?



Parasites enter the body through skin, food, water, dust, from pets or other people. They can settle in all organs and tissues, and also live under the skin or on the skin. Parasites absorb vital elements, secrete toxins, destroy blood cells.



Parasites trigger many chronic diseases of the digestive tract, bronchopulmonary system, and skin. Any health-improving measures (proper nutrition, physical exercises, hardening procedures) without ridding the body of parasites will not give a positive effect.



Parasitic diseases are dangerous because they are perfectly disguised as other diseases. Symptoms of parasite infection include chronic fatigue, irritability and anxiety, hyperactivity in children, anemia, brittle nails and hair, headaches, appetite disorders, decreased immunity.



Treatment for parasites, especially in advanced cases, is very aggressive towards the body. Antiparasitic drugs have a negative effect on the liver and gastrointestinal tract, killing beneficial microflora. Therefore, it is much easier to prevent the entry and multiplication of parasites in the body.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Bacteria, viruses, fungi, parasites or their toxins are the main cause of poor health and one of the leading causes of death. Parasites weaken the immune system, increase the degree of intoxication, affect organs, change the pH of the internal environment, and create an imbalance in the body. The basic anti-parasitic programs on the DeVita AP device affect the majority of known pathogens, quickly, efficiently and without side effects, removing the parasitic load and signs of inflammation.



“Drainage - General” automatic program is used to cleanse the body after applying basic antiparasitic programs, reduces parasitic intoxication and improves overall well-being.

**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

“Healthy Day” automatic program improves blood circulation in organs, relieves from stress, which helps to increase immunity.



“Complete Cleansing” automatic program affects organs taking into account their daily activity, gently cleansing the body of accumulated toxins and harmonizing the work of all organs and systems.

**DeVita Energy** activates self-regulation mechanisms, helps in the acceleration of the recovery processes and the distribution of energy balance. Allows the body to quickly get in shape after cleansing from parasites, maintains vigor and good mood.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- inflammatory processes, the appearance of allergies, the deterioration of the skin, hair and nails;
- a change in appetite, increased fatigue, impaired bowel movements;
- in the presence of unexplained pain;
- to maintain health in high-risk groups: pet owners, travelers.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir pH BALANCE and DeLixir DETOX smart cell nutrition supplements.

## APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Intestine without problems	Daily in the afternoon.
	4	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	5			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	6			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
	7			Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Intestine without problems	Daily in the afternoon.
	4	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night.
	5		Special programs	No bacteria	Daily in the evening after DeVita Ritm programs, 14 days.
6	No herpes simplex, No herpes, No hepatitis A			If necessary, according to guidelines.	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed, at night.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Intestine without problems	Daily in the afternoon.
4	DeVita AP	Special programs	No nematodes and ascarides, No parasites- General, No parasites- Digestive track	Daily in the evening after DeVita Ritm programs, alternating each program every 1-2 weeks.	

## PREVENTIVE MEASURES FOR PARASITIC INFECTIONS



Thoroughly heat food and water before use. Boiling is preferable to frying. Try not to eat meat of wild animals, river fish, as well as dried and slightly salted fish.



Before you travel to tropical countries, take the necessary vaccinations. Be wary of local food.



Try not to walk barefoot on the beach - nematodes live in the sand.



After contact with stray animals, wash your hands with alkaline soap (household) - a soft toilet soap may not completely remove the parasite larvae.



Wear gloves when working in the garden.



If you travel often, walk in nature, work in the garden, have contact with pets, carry out drug prophylaxis of the parasites once a year.



Wash fruits and vegetables with soap and brush for at least 1 minute.



If you have problems with the gastrointestinal tract, muscle pain, general malaise, consult a doctor.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Headaches, nausea, lingering cough, feeling tired, weakness, impaired bowel movements.

After the first use  
of DeVita



Improved general condition of the organism.

After a month of  
using DeVita




Energy boost, disappearance of pain, normalization of all systems.



Use the scheme of the "STOP PARASITES" complex together with the DeLixir pH BALANCE and DeLixir DETOX smart cell nutrition supplements, and your results will become noticeable much faster.



# THEMED COMPLEXES

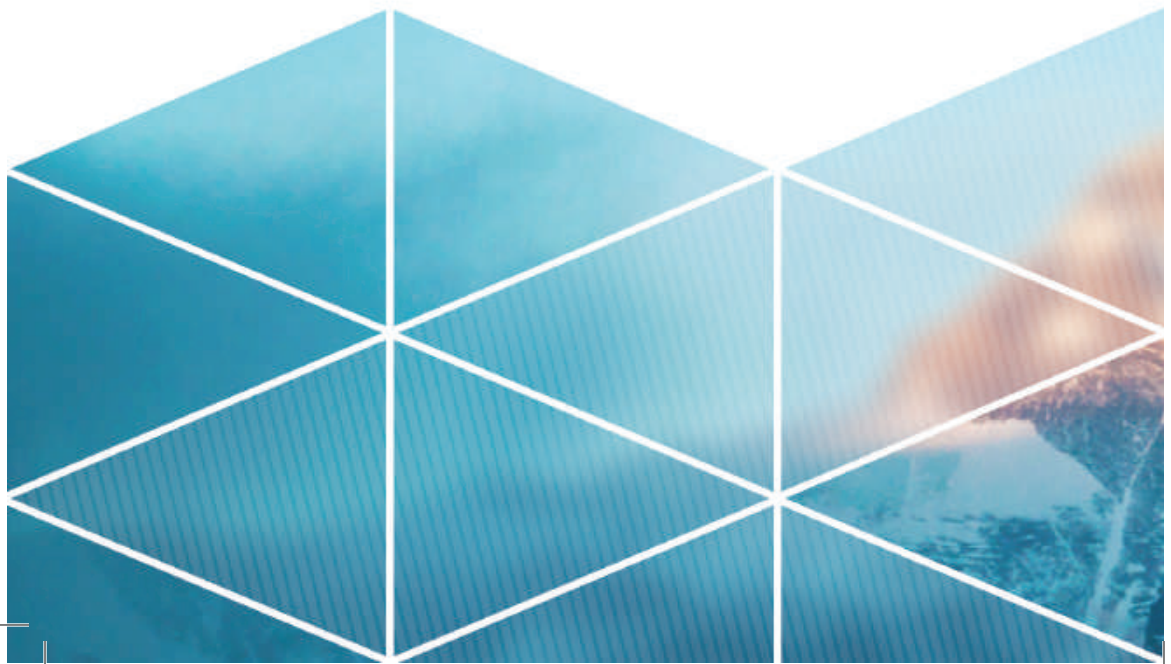





## COMPLEX HEART AND CIRCULATION

### HOW DANGEROUS ARE HEART AND VASCULAR DISEASES?



The main function of the cardiovascular system (CVS) is to provide organs with nutrients, oxygen and energy. Also, decay products and toxins «leave» the blood from organs. The cardiovascular system works constantly, without interruptions, and the slightest malfunctions in its work lead to problems throughout the body.



CVS diseases occupy a leading position among other diseases. Risk factors are stress, negative habits, inactivity, which are seen in the life of almost every person.

CVS diseases are also caused by birth defects, injuries, inflammations, metabolic disorders.



The most common diseases are coronary artery disease, hypertension, arrhythmia, and atherosclerosis. They can lead to heart attacks, strokes and hypertensive crises - the main causes of death. In developed countries, these problems appear not only in older age groups - cases of heart attack at 30-35 are known.



Problems of the cardiovascular system are dangerous because they pass asymptotically for a long time. And their first signs are similar to ordinary fatigue - shortness of breath, dizziness, nausea, migraines. Often a person seeks help when the problems with the heart and blood vessels are already in a developed stage.

## MECHANISM OF ACTION




**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Parasites, entering the human body, exert a toxic load on the cardiovascular system and impair its functioning. Microbes destroy the heart and blood vessels, migrating along the vascular system and settling in the heart muscle and its surrounding membranes. DeVita AP basic antiparasitic programs the parasitic most known pathogens, quickly, effectively and without side effects, removing the parasitic load.

A special program “No schistosomes (blood-flukes) & trematods” cleans the blood of parasites that affect the condition of blood vessels and the human heart.

The “Drainage, General” program is used to cleanse the body after applying basic antiparasitic programs, removing toxins and improving overall well-being.

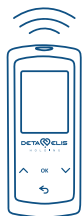

**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

“Healthy Day” automatic program has a beneficial effect on the functioning of the cardiovascular system, helping to open even small capillaries. Improves blood microcirculation, nutrition of tissues and cells. It activates the process of body self-regulation and increases the efficiency of the immune system. The automatic program “Complete Cleansing” helps in the restoration of the excretory system and in harmonizes the work of the body.

Special programs “Healthy heart” and “Cardiac rhythm normalization” restore the heart rhythm, reduce the influence of negative factors, provide normal blood supply to the heart.

**DeVita Energy** relieves from psycho-emotional stress, which is one of the causes of cardiovascular disease. It starts the processes of self-regulation of the body and restores energy.



## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- with a violation of the heart rhythm (arrhythmia)
- with problems with veins and arteries
- with circulatory failure
- with general problems of the cardiovascular system.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir MULTI-ENERGY, DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition.

## APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Healthy heart, Cardiac rhythm normalization	Daily afternoon, if needed.
	4	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	5			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	6			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
	7			Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Healthy heart, Cardiac rhythm normalization	Daily afternoon, if needed.
	4	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night.
	5		Special programs	No schistosomes (blood-flukes) & trematods	Daily in the evening after DeVita Ritm programs for 1-2 weeks .
6	No chlamydia			Daily in the evening after DeVita Ritm programs for 3-4 weeks.	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed for the night.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Healthy heart, Cardiac rhythm normalization	Daily afternoon, if needed.
4	DeVita AP	Special programs	No nematodes and ascarides, No parasites- General, No parasites- Digestive track	Daily in the evening after the DeVita Ritm programs, alternating 1 program for 1-2 weeks.	

## PREVENTION MEASURES FOR CARDIOVASCULAR SYSTEM PROBLEMS



Reduce your salt intake. Salt retains water, which increases blood volume and leads to increased pressure.



Watch your blood pressure. Normal indicators are 120/80 mm Hg.



Once a year do a cardiogram.



Watch your blood glucose and cholesterol levels. Excess sugar and cholesterol cause the creation of deposits on the walls of blood vessels and clog them.



Treat tonsillitis and other bacterial infections in time. Chronic inflammation puts a strain on the heart and causes intoxication.



Eat boiled or steamed fish and poultry, avocados, linseed oil, nuts, pumpkin, broccoli, pomegranate, garlic and onions, ginger, turmeric, oat and corn bran, decoctions of hawthorn and rose hips, dried apricots and figs.



Limit refined oil, fatty and fried, sugar and salt.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Fatigue when standing up and walking fast. General fatigue, memory and efficiency loss.

After the first  
use of DeVita



Increased overall body tone,  
heart rate alignment.

After a month of  
using DeVita



Increased endurance, improved  
vascular and memory status,  
lack of shortness of breath.



Use the "HEART AND CIRCULATION" Complex scheme together with DeLixir MULTI-ENERGY, DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements, and your results will be noticeable much faster.



## COMPLEX STOP STRESS

### WHY IS STRESS DANGEROUS?



Stress is the organism's response to tension, negative emotions or excitement. During periods of severe or prolonged stress, the immune, endocrine, autonomic nervous systems suffer. Stressful situations are the cause of many physical and mental illnesses. The danger of stress is that its effects can accumulate gradually.



With severe or prolonged stress:  
The concentration of attention worsens, migraines appear;  
Problems with the cardiovascular system arise;  
Peptic ulcer and gastritis worsen;  
Immunity decreases  
In especially severe cases, degeneration of the cells of the spinal cord and brain occurs.



It is hard to work under stress and engage in habitual activities. Irritability, insomnia, and neurosis appear. Physical stamina is reduced. It becomes difficult to concentrate on the usual things – whether it is at work or at home. Relations with loved ones also suffer. An addiction to alcohol and/or drugs may appear.



Severe stressful situations require serious treatment – from psychotherapy sessions to taking medications that can be addictive. The way out of severe stress can last several years. It is much more effective to notice the first signs of stress in time or to prevent stressful situations.


## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Stress reduces the body's defenses, which is a favorable factor for infestation with parasites. They negatively affect the nervous system, causing fatigue, irritability and impair sleep quality. DeVita AP basic antiparasitic programs work against most known pathogens, quickly, efficiently and without side effects, removing the parasitic load.

"Drainage – General" program is used to cleanse the body after using basic antiparasitic programs, removing toxins and improving overall well-being.

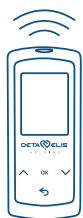

**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

The automatic program "Healthy Day" increases immunity and harmonizes the work of all body systems. Helps with acclimatization and stressful situations. The automatic "Complete Cleansing» helps in the restoration of the excretory system and in the harmonization of the body.

Special programs "Healthy sleep", "Good mood" restore the nervous system, relieve nervous tension, improve the process of falling asleep and the quality of sleep.



**DeVita Energy** relieves stress and chronic fatigue. Restores biorhythms. Launches the functions of regulating the body, helping a person to be in good shape, maintain vigor and good mood throughout the day.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating "junk food" (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- in case of problems with sleep - insomnia, worsening sleep quality;
- when you feel tired in the morning, lethargy and weakness during the day;
- with nervous strain, severe excitement, frequent manifestations of irritability, apathy.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir MULTI-ENERGY and DeLixir DETOX smart cell nutrition supplements.

## APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Good mood, Healthy sleep	If necessary, 2-3 hours before bedtime.
	4	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	5			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	6			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
	7			Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Good mood, Healthy sleep	If necessary, 2-3 hours before bedtime.
	4	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night.

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed for the night.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Good mood, Healthy sleep	If necessary, 2-3 hours before bedtime.
4	DeVita AP	Special programs	No nematodes and ascarides, No parasites - General, No parasites - Digestive track	Daily in the evening after the DeVita Ritm programs, alternating 1 program for 1-2 weeks.	

## STRESS PREVENTION MEASURES



Do not suppress spontaneous emotions. Suppressed anger or anxiety leads to stress.



Eat right. Add oranges and bananas, oatmeal, sea fish, herbal tea, avocado, nuts, fresh vegetables, dark chocolate to your diet.



Go to bed no later than 23:00 p.m. Sleep for at least 7-8 hours.



Do meditation or special breathing exercises.



Reduce your coffee, tea, and other caffeinated drinks.



Exercise. Exercising causes the production of endorphins - hormones of good mood, and also the adrenaline levels get higher.



Non-stressful monotonous work or creativity helps with stress - knitting, coloring pictures, etc.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Aggressiveness,  
fatigue, apathy,  
insomnia.

After first use  
of DeVita



Tranquility,  
mood improvement,  
burst of energy.

After a month of  
using DeVita



Positive mood,  
energy,  
poise.



Use the recovery scheme of the complex "STOP STRESS" together with the DeLixir MULTI-ENERGY and DeLixir pH BALANCE smart cell nutrition supplements, and your results will be noticeable much faster.





## COMPLEX STOP INSOMNIA

### WHAT DANGER DOES INSOMNIA HAVE?



The mechanism of sleep is not fully understood, but its primary importance has been proven. During sleep, recovery processes, the structuring of information received during the day, and physical relaxation occur. It is during this period that information is assimilated, tissues of the body are updated and restoration of immunity occurs.



Without a full night sleep, the body does not have time to recover and malfunctions. First of all, the nervous system suffers. Then irritability, aggressiveness, lethargy, apathy and reduced performance occur.

The mood drops sharply, a person becomes uncommunicative. This leads to numerous conflicts at work and at home.



Then physiological problems appear. A sharp decrease in immunity occurs, since immune cells - T-lymphocytes - are synthesized during sleep. Body weight increases due to the production of the stress hormone - cortisol. Chronic lack of sleep leads to frequent illnesses, fatigue and poor health.



The lack of a full and healthy sleep significantly affects the appearance. Collagen production decreases in the skin of the face, wrinkles appear, and elasticity decreases. The complexion acquires an earthy hue, swelling, dark circles and wrinkles appear under the eyes. Irritations or rashes may occur.


## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Microbes and larvae of large parasites can inhabit the organs and tissues of the body, including the nervous system. Their presence and the toxins secreted by them disrupt the system, leading to rapid fatigue, a state of nervous excitement and, subsequently, sleep disturbances. DeVita AP basic antiparasitic programs affect most known pathogenic organisms, quickly, efficiently and without side effects, removing parasitic load.

“Drainage - General” program is used to cleanse the body after applying basic antiparasitic programs, removing toxins and, as a result, reducing nervous system excitement and improving overall well-being.

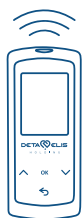

**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

“Healthy Day” automatic program improves the functioning of the excretory system, activates the circulatory and immunity systems. Helps in coping with physical, mental and psychological stress, creating favourable conditions for healthy sleep. The automatic program “Complete Cleansing” helps in the restoration of the functions of the excretory system and harmonizes the body.

“Healthy Sleep” special program reduces nervous excitement, promotes fast falling asleep and good quality sleep.



**DeVita Energy** relieves from psycho-emotional stress, stabilizes the nervous system, helps in restoring energy. It positively affects the quality of sleep, gives vigor when waking up.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water for at least 1-1.5 l / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- problems with falling asleep;
- sleep disorders (light sleep, waking up in the middle of the night, etc.);
- state of nervous excitement.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir MULTI-ENERGY and DeLixir DETOX smart cell nutrition supplements.

## APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Healthy sleep	If necessary, 2 hours before bedtime.
	4	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	5			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	6			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
	7			Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating every other day.
	3		Special programs	Healthy sleep	If necessary, 2 hours before bedtime.
4	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night.	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed for the night.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Healthy sleep	If necessary, 2 hours before bedtime.
4	DeVita AP	Special programs	No nematodes and ascarides, No parasites- General, No parasites- Digestive track	Daily in the evening after the DeVita Ritm programs, alternating 1 program for 1-2 weeks.	

## PREVENTIVE MEASURES FOR SLEEP DISORDERS



Maintain a comfortable temperature for sleeping in the bedroom - 19-20 ° C. Walk before bedtime and ventilate the room, and in the warm period, sleep with the window open.



Create a ritual of going to bed: take a hot bath, listen to music, relax. Try not to do work or daily routine in the bedroom.



Try not to sleep during the day. Or at least reduce your daytime sleep to 20-30 minutes.



The key to a good sleep is complete silence and darkness. Use blackout curtains and for protection against extraneous sounds - ear plugs.



A way of falling asleep easier is by using "white noise" sounds - monotonous sound of the radio on an "empty" frequency. Other examples are the sounds of waterfall, fan, wind, rain etc.



Exercise and eat no later than 2 hours before bedtime.



Remove flowers and other plants from the bedroom: they produce oxygen only in the daytime, and at night, on the contrary, they absorb it.



If possible, use blue and purple tones in the design of the bedroom.



Include in your evening menu products with a calming effect: honey, milk, bananas, pumpkin, date fruits.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements

Before  
use



Difficulty in falling asleep, drowsiness, restlessness, apathy.

After the first use  
of DeVita



Emotional and energetic recovery, improvement of general condition.

After a month  
of using DeVita



Cheerfulness, improved efficiency, easy awakening, improved communication.



Use the scheme of the "STOP INSOMNIA" complex together with the DeLixir MULTI-ENERGY and DeLixir DETOX smart cell nutrition supplements, and your results will become noticeable much faster.



## COMPLEX WOMEN'S HEALTH

### WHY ARE WOMEN'S HEALTH PROBLEMS DANGEROUS?



Women are more adapted to long-term stress than men. But the female body has its own resources, and in case of serious prolonged loads, female health is weakened. Ecology has a huge impact on this matter, as every 5th female resident of any city is faced with «female» problems.



The harbingers of the failure of the female body are: fatigue, bad mood, irritability, constant prolonged fatigue, unwillingness to get up in the morning and loss of interest in intimacy with a man. If you ignore these changes, they can lead to serious health problems.



Failures in the female genitourinary system lead not only to disruption of the cycle and discomfort, but also to problems with pregnancy, inflammation and even tumors. Common female problems include an irregular cycle, cervical erosion, uterine inflammation and inflammation of fallopian tubes, cystitis, bacterial and fungal infections.



Problems with female health lead to negative changes in appearance. The skin becomes more oily and body hair growth is increased. Inflammation and acne appear. Failure in the production of female hormones can lead to rapid weight gain and hair loss.


## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibratory frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Parasites and microbes cause inflammation of the organs of the female genitourinary system and even lead to infertility and chronic kidney diseases. DeVita AP basic antiparasitic programs affect most known pathogens, quickly, effectively and without side effects, removing the parasitic load from the organism.

“Drainage - General” automatic program is used to cleanse the body after applying basic antiparasitic programs, remove toxins, and has an antiseptic effect.

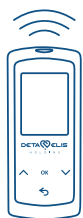

**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

The automatic program “Healthy Day” harmonizes the work of the body, the energy of the chakras, activates the circulatory system, digestion processes, the excretory system, normalizes the functions of the genitourinary system and improves the quality of intimate life. The automatic program “Complete Cleansing” contributes to the restoration of the main excretory system and the harmonization of the body.

Special programs “Female sexuality”, “Women’s health”, harmonize sexual desire, regulate the reproductive system, facilitate the condition of women during menstruation and menopause.



**DeVita Energy** relieves from psycho-emotional stress, restores and strengthens physiological indicators. It gives energy and vitality, fights against sexually transmitted diseases, removes psychological burden.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- violation of the menstrual cycle;
- severe menstrual pain;
- to harmonize the work of the body during menopause;
- decrease in sexual desire.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir MULTI-ENERGY, DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements.

## APPLICATION SCHEMES









SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Women's health	Daily in the afternoon.
	4	Female sexuality		If necessary, according to guidelines.	
	5	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	6			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	7			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
	8			Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Women's health	Daily in the afternoon.
	4	Female sexuality		If necessary, according to guidelines.	
	5	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night.
6	Special programs		No chlamydia	Daily in the evening after DeVita Ritm programs.	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed, at night.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Women's health	Daily in the afternoon.
	4			Female sexuality	If necessary, according to guidelines.
5	DeVita AP	Special programs	No nematodes and ascarides, No parasites- General, No parasites- Digestive track	Daily in the evening after the DeVita Ritm programs, alternating 1 program every 1-2 weeks.	

## PREVENTIVE MEASURES FOR DISEASES OF THE GENITOURINARY SYSTEM

-  Visit your gynecologist and urologist at least once per year.
-  Observe hygiene. Use special cleansers - regular soap disrupts the microflora of the genitourinary system.
-  Take hormonal drugs and contraceptives only if they are prescribed by your doctor. Do not take medicines, especially antibiotics, for longer than 7 days without special instructions.
-  Have a regular sex life, be sure to use means of protection against sexually transmitted diseases.
-  Protect the genitourinary system from excessive heat or cold.
-  Do not lift weights and control physical activity to avoid prolapse of the uterus and other problems.
-  Walk at least 20 minutes a day to avoid stagnation in your pelvic organs.
-  Exclude sports and sexual intercourse during menstruation.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Bad mood, decreased libido, discomfort and cycle disruptions.

After the first use  
of DeVita



Burst of energy, decrease of pain, elimination of irritability and apathy.

After a month  
of using DeVita



Improved mood, normalization of the reproductive system, increased libido.



Use the scheme of the "WOMEN'S HEALTH" complex scheme together with the DeLixir MULTI-ENERGY, DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements, and your results will be noticeable much faster.





## COMPLEX MEN'S HEALTH

### WHY ARE MALE HEALTH PROBLEMS DANGEROUS?



Even young men may have problems with erection due to unhealthy lifestyle and urological diseases, the most common of which is prostatitis. Another dangerous disease is prostate adenoma, a benign tumor of the prostate gland.

Every second male over age of 50, suffers from prostate adenoma.



Acute prostatitis is difficult, with increased body temperature up to 39-40 ° C. Urination is difficult and very painful. If untreated, prostatitis becomes chronic. It can affect the bladder and kidneys, cause inflammation of the prostate and sepsis. And the source of infection can not always be determined.



Prostate adenoma occurs due to a combination of reasons - age, hormonal disorders, obesity, stress, etc. Prostate adenoma is dangerous and complications cause: acute urinary retention, urinary incontinence, problems in sexual life, and infections of the genitourinary system. Adenoma often requires surgical treatment.




Due to prostatitis and prostate adenoma, not only the urination process is disturbed, but also the potency. Also, the complications of these diseases can lead to infertility. Often men are faced with psychological disorders - this leads to depression, neurosis, chronic stress and other neurological diseases.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Parasites cause inflammation of the genitourinary system. Diseases proceed unnoticed by humans, having a serious impact on the function of the urinary and reproductive systems. DeVita AP basic antiparasitic programs affect most known pathogens, quickly, effectively and without side effects, removing the parasitic load from the organism.

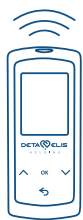

“Drainage - General” automatic program is used to cleanse the body after applying basic antiparasitic programs, remove toxins, and has an antiseptic effect.

**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

“Healthy Day” automatic program harmonizes the work of the organism, normalizes the functions of the genitourinary system and improves the quality of intimate life. The automatic program “Complete Cleansing” helps in main excretory systems restoration and harmonizes the organism.



Special programs “Male sexuality”, “Men’s Health”, increase self-confidence and sexual desire, regulate sexual activity, ensure the normal functioning of the genitourinary system.

**DeVita Energy** relieves from psycho-emotional stress, restores and strengthens physiological indicators. Gives energy, reduces fatigue accumulated throughout the day.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the device. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- violation of potency;
- psychological problems associated with the “neurosis of expectation of failure”;
- problems with the bladder and prostate.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir MULTI-ENERGY, DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements.

## APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Men's health	Daily in the afternoon.
	4	Male sexuality		If necessary, according to guidelines.	
	5	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	6			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	7			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
	8			Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Men's health	Daily in the afternoon.
	4	Male sexuality		If necessary, according to guidelines.	
	5	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night.
	6		Special programs	No chlamydia	Daily in the evening after DeVita Ritm programs.

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed, at night.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Men's health	Daily in the afternoon.
		Male sexuality		If necessary, according to guidelines.	
4	DeVita AP	Special programs	No nematodes and ascarides, No parasites- General, No parasites- Digestive track	Daily in the evening after the DeVita Ritm programs, alternating 1 program every 1-2 weeks.	

## PREVENTIVE MEASURES FOR UROGENITAL PROBLEMS



Maintain a regular and stable sex life. Excessive activity and prolonged abstinence lead to disruption of the genitourinary system.



Exercise at least 20 minutes a day.



Avoid hypothermia, especially in the groin area.



When changing your regular partner, be sure to use condoms. Sexually transmitted infections can cause prostatitis.



Keep a track of your bowel movements. Frequent constipation leads to prostatitis..



4 times a year, have a massage of the lower back area.



Have a routine examination by a urologist 2 times a year.



Diversify your diet with fresh fruits, vegetables, herbs, seafood, nuts and seeds.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Pain in the lower abdomen,  
frequent urges to urinate,  
sexual problems.

After the first use  
of DeVita



Pain relief,  
mood improvement,  
surge of energy.

After a month of  
using DeVita



Stable genitourinary system  
functioning, stable erection,  
self-confidence.



Use the scheme of the "MEN'S HEALTH" complex in conjunction with DeLixir MULTI-ENERGY, DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements, and your results will be noticeable much faster.



## COMPLEX MOVEMENT WITHOUT PAIN

### WHY SPINE AND JOINT PROBLEMS ARE DANGEROUS?



Problems with spine and joints are a fee that's humans have to pay for upright posture. Having only two points of support, the human body is forced to constantly maintain balance. Strictly speaking, each person is fated with problems with the spine. Because of this, the joints experience increased stress.



We exacerbate these problems by leading a sedentary lifestyle, lifting weights, or not paying attention to posture. A huge role is played by shoes, which are designed to soften the load on the spine, but do not always fulfill this function. Meanwhile, over 70% of people have spinal curvature.



The most common problems are impaired mobility, severe pain and inflammation. And this indirectly affects the appearance, causing obesity, varicose veins, problems with gastrointestinal tract and general fatigue. In especially severe cases, problems with the spine and joints lead to restricted movement and disability.



Diseases of the back and joints can go unnoticed for a long time, but at the most unexpected moment declare about themselves with severe pain. Treatment can, basically, only relieve symptoms - reduce pain, improve mobility and blood circulation. Therefore, it is important to take care of the health of the spine and joints on time.

## MECHANISM OF ACTION




**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

The parasite larvae that enter the body settle on the membrane lining the joint cavity, or in the periarticular tissues, causing inflammation. Parasites also affect the spine. DeVita AP basic antiparasitic programs affect most known pathogens, quickly, effectively and without side effects, removing the parasitic load from the organism.

“Drainage - General” automatic program is used to cleanse the body after applying basic antiparasitic programs, reduce the parasitic intoxication and improve overall well-being.

Special programs “No bacteria” and “Drainage Antitox” help to reduce toxic load in the presence of specific groups of microorganisms and salts of heavy metals in the body that have a negative effect on the musculoskeletal system.

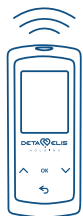

**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

“Healthy Day” automatic program improves blood circulation, gently reduces the toxic load, which enhances immunity, improves the outflow of metabolic products from damaged areas. “Complete Cleansing” automatic program helps in main excretory system restoration and harmonizes the organism.

The programs “Recovery from injuries”, “Joints without inflammation”, “Joints without pain” relieve from inflammation, pain, swelling, activate metabolic processes in damaged tissues, normalize muscle tone, improve the functioning of the ligamentous apparatus.



**DeVita Energy** triggers self-regulation mechanisms, restores physiological processes in organs and tissues. Gives vigor and energy for movement.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- pain in the joints and spine;
- injuries and inflammation of the muscles and joints;
- pinched nerves.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir MULTI-ENERGY, DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements.

## APPLICATION SCHEMES









SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Joints without pain and / or Joints without inflammation and / or Recovery from injuries	Daily in the afternoon, if necessary.
	4	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	5			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	6			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
	7			Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Joints without pain and / or Joints without inflammation and / or Recovery from injuries	Daily in the afternoon, if necessary.
	4	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night
5	Special programs		No bacteria, Drainage Antitox	Daily in the evening after DeVita Ritm programs for 14 days.	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed, at night.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Joints without pain and / or Joints without inflammation and / or Recovery from injuries	Daily in the afternoon, if necessary.
4	DeVita AP	Special programs	No nematodes and ascarides, No parasites- General, No parasites- Digestive track	Daily in the evening after DeVita Ritm programs, alternating 1 program every 1-2 weeks.	

## PREVENTIVE MEASURES FOR SPINAL AND JOINT PROBLEMS

-  Use vitamins and minerals in spring and autumn.
-  Watch your posture. Lean your back against a flat wall - your nape, shoulder blades, buttocks and feet should touch the wall. If necessary, consult a doctor and purchase a corrective corset.
-  Sleep on an orthopedic mattress and pillow. In extreme cases - on a dense hard surface. Fall asleep lying on your back.
-  Exercise and strengthen your back muscles and abs. Developed abdominal muscles relieve pressure on the lower back.
-  When doing sports and lifting weights, follow the correct technique - do not round your back, strain your abs.
-  Eat sour-milk products rich in calcium, as well as products with gelatin content - meat with gelatin and jelly.
-  Consult a doctor if you want to try a new physical activity - yoga, step aerobics, strength training, etc.
-  Be wary of osteopathy, manual therapy, mechanical massage.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Pain in the back and joints,  
fatigue, headache.

After the first use  
of DeVita



Decreased pain  
and tiredness.

After a month of  
using DeVita



Fast recovery,  
ease of motion,  
no pain.



Use the scheme of the "MOVEMENT WITHOUT PAIN" complex together with the DeLixir MULTI-ENERGY, DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements, and your results will be noticeable much faster.





## COMPLEX CLEAN SKIN

### WHAT IS CAUSING SKIN PROBLEMS?



The condition of the skin is an indicator of human health. Failures in the work of some body systems: digestive, endocrine, immune, circulatory, are reflected on the skin of the face and body. Also, the condition of the skin and hair is affected by ecology, stress, the use of low-quality cosmetics.



The most common facial skin problems are acne and blackheads, age spots, and rosacea. Basically, these problems have age-related causes associated with hormone imbalance. However, facial skin conditions may worsen at an earlier age. The reason for this may be improper skin care or the selection of poor-quality cosmetics.



Common scalp problems are dandruff and seborrhea, which cause excessive oiliness or dryness of the skin. This affects the health of the hair - they become thin, brittle and weak, begin to fall out. Problems with the scalp can be dermatological, they are also often caused by microorganisms – bacteria and fungi.



Body skin problems are similar to scalp problems. Those are dryness in the winter season, acne and blackheads, especially in the back, low skin tone and sagging. These problems are caused by both external pathogens and internal malfunctions in the body, as well as improper care, unbalanced diet and daily routine.


## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Viral, bacterial, fungal and other types of infections cause problems with the skin, hair, and nails. Dandruff, the formation of crusts, scales, pustules, allergic skin reactions often indicate the presence of pathogenic organisms. The DeVita AP basic antiparasitic programs affect the majority of known pathogens, quickly, efficiently and without side effects, removing the parasitic load and signs of inflammation.

“Drainage - General” automatic program is used to cleanse the body after applying basic antiparasitic programs, reduces parasitic intoxication and improves both overall well-being and skin condition.

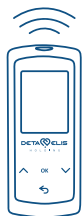

**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

“Healthy Day” automatic program improves blood circulation, normalizes the functioning of the gastrointestinal tract. Gently reduces toxic load, improving the nutrition of skin cells and hair follicles, normalizing cellular metabolism and activating regeneration processes. The automatic program “Complete Cleansing” contributes to the restoration and harmonization of the body.

The programs «Beautiful Skin», «Beautiful Hair and Nails» relieve from inflammation, activate metabolic processes, improve skin tone, strengthen hair and nails, slow down the aging process.



**DeVita Energy** launches self-regulation mechanisms, restores physiological processes in organs and tissues, enhancing the action of the DeVita AP and DeVita Ritm programs.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) – this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- inflammatory processes on the skin, acne;
- deterioration of hair and nails;
- hair loss, dandruff.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements.

## APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Beautiful skin and / or Beautiful hair and nails	Daily in the afternoon, if necessary.
	4	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	5			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	6			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
	7			Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Beautiful skin and / or Beautiful hair and nails	Daily in the afternoon, if necessary.
	4	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night.
5	Special programs		No fungi, No herpes zoster, No herpes simplex	If necessary, according to guidelines.	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed, at night.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Beautiful skin and / or Beautiful hair and nails	Daily in the afternoon, if necessary.
4	DeVita AP	Special programs	No nematodes and ascarides, No parasites- General, No parasites- Digestive track	Daily in the evening after DeVita Ritm programs, alternating each program every 1-2 weeks.	

## PREVENTIVE MEASURES FOR PROBLEMS WITH FACE, BODY AND HEAD SKIN



With pronounced problems that arise for no apparent reason, conduct an examination of the body for internal pathologies.



Do not use other people's cosmetic and hygiene items, do not try on other people's hats.



Try to avoid prolonged sun exposure, limit your visits to the solarium. Use cosmetics with UV filters.



Monitor the balance of hormones, with a predisposition to hormonal disruptions, undergo an examination at least once a year.



Use only high-quality cosmetics, monitor their expiration date. Thoroughly clean your face with a special cleanser in the morning and evening.



Do not perform scrubs and mechanical cleanings frequently. Frequent exfoliation causes irritation and dryness.



In the cold period, be sure to moisturize the skin after a shower, use thermal water for the face every few hours.



Red fish, eggs, greens, nuts and legumes are useful for good skin and hair condition.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements

Before  
use



Acne, dandruff,  
age spots,  
sagging skin.

After the first use  
of DeVita



Improved skin tone,  
reduced rashes,  
regulated oil balance.

After a month of  
using DeVita



Dark spots and inflammation  
becomes less noticeable,  
dandruff disappears.



Use the scheme of the "CLEAN SKIN" complex with the DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements, and your results will be noticeable much faster.



## COMPLEX SPORT PERFORMERS

### HOW DO TRAINING AND SPORTS AFFECT THE ORGANISM?



Physical activity seriously changes the metabolism and vital functions of the body. Athletes need more nutrients and vitamins, relaxation and full sleep. Bones and joints experience increased stress. Metabolism is accelerated, the need for fluid intake increases.



When doing fitness or sports, we want to get high results in a short time. However, the body's resources are limited. Due to active metabolism during training, a large expenditure of energy occurs, the accumulation of decay products, which are removed only after 24 hours, increases the toxic load on the liver and kidneys.



Intense and active workouts cause post-workout muscle pain and a feeling of fatigue. In addition, injuries are common in sports. For best results, you need to resume training only after a full recovery, but this is fraught with a significant loss of time. And in sports, every day training is important.



Each training or competition is a great stress for the body. The accumulated stress negatively affects the general state of health and can reduce working capacity, and in the future - cause heart problems. The success of training is influenced by psychological factors - mood, personal relationships.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Parasites and microorganisms can significantly undermine health. This affects the quality of training and recovery. DeVita AP basic antiparasitic programs affect the majority of known pathogens, quickly, efficiently and without side effects, removing the parasitic load and signs of inflammation.

"Drainage - General" automatic program is used to cleanse the body after applying basic antiparasitic programs, reduces parasitic intoxication and improves overall well-being.



**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

"Healthy Day" automatic program improves blood circulation, facilitating quick recovery after training. The work of the digestive tract is normalized, due to which the load on the liver is reduced, the elimination of toxins is accelerated. "Complete Cleansing" automatic program helps to restore and harmonize the body.



The "Warm-up" and "Workout" programs help to increase strength indicators, allow you to train longer and more intensively, give energy and vigor. The programs "Recovery from injuries", "Joints without inflammation", "Joints without pain" relieve from inflammation and pain, normalize muscle tone, improve the functioning of the ligamentous apparatus, and help recover after an injury.

**DeVita Energy** activates the mechanisms of self-regulation, restores the physiological processes in organs and tissues. Due to this, recovery is much faster.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating "junk food" (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- intense physical exertion or sports training;
- pain in the joints, muscles and spine after exertion;
- injuries and inflammation of the organs of the musculoskeletal system;
- for the best recovery after training.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir MULTI-ENERGY, DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements.

## APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Warm-up, Workout	Accordingly, before and during training.
				Joints without pain and / or Joints without inflammation and / or Recovery from injuries	Daily in the afternoon, if necessary.
	4	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	5			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	6			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
7	Mode 4 automatic			Daily at night during the 4 <sup>th</sup> week.	

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic	Daily in the morning.
	3		Special programs	Warm-up, Workout	Accordingly, before and during training.
4	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night.	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed, at night.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Warm-up, Workout	Accordingly, before and during training.
	4			Joints without pain and / or Joints without inflammation and / or Recovery from injuries	Daily in the afternoon, if necessary.
5	DeVita AP	Special programs	No nematodes and ascarides, No parasites- General, No parasites- Digestive track	Daily in the evening after DeVita Ritm programs, alternating each program every 1-2 weeks.	

## PREVENTIVE MEASURES AND RECOMMENDATIONS FOR EFFECTIVE TRAINING



Start your workout with more complex exercises for large muscle groups, and end with isolated or cardio exercises.



Watch your pulse. The maximum heart rate is calculated from the formula  $214 \text{ minus age times } 0.8$ . The pulse during cardio training should be 50-59% of the maximum, with light fitness - 60-69%, during endurance training - 70-79%, during anaerobic training - 80-89%.



Eat right. In terms of quality composition, proper nutrition for an athlete should be close to the formula: 30% - proteins, 60% - carbohydrates, 10% - fats.



Use natural energy sources - ginseng extract, guarana, succinic acid.



Variety of exercises and training methods - muscles quickly get used to uniformity.



If possible, train at the same time. Choose it depending on your daily routine and activity level.



If you have severe muscle pain, take a hot bath and do self-massage.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements

Before  
use



Fatigue, muscle tension  
and pain,  
long recovery time.

After the first use  
of DeVita



Burst of energy,  
reduction of pain.

After a month of  
using DeVita



Increased efficiency,  
quick recovery, improved  
stamina and strength.



Use the scheme of the "SPORT PERFORMERS" complex together with the DeLixir MULTI-ENERGY, DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements, and your results will be noticeable much faster.





## COMPLEX STOP ALLERGY

### HOW DANGEROUS ARE ALLERGY MANIFESTATIONS?



Allergy is a disrupted mechanism for the implementation of the body's immune response. The body mistakenly takes a substance as a threat to life and forces the immune system to produce antibodies and various chemicals, one of which is histamine, which causes allergy symptoms.



In addition to food, wool and pollen, the role of an allergen can be played by chemical substances, some metals, alcohol, cold and even sunlight. There is also an allergy to substances formed inside the body itself. Often the symptoms are similar to the common cold or usual fatigue, making diagnosis difficult.



Allergic reactions may be the cause of development of bronchial asthma, hay fever, allergic dermatitis. The most serious complications are convulsions, acute vascular insufficiency, stomach cramps. The most severe allergic reactions - anaphylactic shock and Quincke's edema - lead to instant death.



According to statistics from the World Health Organization, about 40% of the world's population suffer from allergies. Diagnosing and identifying allergens is a complex and time-consuming process, and treatment is reduced to reduce symptoms and alleviate the general condition of the body, but does not eliminate the allergy in general.

## MECHANISM OF ACTION




**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Bacteria, fungi, viruses, helminths are biological allergens. The basic antiparasitic programs of DeVita AP affect most microorganisms, quickly, efficiently and without side effects, removing the pathogenic load.

Parasitic burdening of the liver and blood reduces their functions, enhances intoxication and allergic reactions. Special programs “No schistosomes (blood-flukes) & trematods”, “No allergies”, “No herpes” relieve the pathogenic load.

The program “Drainage – General” is used to cleanse the body after applying antiparasitic programs, the removal of toxins, which can also be allergens.

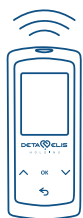

**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

“Healthy Day” automatic program has a beneficial effect on the work of the whole organism, increasing its resistance to foreign agents. Strengthens the immune system, increases vitality. Special programs “Without allergens” and “Anti-cold” relieve from itching, reduce lacrimation, swelling of the mucous membranes, and activate immune defense in case of allergies. The automatic program “Complete Cleansing” helps in the restoration of the basic functions of the excretory system and in the harmonization of the body.



**DeVita Energy** harmonizes the work of the body, energizes, activates the process of elimination of toxins. Improves the general condition when allergic reactions manifest, reduces lethargy and drowsiness arising when taking antihistamines.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks – they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) – this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

Allergy is a consequence of intoxication of the body. To minimize its manifestations, it is recommended to carry out cleansing measures within 6 months. In the first three months it is necessary to use the devices according to schemes 1-3, in the next 3 months – alternate schemes 2 and 3. When cleansing and restoring the body, minor exacerbations are possible, in this case, a break in the use of devices is necessary.

In acute and severe conditions, special programs can be included several times a day.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices together with DeLixir DETOX and DeLixir pH BALANCE smart cell nutrition supplements (if there is no food allergy).

## APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Without allergens	Daily in the morning, alternating with each other every other day.
	4	DeVita AP	Auto-programs	Anti-Cold, Free breath	If necessary, according to guidelines.
	5			Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	6			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	7			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
	8			Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.
	9		No allergies	In the evenings after the DeVita Ritm programs.	

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Without allergens	Daily afternoon.
	4	DeVita AP	Auto-programs	Anti-Cold, Free breath	If necessary, according to guidelines.
	5			Cleansing-automatic	Daily at night.
	6			Special programs	No allergies
7			No schistosomes (blood-flukes) & trematods	In the evenings after the DeVita Ritm programs, 14 days.	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If necessary, at night.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Without allergens	Daily afternoon.
	4	DeVita AP	Special programs	Anti-Cold, Free breath	If necessary, according to guidelines.
	5			No allergies	In the evenings after the DeVita Ritm programs.
6	No nematodes and ascarides, No parasites- General, No parasites- Digestive track.			Daily in the evening after the DeVita Ritm programs, alternating 1 program for 1-2 weeks.	

## PREVENTION MEASURES FOR ALLERGIC REACTIONS



If you already have an allergy to any substance, go through a full allergen test. Often an allergy-prone organism cannot tolerate several foods.



Check the gastrointestinal tract. The presence of helminths, dysbiosis and diarrhea problems exacerbate the course of allergic reactions.



If you are allergic to dust, store bedding, pillows, and books in tightly closed cabinets. Wash clothes weekly at a temperature of at least 60 °C.



When hay fever starts, perform frequent wet cleaning of your personal and living space every day. Take off your outdoor shoes and clothes at the doorstep - pollen and plant seeds accumulate on them.



For food allergies, avoid foods with strong allergenic properties: honey, peanuts, fish and seafood, spices, chocolate, citrus fruits, tomatoes, strawberries.



Use the medicine only as directed by your doctor. Keep track of the expiration date of drugs - expired drugs emit volatile allergenic substances.



If you are allergic to mold, do not close the doors to the bathroom and the kitchen. Throw away moldy foods. Check the flower pots - mold may form in them.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Runny nose, itching, swelling, skin rashes, lacrimation.

After first use  
of DeVita



Reduced itching, relieve from swelling, better general condition.

After a month  
of using DeVita




Reduced allergy manifestation, reduced rashes, breathing relief.



Use the recovery scheme of the "STOP ALLERGY" complex together with DeLixir DETOX and DeLixir pH BALANCE smart cell nutrition supplements (if there is no food allergy), and your results will become noticeable much faster.



# PARTICULAR COMPLEXES

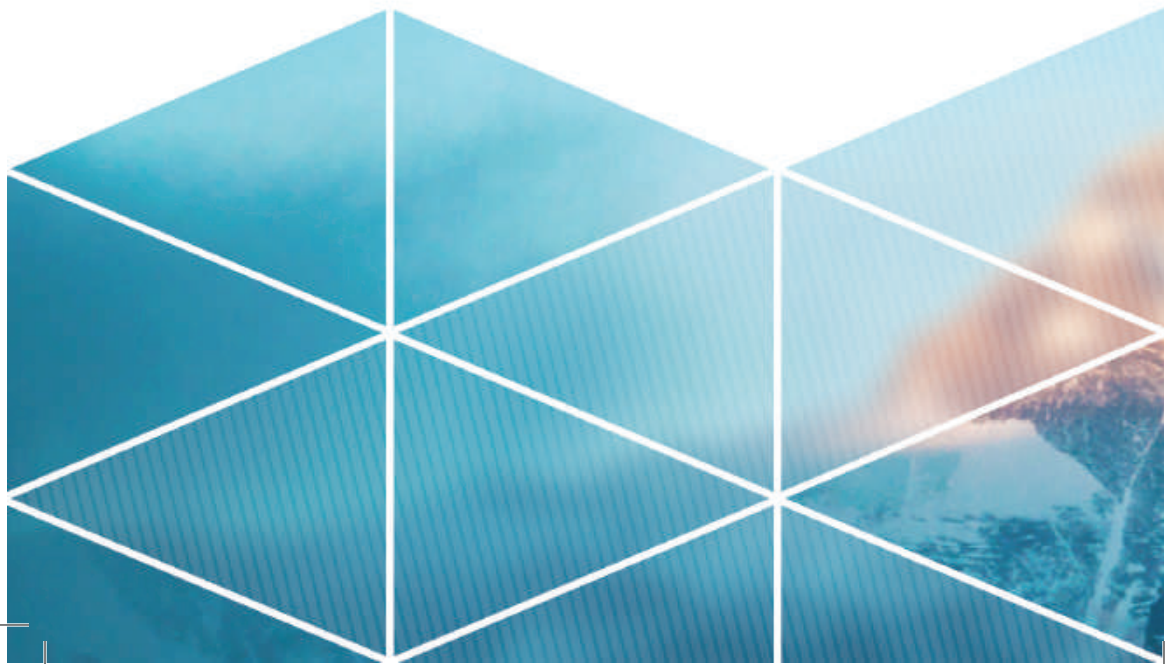



## COMPLEX WITHOUT COLD

### WHY ARE DIGESTIVE DISEASES DANGEROUS?



The immune system protects the body, participates in its life activity and in adaptation processes to changing conditions.

Immune cells should always be on guard, because the body is attacked by pathogenic microbes every second.

To some of their species, long-term immunity is developed, but not to colds and flu viruses.



Frequent catarrhal diseases lead to disorders

in the work of organs and an even greater decrease in immunity. Microbes freely enter the body and destroy it.

As a result, the common cold results in serious complications. They require outpatient treatment - from 2 weeks to several months of bed rest.



The most common complications are inflammations:

sinuses of the nose (sinusitis, frontal sinusitis);

middle ear (otitis media, labyrinthitis);

respiratory organs (bronchitis, pneumonia, pleurisy).

Every complication affects major organ systems and causes severe intoxication.



Medicines against flu and colds also cause negative reactions.

Antipyretic drugs reduce immunity and contribute to inflammation. Vasoconstrictor nasal drops affect the mucous membrane and small vessels, causing chronic edema. Antitussive drugs produce sputum and respiratory infections.


## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

DeVita AP's basic antiparasitic programs affect most known pathogens. The use of the programs helps in coping with the disease quickly, efficiently and without side effects. In case of a long and severe course of the disease, combine DeVita AP programs from schemes 2 and 3.

The "Drainage – General" program is used to cleanse the body after basic antiparasitic programs, remove toxins, strengthen immunity and recover from the disease. As a result, overall well-being improves.

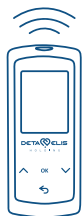

**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

Auto program "Healthy Day" improves local and general blood circulation, the functioning of the digestive system, which is important for an imperfect diet. Helps increase immunity and harmonize the work of all systems. Helps with acclimatization and stressful situations. The special program "Free Breath" restores the bronchopulmonary system and reduces coughing. A special program «Anti-cold» relieves inflammation of the upper respiratory tract.



**DeVita Energy** relieves stress and chronic fatigue. It neutralizes the effects of negative environmental factors and domestic radiation, helps with weather sensitivity. Launches the body's regulation functions, helping a person to be in good shape, maintain vigor and good mood throughout the day.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - the impact of the devices. Remove or minimize coffee in your diet. Avoid eating "junk food" (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir MULTI-ENERGY and DeLixir pH BALANCE smart cell nutrition supplements.



## APPLICATION SCHEMES



USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

SCHEME №1: in the season of colds and flu, in close contact with sick people, with overcooling

SCHEME №1	DEVICE		PROGRAM	USAGE
	1	DeVita Energy	Mode program	1-2 times in the morning.
	2	DeVita Ritm	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3	DeVita AP	Antivirus + Without Toxins	1-2 times per week.

SCHEME №2 : at the first sign of bacterial and viral infections (SARS, ARI, tonsillitis, rhinitis, etc.), except for flu.

SCHEME №2	DEVICE		PROGRAM	USAGE
	1	DeVita Energy	Mode program	1-2 times in the morning.
	2	DeVita Ritm	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Anti-Cold, Free breath	Daily afternoon, 5-7 days.
	4	DeVita AP	Antivirus + Without Toxins	First 2 days, after DeVita Ritm programs.
	5		Mode 2 automatic + Anti-Cold	Daily, alternating with each other for 5-7 days.

SCHEME №3: at the first signs of flu.

SCHEME №3	DEVICE		PROGRAM	USAGE
	1	DeVita Energy	Mode program	1-2 times in the morning.
	2	DeVita Ritm	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Anti-Cold, Free breath	Daily afternoon, 5-7 days.
	4	DeVita AP	Antivirus + Without Toxins	First 2 days, after DeVita Ritm programs.
	5		Anti-Flu + Against flu - General + Without toxins	Daily, 5-7 days.
6	Mode 2 automatic		A week after influenza programs, daily for 5-7 days.	

## PREVENTIVE MEASURES DURING THE EPIDEMIC OF FLU AND COLD



Perform frequent wet cleaning of your personal and living space.



Ventilate the room every day.



While being outside, keep your legs, head and neck warm and dry.



Change your disposable hygiene mask every 3 hours.



Use disposable wipes to cover your nose and mouth when you sneeze or cough, and throw them away immediately after use.



Eat garlic, red bell pepper, cranberries, rose hips, carrots, sea fish and seafood, yogurt without additives, kefir, honey, broccoli.



In addition to juices, tea and other liquid products, be sure to drink plenty of clean water at least 1-1.5 l / day.



Moisten the mucous membranes with saline (1 tsp salt in 1 liter of water). Spray or drip it into your nose every 1-2 hours.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Fatigue in the morning, complications after illness, lethargy and weakness during the day.

After first use  
of DeVita



A burst of energy, improvement in well-being and mood.

After a month  
of using DeVita



Easy rise in the morning, vigor, increased adaptive reserves of the body, reduced morbidity.



Use the recovery scheme of the complex "WITHOUT COLD" together with the DeLixir MULTI-ENERGY and DeLixir pH BALANCE smart cell nutrition supplements, and your results will be noticeable much faster.



## COMPLEX STOP CORONAVIRUS- PREVENTION

### BRIEFLY ABOUT CORONAVIRUS



Coronaviruses are a family of viruses that includes 40 types of RNA-viruses from the group of coronaviruses of animal origin, that infect humans and animals. The name is associated with the structure of the virus, the structures on the virus itself, which resemble a crown, thanks to which the virus enters the cell of a person or animal.



Sources of coronavirus infections can be an infected person or animal. The main transmission routes: airborne, dust, fecal-oral, contact. The infection rate increases during winter and spring.



According to WHO, the virus spreads only within a radius of 1 meter around an infected person. In air, the virus retains its activity for 8-10 hours, in water – up to 9 days. In the most favorable conditions, but outside the human or animal body – 48 hours.



The most effective means to prevent infection are to treat your hands with alcohol-containing antiseptics and thoroughly wash your hands with soap. Do not touch with dirty hands your eyes, nose, and mouth. Masks are ineffective because they get wet from breathing and the virus penetrates through the moisture in the mask.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist discoverer Dr. Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency; it is the increase in amplitude which leads to the death of the microorganism without harming other cells.

DeVita AP antimicrobial programs help the body's efforts to destroy respiratory viruses and bacteria. The programs quickly and effectively help to cope with the infection, as well as help the body recover from the disease.

Special programs "Without cough" and "Without pneumonia" are used to reduce inflammation and facilitate the discharge of sputum from the respiratory tract.



**DeVita Ritm** – here the action is based on the studies of the scientists Dr. Reinhold Voll and A.F. Cramer, which found that each organ has a healthy frequency of oscillation. DeVita Ritm gently affects the body, pointing out the frequency of the vibrations of healthy organs. As a result, the body is capable of self-healing.

The automatic programs "Healthy Day - automatic" and "Health - automatic" have a beneficial effect on the work of the whole organism, increasing its resistance to foreign agents. They strengthen the defenses, increase vitality. Using the "Anti cold - automatic" automatic program helps to restore immunity and increase the energy potential of the respiratory system.



**DeVita Energy** harmonizes the function of the body, energizes, activates the removal of toxins, starts the processes of self-regulation.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING THE DEVITA DEVICES

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the effectiveness of the devices. Remove or minimize coffee in your diet.

Avoid eating "fast food" (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING CASES:

- if you have to communicate with sick people;
- if your work is connected with a large number of people;
- if you have been in countries with an unfavorable situation for this disease;
- if you often get colds;
- if you want to protect yourself from infection.



### COMBINE THE DEVICES WITH THE DELIXIR SMART CELL NUTRITION SUPPLEMENTS

For best results, it is recommended to use the DeVita devices in conjunction with the smart cell nutrition supplements DeLixir pH Balance, DeLixir Top Form, DeLixir Multi-Energy.

# APPLICATION SCHEMES

SCHEME №1 For infection prevention:

SCHEMA №1	DEVICE	PROGRAM		APPLICATION
	1 DeVita Energy	Basic program		1-2 times, in the morning.
	2 DeVita Ritm	Automatic programs	Healthy Day - automatic	1 <sup>st</sup> day, in the morning.
			Health - automatic	2 <sup>nd</sup> day, in the morning.
			Immune protection	Use after the automatic program "Healthy Day – automatic" on the DeVita Ritm device.
	4 DeVita AP	Manual programs	Without cold - automatic	2-3 times a week, at night.
			Stop Coronavirus- Prevention- New	2-3 times a week after the use of the DeVita Ritm device.
	7 DeLixir		Multi-Energy	15 ml, in the morning, after breakfast.
	8 DeLixir		pH Balance	15 ml, in the evening, after dinner.

SCHEME № 2 During the infection, at its first symptoms:

SCHEMA №2	DEVICE	PROGRAM		APPLICATION		
	1 DeVita Energy	Basic program		1-2 times, in the morning.		
	2 DeVita Ritm	Automatic programs	Healthy Day - automatic	Alternate every second day for 2 weeks. Use in the morning.		
			Health - automatic			
			4 DeVita Ritm	Automatic programs	Anti cold - automatic	Use after the automatic program "Healthy Day – automatic" on the DeVita Ritm device (first week) and after the automatic program "Health - automatic" (second week).
	5 DeVita AP	Automatic programs	Without cold - automatic	Alternate every second day during the first week. Use at night.		
			Healthy breathing - automatic			
			7 DeVita AP	Automatic programs	Restoration 2 - automatic	Use after the automatic program "Healthy Day – automatic" on the DeVita Ritm device (first week).
			8 DeVita AP	Automatic programs	Mode 2 - automatic	Alternate every second day during the second week. Use at night.
			9 DeVita AP	Automatic programs	Healthy breathing - automatic	
	10 DeVita AP	Manual program	Without cough, Without pneumonia	If necessary, use after the automatic program on the DeVita Ritm device "Health - automatic" (second week).		
	11 DeLixir		Multi Energy or Top Form	15 ml, in the morning, after breakfast.		
12 DeLixir				pH Balance	15 ml, in the evening, after dinner.	

SCHEME № 3 Maintaining health after infection (2 weeks from start):

SCHEMA №3	DEVICE	PROGRAM		APPLICATION
	1 DeVita Energy	Basic program		1-2 times, in the morning.
	2 DeVita Ritm	Automatic programs	Healthy Day - automatic	Alternate every second day for 2 weeks. Use in the morning.
			Health - automatic	
			4 DeVita Ritm	Automatic programs
	5 DeVita AP	Special programs	Healthy heart	Use after the program "Anti cold - automatic" on the DeVita Ritm device.
			6 DeVita AP	Special programs
	7 DeVita AP	Automatic programs	Restoration 2 - automatic	Alternate every second day for 2 weeks. Use at night.
	8 DeVita AP	Manual program	Mode 4 - automatic	
	9 DeLixir		Multi Energy or Top Form	15 ml, in the morning, after breakfast.
10 DeLixir				pH Balance

## DISEASE PREVENTION RECOMMENDATIONS



Avoid visiting public places: shopping centers, sports and entertainment events, public transport during rush hour.



Avoid close contact and staying in the same room with people who have visible signs of SARS (coughing, sneezing, nasal discharge).



Wash your hands with soap and water thoroughly after returning home, or after contact with strangers. Do not touch your eyes, nose, mouth with dirty hands.



Sterilize gadgets, office equipment and the surfaces you touch.



Limit tight hugs and handshakes when greeting.



Use only personal hygiene items (towel, toothbrush).



Regularly carry out wet cleaning with disinfectants and frequent ventilation of the living spaces.



You will get the most effective result by combining preventive measures with the use of DeVita-technologies and DeLixir cell nutritionstrition

Before  
use



Dry cough, fever up to 39°C and weakness throughout the day.

After the first use  
of DeVita



Improvement of well-being, temperature drop and burst of energy.

After a month of  
using DeVita



Increase of adaptive reserves of the organism, disease reduction, vigor.




Use the recovery scheme of the "Stop Coronavirus- Prevention" Complex with the smart cell nutrition supplements DeLixir MULTI-ENERGY, DeLixir pH BALANCE and DeLixir TOP FORM for much faster results.



# SPECIAL COMPLEXES

(ARE USED ONLY AFTER ONE OF THE BASIC COMPLEXES)

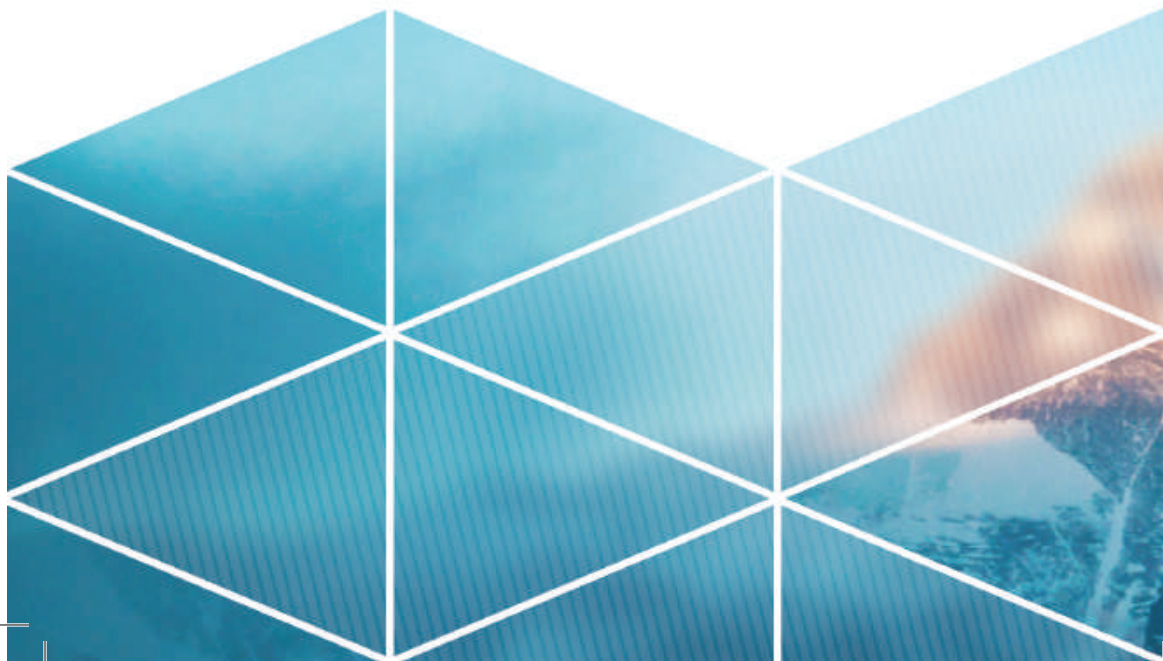





## COMPLEX BODY FIGURE

### WHY IS EXCESS WEIGHT DANGEROUS?



According to WHO, more than 3 billion people on Earth are already overweight, and this number is constantly increasing. Children also begin to suffer from problems with body weight - about 40 million children under the age of 15 weigh much more than the norm. Doctors sound the alarm - being overweight leads to many health problems.



With an increase in body weight, the load on the musculoskeletal system increases. Adipose tissue forms between the internal organs, disrupting their work. Vascular problems and shortness of breath appear. Due to the increased load on the legs, varicose veins develop. All this reduces the immune defense, and overweight people are more likely to get sick.



The main factor in gaining excess weight is unhealthy diet and lack of activity. This is due to the accelerating rhythm of life. We travel by car and public transport, walk a little, have a bite on the run. After work, we feel tired, and we have absolutely no strength and desire to play sports.



For most overweight people, losing weight is a very difficult task. Weight loss is a complex process that includes changing the diet, including physical activity in the daily regimen and the use of special supplements or drinks that help balance the metabolism.

## MECHANISM OF ACTION




**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Parasitic forms that live in the human body cause malfunctions in the work of organs and, as a result, metabolic disorders. And this is one of the main factors in weight gain. DeVita AP antiparasitic programs affect most parasites that negatively affect metabolism, destroying them and effectively removing parasitic load.

“Weight Reduction” automatic program reduces the parasitic load and destroys viruses that cause excess weight. In addition, it normalizes the functioning of the autonomic system, maintains the tone of the body and improves the psychoemotional state.

“Reduce Points” is used to cleanse the body of toxins, normalize lymphatic drainage, maintain immunity and increase energy.

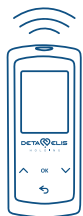

**DeVita Ritm** gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

“High Tonus” automatic program is necessary to relieve toxic load from the liver, lymphatic system and kidneys, activate their work, and improve blood circulation. This helps to increase energy and improve the tone of the body. In addition, the automatic program launches natural rejuvenation processes.

“Healthy Metabolism” automatic program harmonizes the work of the body, regulates the activity of the autonomic system, stabilizes the appetite, and activates the digestive and excretory organs. This normalizes metabolic processes and the psychoemotional state, helps to reduce body weight and improve the quality of life.



**DeVita Energy** activates self-regulation mechanisms, helps accelerate the recovery process and the distribution of energy balance. Allows the body to maintain shape in the process of cleansing and losing weight, maintain vitality and good mood.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the effects of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water for at least 1-1.5 l / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- to reduce your weight;
- to improve the figure for the summer period or an important event;
- for athletes: to “dry” the body before the competition.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir ULTRA SLIM+ smart cell nutrition supplement.

## APPLICATION SCHEMES

SCHEME	DEVICE		PROGRAM		USAGE
	1	DeVita Ritm	Auto-program	High Tonus- automatic	1 <sup>st</sup> day starting at 12:00
2	Auto-program		Healthy Metabolism- automatic	2 <sup>nd</sup> day starting at 12:00	
3	DeVita AP	Auto-program	Weight Reduction- automatic	1 <sup>st</sup> day after the usage of DeVita Ritm or at night	
4		Auto-program	Reduce points- automatic	2 <sup>nd</sup> day after the usage of DeVita Ritm or at night	



It is advisable to use the automatic programs of DeVita AP and DeVita Ritm devices after the general complex of body cleansing, and later, in addition to the DeVita Energy device and the DeVita Ritm device "Healthy Day - automatic" and "Health-automatic" automatic programs.

Attention: when using this program scheme, increased sweating and urination are possible.

## PREVENTION MEASURES FOR INCREASED BODY WEIGHT



Eat 5 times a day small portions.



The last meal should be 2 hours before bedtime.



Drink plenty of water. Water reduces appetite and helps to involve deferred fats in the metabolism, as well as remove toxins.



Eat fruits, they allow you to effectively lose weight.



Vitamins also accelerate the metabolism in cells.



Anti-cellulite massage of problematic areas accelerates blood circulation, the influx to the site of exposure increases, which significantly speeds up the metabolism.



Exercise and physical activity allow you to effectively lose weight and adjust your figure.



Observe your sleep patterns. This contributes to the renewal of brain cells, which, in turn, leads to an acceleration of the metabolism.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements

Before  
use



Increased weight, bad mood,  
heaviness and stiffness in  
movements.

After 1st week of  
using DeVita



Significant improvement  
in immunity, increased overall  
body tone.

After a month of  
using DeVita



Well-being, good mood,  
lightness and ease  
of movement.



Use the scheme of the "BODY FIGURE" complex together with DeLixir ULTRA SLIM+ smart cell nutrition supplement, and your results will be noticeable much faster.



## COMPLEX HEALTHY CELL

### WHY CARCINOGENS ARE DANGEROUS?



Carcinogens are certain factors under the influence of which a person is more likely to develop tumor processes. Carcinogens include organic (food), inorganic (household chemicals) substances, and ionizing radiation.

The development of the process is influenced by bad habits and the ecology of the area.



Smoking is the most significant cause of cancer. Nutrition is another reason for the development of pathological processes, especially in the digestive tract. Civilization costs lead to the fact that food is contaminated with harmful environmental substances and chemicals of agriculture, including carcinogens.



Moderate doses of solar radiation has a beneficial effect on the human body.

But the desire to tan better leads to overheating, skin burns and thermal shock. Such abuse of sunbathing can lead to serious disorders in the body.



Infectious agents (some viruses, bacteria, fungi) with prolonged exposure to the body and certain conditions are also risk factors for the development of such diseases. This was proven back in 1910 and confirmed in the 1980s.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Prolonged presence of carcinogenic factors in the body: biological, chemical, physical, can lead to severe intoxication and uncontrolled division and proliferation of atypical cells in any human organs. The complex antiparasitic programs on the DeVita AP device help in toxic load and inflammatory reactions reduction, neutralize infectious pathogens.



“Cleansing” automatic program is used to remove toxins and reduce the parasitic load on the excretory system after applying the complex of antiparasitic programs.

**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.



“Healthy Day” and “Health” automatic programs have a beneficial effect on the work of the whole organism, increasing its resistance to foreign agents. Strengthens defenses, increases vitality. The use of special programs helps to reduce the load on specific organs and systems, restores their energy potential. The automatic program “Complete Cleansing” helps in main excretory system restoration and harmonizes the body.

**DeVita Energy** harmonizes the body’s work, energizes, activates the removal of toxins, activates the processes of self-regulation.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- if you have a bad diet;
- if you have bad habits;
- if you get sick often;
- if you live in an unfavourable ecological region;
- if you are under constant stress.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir pH BALANCE and DeLixir DETOX smart cell nutrition supplements.

## APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning
	2	DeVita Ritm	Auto-programs	Healthy Day automatic, Health-automatic	Daily in the morning, alternating with each other every other day
	3		Special programs	Digestion normalization and / or Intestine without problems	Daily in the afternoon, if needed
	4	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	5			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	6			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
	7			Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning
	2	DeVita Ritm	Auto-programs	Healthy Day automatic, Health-automatic	Daily in the morning, alternating with each other every other day
	3		Special programs	Digestion normalization and / or Intestine without problems	Daily in the afternoon, if needed
	4	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night
	5		Special programs	No schistosomes (blood-flukes) & trematods	Daily in the evening after DeVita Ritm programs for 1 <sup>st</sup> -2 <sup>nd</sup> weeks
6	Drainage Antitox			Daily in the evening after DeVita Ritm programs during 3 <sup>rd</sup> -4 <sup>th</sup> weeks	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning
	2	DeVita Ritm	Auto-programs	Healthy Day automatic, Health-automatic	Alternate every other day. Turn on at 8:00a.m.
	3		Special programs	Organ Correction Programs	After the DeVita Ritm program "Healthy Day"
	4	DeVita AP	Auto-programs	Cell Protection 1 - automatic Cell Protection 2 - automatic	Alternate every other day within 2 weeks. Turn on at 5:00p.m.
5	Cell Protection 3 - automatic Cell Protection 4 - automatic			Alternate every other day within 2 weeks. Turn on at 5:00p.m.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Ritm	Auto-programs	Healthy Day - automatic	1 <sup>st</sup> day at 8:00 a.m.
	2		Special programs	Organ Correction Programs	After the DeVita Ritm program "Healthy Day"
	3	DeVita AP	Auto-programs	Cleansing - automatic	1 <sup>st</sup> day at night
4	Special programs		Without Papilloma - General	2 <sup>nd</sup> day in the evening and at night	

Programs of 3<sup>rd</sup> and 4<sup>th</sup> months to repeat for another 2 cycles!

## PREVENTIVE MEASURES FOR NEOPLASMS DEVELOPMENT



Exclude sweets, white bread, convenience foods, margarine, smoked meats from ration.



Give up bad habits (alcohol, tobacco).



Timely treat inflammatory infections.



When staying in the open sun, use sunscreen, a hat, clothes that cover open areas of the body as much as possible.



Move more. This contributes to the physiological course of all processes in the body.



Strengthen immunity with regular walks in the fresh air and normalization of the psycho-emotional state.



Have regular physical examinations.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Weakness, fatigue,  
decreased immunity.

After a week of  
using DeVita



Increased  
vitality.

After a month of  
using DeVita



Strengthened immunity,  
improved performance.



Use the scheme of the "HEALTHY CELL" complex together with DeLixir pH BALANCE and DeLixir DETOX smart cell nutrition supplements, and your results will become noticeable much faster.





## COMPLEX LYME DISEASE

### WHAT DANGERS DO TICK BITES CARRY?



Ticks are one of the most diverse and ancient groups of arthropods on Earth. As a rule, mites are free-living, but parasitic blood-sucking species are also found. Among them, the most famous are the ixodid ticks (Ixodoidea). They carry the causative agents of human diseases: tick-borne encephalitis, tick-borne borreliosis (Lyme disease), etc.



Ticks are active in the period: late March, early April, increasing numbers by mid-May. The maximum number is preserved until the end of June, and sometimes until September. Ticks prefer moderately shaded and moistened deciduous and mixed forests. Ticks lie in wait for their prey, sitting on the ends of blades of grass, sticks and twigs sticking up.



Tick saliva is able to dissolve the walls of blood vessels and reduce immunity. Infectious agents are also transmitted through saliva, which pose a danger to human health. Agent transfer occurs in the first minutes after being bitten. Tick-borne encephalitis is the most dangerous of tick-borne infections. Tick-borne borreliosis is the most common.



Infectious agents caught by a tick bite infect the nervous system, musculoskeletal system, and the cardiovascular system. They lead to intoxication and further damage to the liver and kidneys, a general decrease in immunity and a state of constant fatigue.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Parasites that enter the body cause severe intoxication with damage to the nervous, cardiovascular systems and musculoskeletal system. The DeVita AP complex of antiparasitic programs acts on specific pathogens and associated parasites, quickly, effectively and without side effects, removing pathogenic toxic load and signs of inflammation.



“Cleansing” automatic program is used to remove toxins and reduce the parasitic load on the excretory system after applying a complex of antiparasitic programs.

**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

The automatic programs “Healthy Day” and “Health” have a beneficial effect on the work of the whole organism, increasing its resistance to foreign agents. Strengthen defenses, increase vitality. Special programs “Nervous system- general”, “Healthy joints” and “Healthy heart” help reduce the load on these organs and restore their energy potential. “Complete Cleansing” automatic program helps in main excretory system restoration and harmonizes the body.



**DeVita Energy** harmonizes the body’s work, energizes, activates the processes of toxin elimination. Improves the general condition of intoxication resulting from infection after a bite. Reduces lethargy and drowsiness.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- immediately after a tick bite;
- if the tick bite was earlier;
- with problems with the cardiovascular system, joints, kidneys, liver, resulting from a tick bite.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir MULTI-ENERGY, DeLixir pH BALANCE and DeLixir DETOX smart cell nutrition supplements.

## APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Digestion normalization and / or Intestine without problems	Daily in the afternoon, if needed.
	4	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	5			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	6			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
	7			Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Digestion normalization and / or Intestine without problems	Daily in the afternoon, if needed.
	4	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night.
	5		Special programs	No schistosomes (blood-flukes) & trematods	Daily in the evening after DeVita Ritm programs for 1-2 weeks.
6	Drainage Antitox			Daily in the evening after DeVita Ritm programs for 3-4 weeks.	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Health Day automatic, Health - automatic	Alternate every other day. Turn on time 8:00.
	3		Special programs	Nervous System – General	After the DeVita Ritm program “Healthy Day”.
	4	DeVita AP	Auto-programs	Lyme 1 - automatic, Lyme 2 - automatic	Alternate every other day for 2 weeks. Turn on time 17:00.
5	Lyme 3 - automatic, Lyme 4 - automatic			Alternate every other day for 2 weeks. Turn on time 17:00.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Ritm	Auto-programs	Healthy Day - automatic	1 <sup>st</sup> day from 8:00 a.m.
	2		Complete Cleansing - automatic	2 <sup>nd</sup> day from 8:00 a.m.	
	3		Special programs	Nervous System – General	After the DeVita Ritm program “Healthy Day”.
	4	DeVita AP	Auto-programs	Cleansing-automatic	1 <sup>st</sup> day at night.
5	Nerve Regeneration - automatic			2 <sup>nd</sup> day at night.	

Programs of 3rd and 4th months are to repeat for 2 more cycles!

## PREVENTIVE MEASURES FOR TICK BITES



When in tick habitats, avoid dark tones in clothing, as ticks are harder to spot on a dark background.



Pull outerwear into your pants and pants into your socks. If there is no hood, put on a hat.



Every 15 minutes, inspect your clothes, and if possible, inspect the head and body, especially above the waist, ticks are most often found there.



Use chemical protection against ticks. Treat exposed body parts, clothing, sleeping bag, tent and other items.



If you most often spend time in a summer cottage, where ticks live, then this territory should be treated with a special insecticidal agent that kills ticks.



Arriving home, remove clothing outside the living room and carefully inspect it, paying particular attention to folds, seams, pockets.



If you find a sucking tick, do not remove it yourself. Contact the regional clinic or district emergency room.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Headache, weakness,  
lack of appetite.

After a week of  
using DeVita



Reduced pain, increased  
overall tone.

After a month of  
using DeVita



Normalization of the functional  
and psychological state.



Use the scheme of the "LYME DISEASE" complex together with DeLixir pH BALANCE, DeLixir DETOX and DeLixir MULTI-ENERGY smart cell nutrition supplements, and your results will be noticeable much faster.



## COMPLEX HEALTHY TEETH

### WHAT ARE THE DANGERS OF ORAL CAVITY PROBLEMS?



The human mouth is the main gateway to the body. Through it, not only vital substances enter the body, but also pathogenic forms. Therefore, it is necessary to respond in a timely manner to various changes and strange symptoms, because they can indicate not only problems with teeth, but also other organs.



The most common problems are bad breath, tooth decay, sores on the mucous membrane, plaque on the tongue. These manifestations may indicate the presence of viral (herpes), bacterial (streptococcus) and fungal (candida) infections.



Violation of the blood microcirculation, the metabolism can lead to bleeding and inflammation of the gums, which may cause loosening and loss of teeth. The same reasons often cause inflammation of the lips and salivary glands.



The condition of the oral cavity (teeth, mucous membrane, gums, tongue) is an indicator of the work of many internal organs. It is affected by: prolonged use of medications, immune failure, inflammatory processes, unbalanced nutrition, bad habits, vitamin deficiency and other factors.


## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Failure to comply with hygiene rules, bad habits, the presence of chronic inflammatory processes in the body creates a favorable environment for the development of pathogenic organisms in the oral cavity, and as a result, the occurrence of problems with teeth, gums and other organs of the oral cavity. DeVita AP antiparasitic programs help reduce pathogenic infectious load, while Special programs reduce inflammatory reactions and help recover faster after dental interventions.

“Cleansing” automatic program is used to remove toxins and reduce the parasitic load on the excretory system after applying the complex of antiparasitic programs.

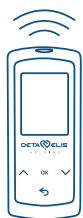

**DeVita Ritm** –In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

The automatic programs “Healthy Day” and “Health” have a beneficial effect on the work of the whole organism, increasing its resistance to foreign agents. Strengthen defenses, increase vitality. The use of special programs improves blood microcirculation, strengthens gums, reduces bleeding, and restores mucous membranes. The automatic program “Complete Cleansing” normalizes the functioning of the main excretory systems and the coordination of the work of the whole organism.

**DeVita Energy** harmonizes the body’s work, energizes, activates the removal of toxins, starts the processes of self-regulation.



## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- if you have a bad diet;
- you have bad habits;
- you are prone to frequent colds;
- you have bleeding gums, loose teeth, damage to the oral mucosa;
- you have planned dental procedures (treatment, removal, implantation, etc.).



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir MULTI-ENERGY, DeLixir pH BALANCE smart cell nutrition supplements.

## APPLICATION SCHEMES

Scheme 1. Use after dental interventions (removal, implantation, treatment, etc.). Duration of use 2 weeks:

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
2	DeVita Ritm	1 <sup>st</sup> day	Auto-program	Healthy Day - automatic	1 <sup>st</sup> day. Turn on at 8:00 a.m.
			Special program	Recovery from injuries	
3		2 <sup>nd</sup> day	Auto-program	Health - automatic	2 <sup>nd</sup> day. Turn on at 8:00 a.m.
			Special program	Strong teeth	
4	DeVita AP	Auto-programs	Restoration 1 automatic		1 <sup>st</sup> day in the evening after DeVita Ritm programs.
Mode 2 automatic			1 <sup>st</sup> day. At night.		
Restoration 2 automatic			2 <sup>nd</sup> day. At night.		

In the presence of abscesses, fistulas add the program "Without fistula" daily.

Scheme 2. Use in chronic problems of the oral cavity:

The first month according to the scheme of 1 month of the complex "Digestion without problems".

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE
	1	DeVita Energy	Mode program	
2	DeVita Ritm	Auto-programs	Healthy Day - automatic, Health - automatic	In the daytime from 8:00. Alternating every other day.
3		Special programs	Strong teeth	After the DeVita Ritm program "Healthy Day".
4	DeVita AP	Auto-programs	Healthy Gums - automatic Cleansing - automatic	At night. Alternating every other day.

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE
	1	DeVita Ritm	Auto-programs	Healthy Day - automatic, Complete Cleansing - automatic
2	Special programs		Strong teeth	After the DeVita Ritm program "Healthy Day".
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed, at night.

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE
	1	DeVita Ritm	Auto-programs	Healthy Day - automatic, Health - automatic
2	Special programs		Strong teeth	After the DeVita Ritm program "Healthy Day".
3	DeVita AP	Auto-programs	No parasites automatic	At night, daily for 1 week.

## PREVENTIVE MEASURES FOR DENTAL PROBLEMS



Brush your teeth twice a day: morning and evening.



Use a damp brush to clean the gums, tongue, and inside of the cheeks..



Timely treat inflammatory infections of the ENT and other internal organs.



Use dental floss to clean the interdental spaces.



Give up bad habits.



Strengthen immunity, normalize the psycho-emotional state, eat right.



Check with your dentist regularly.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Soreness and bleeding  
gums, decreased appetite.

After the first use  
of DeVita



Reduced soreness, improved  
overall well-being.

After a month of  
using DeVita



Reduced bleeding gums,  
strengthened immunity.



Use the scheme of the "HEALTHY TEETH" complex together with the DeLixir pH BALANCE and DeLixir MULTI-ENERGY smart cell nutrition supplements, and your results will be noticeable much faster.



## LIST OF DEVITA AP MANUAL PROGRAMS

Nº p/p	The name of the program	Program description
1	No staphylococcus	<p>According to statistics, 20% of the population are regular carriers of staphylococcus. Staphylococcus is a fairly common causative agent of infectious diseases affecting the skin, soft tissues, upper respiratory tract, and also the musculoskeletal system. The main cause of infection is often the hospital and doctors who become chronic carriers of bacteria and people who have taken antibiotics. Staphylococcus cause secondary diseases with influenza, wound infections, postoperative suppuration, and can also cause food poisoning. The program removes the cause of the infection, which can have a detrimental effect on the human body. It is also effective as a complex treatment for diseases caused by these bacteria.</p> <p><b>Application:</b> daily up to 2 weeks, depending on the stage of the infection.</p>
2	No streptococcus	<p>Streptococcus is a widespread group of microorganisms. They usually live in the skin of a person, in the respiratory tract, intestines. They can cause acute inflammation of the ENT organs, urinary system, organs of the abdominal cavity, skin and soft tissues. Streptococcal infection is more severe than staphylococcal infection, which can lead to the appearance of joint, heart, kidney damage, sinusitis, etc. The program is effective for complex treatment, and separately.</p> <p><b>Application:</b> daily up to 2 weeks.</p>
3	No chlamydia	<p>Chlamydia cause sexually transmitted infections. Chlamydia are often asymptomatic, or with minimal manifestations. The organs of the genitourinary sphere are mainly affected. Women are more susceptible to this infection, children can become infected when passing through the birth canal of the mother. Chlamydia can cause serious complications, therefore it is very important to take timely measures to prevent the disease, observing hygiene rules and using DeVita technologies.</p> <p><b>Application:</b> daily for 2-3 months.</p>
4	No helicobacter	<p>Helicobacter Piloni is a microorganism that lives in the human stomach and duodenum and provokes the development of many dangerous diseases (gastritis, ulcers, hepatitis and many others). According to the WHO, helicobacteriosis infects 60% of the world's population. Infection occurs very easily and is often considered a family disease, because if one member of a family becomes infected, the probability of infection of the rest is 95%. Often the transmission of the parasite occurs as a result of the general use of personal hygiene products, utensils, with a kiss. Symptoms such as heartburn, bloating, pain, and heaviness in the stomach can be caused by this microorganism.</p> <p><b>Application:</b> daily for 2 weeks.</p>
5	No bacteria	<p>Bacterial infections include a huge group of diseases caused by microorganisms - bacteria. The most important in the development of diseases caused by a bacterial infection is the release of toxins by the bacteria, which have an inflammatory effect, cause intoxication of the body and damage internal organs. The "No bacteria" program is used to reduce the severity of intoxication symptoms and alleviate the human condition during bacterial infections of the digestive tract and respiratory system caused by such forms as pneumococcus, listeria, and Klebsiella.</p> <p><b>Application:</b> daily until the disappearance of signs of infection.</p>
6	No hepatitis A	<p>In developing countries and in places where hygiene rules are neglected, the incidence of viral hepatitis A is very high. Since the virus is transmitted through dirty hands, food, water, even through salty (sea) water. If you plan to travel to Asia and Africa, you need to take appropriate safety measures: observe hygiene, drink only bottled water, use the "No hepatitis A" program at the first suspicion of the presence of infection in you or others.</p> <p><b>Application:</b> daily up to 2-3 times a day, at least 1 month.</p>
7	Without toxins	<p>The program promotes the removal of toxins, improves drainage functions, as a result, immunity and vitality are increased. Use is as a detoxification program after parasitic diseases, as well as independently for a month. Unlike other detoxification programs, it does not expend the energy resources of the body.</p> <p><b>Application:</b> once a day after antiparasitic programs. Included in most automatic programs.</p>

8	No herpes zoster	<p>The disease caused by this virus is called herpes zoster. If a person has had chickenpox (one of the forms of the zoster virus), the virus remains in the body in a latent form and hides in the nerve cells. With hypothermia, colds or chronic diseases, the virus comes to the surface of the skin with the formation of characteristic rashes on it and strong pain along any sensory nerve. The program is aimed at the destruction of the virus and the removal of the main signs of its manifestation.</p> <p><b>Application:</b> daily up to 3-4 times a day until the disappearance of undesirable manifestations. At least for 1 month. In conjunction with the "No herpes" program.</p> <p><b>Attention:</b> the program is not used for chickenpox.</p>
9	No herpes simplex	<p>Herpes simplex has been known to mankind since ancient times. Every third person can be a virus carrier. It is localized on the mucous membranes of the lips, mouth, genitals. A decrease in immunity during hypothermia or overheating, during infectious diseases or during the period of menstruation in women contributes to the appearance of unpleasant rashes on the mucous membranes. The virus is transmitted primarily through contact with body fluids of an infected person (blood, saliva, sperm, mucous secretions). The program effectively reduces herpetic burden and prevents the appearance of signs of infection on the mucous membranes.</p> <p><b>Application:</b> at the first sign of a rash every 3 hours, and then 1-2 times a day for 1 month. In conjunction with the program "No herpes».</p>
10	No herpes	<p>Virus carriers are more than 90% of all people on the planet. You can get infected in various ways: through saliva, sex, hygiene items, as well as by airborne droplets. At the beginning, the disease can be asymptomatic, then characteristic signs appear: headache, malaise, fever, herpetic eruptions. The program will help reduce recovery time and reduce inflammatory and pain manifestations.</p> <p><b>Application:</b> 1-2 times a day for a month in conjunction with the "No Herpes Zoster" or "No Herpes Simplex" program.</p>
11	Healthy lymph	<p>Lymph is the liquid tissue of the body. It plays a very important role in maintaining water balance, homeostasis of the microenvironment of cells and capable of leaching pathological microorganisms and toxins from tissues. In case of insufficiency of lymphatic circulation, the lymph does not cope with its main tasks, which leads to the emergence of a pathological process. The program helps to normalize lymphatic outflow, strengthen the walls of lymphatic vessels, improve the condition of lymph nodes, which leads to a decrease in swelling and better discharge of metabolites from the tissues. When using the program, the body's defenses increase.</p> <p><b>Application:</b> if necessary, use 1 time per day or according to the schemes of the complexes (included in the "Cleansing" automatic program).</p>
12	Anti-Cold	<p>The common cold is one of the most common infectious diseases in the world. The causative agents are respiratory viruses (adenoviruses, respiratory syncytial virus) that affect the upper respiratory tract (nose, pharynx, trachea, larynx). Diseases caused by this group of viruses are contagious, because they are transmitted by airborne droplets and are not treated with antibiotics. Complications after a cold - rhinitis, sinusitis, otitis media and more. Very often, an attenuation of the weakened organism can result in the accumulation of a bacterial infection. Therefore, it is advisable to use this program in combination with the antibacterial "No Staphylococcus", "No Streptococcus". The program is aimed at reducing inflammation of the mucous membranes of the nose, throat, bronchi and activation of local immunity.</p> <p><b>Application:</b> at the first sign of a cold the first 3 days - up to 2-3 times a day. Then 1 time per day for 7-14 days.</p>
13	No mould & fungi	<p>Microscopic fungi are ubiquitous. A modern heated apartment or private house is an excellent place for mould and fungi. Favorite mold habitats are bathrooms, libraries, and poorly ventilated living spaces. Fungi can enter the body through food, water, air, contact with old dusty things. Mold fungi cause allergies, not only respiratory, but also cause allergic reactions on the skin. Mold can penetrate the skin barrier wherever we bathe, where we walk barefoot, wherever there are conditions for mold (humidity, temperature, organic residues in the form of dirt or skin fat). Fungi and mold cause diseases of the respiratory organs, gastrointestinal tract, skin with a change in the tissue of these organs.</p> <p>To prevent fungal infections and reduce the occurrence of allergic reactions, it is recommended to use the "No mould &amp; fungi" program. The program is effective against aspergillus.</p> <p><b>Application:</b> 1-3 times a day for 4 weeks. For problems with skin and nails, use the program in conjunction with the «No fungi» program.</p>

14	No candida	<p>The program protects you from infection with such a common fungus as candida. It is part of the natural microflora of the intestine, vagina, and mouth. But with a decrease in immunity, candida begins to multiply actively and harms the body. Diseases are also caused by pathogenic candida, which enter the body through food, through household items, after the use of potent antibiotics, during intercourse with an infected partner. Symptoms are manifested in the form of allergic rashes on the skin, pain and burning in the genitals, dysbiosis. Prevent candida from multiplying using the «No candida» program.</p> <p><b>Application:</b> the program should be applied 1-3 times a day for 2-4 weeks. Use in conjunction with the "No protozoans» program".</p>
15	No fungi	<p>Human fungal diseases belong to the group of infectious diseases caused by pathogenic fungi. The main danger of fungal diseases is that often the areas affected by the fungus are an open road for pathogens of many different infections. Most often, fungi affect the skin, its derivatives (hair, nails), mucous membranes of internal organs. Treatment requires an integrated approach and hygiene. The program actively fights against such fungi as trichophyton, microsporium, and mucor.</p> <p><b>Application:</b> In order to prevent infection with fungal infections, it is recommended to use the program every time after visiting the sauna, pool, beach, communicating with an infected person. If there is a fungal infection in the body, the program is used 1-2 times a day for a month.</p>
16	No protozoans	<p>The program protects you from the negative effects of a large number of protozoa. Every fifth person suffers from them, and becomes infected through dirty hands or dirty food and water. The infected person experiences severe abdominal pain, bloating, constant feeling of fatigue, depression. The liver and pancreas may be affected. Of course, all these signs interfere with your life, and besides, you can infect your loved ones. The program works on the amoeba, lamblia, destroying them. Observing the rules of personal hygiene, as well as the "golden rule" - not to drink tap water and regularly using the program "No protozoans" as a preventive measure, you will be protected from these microorganisms.</p> <p><b>Application:</b> it is recommended to use the program 1-2 times a day for 10-14 days. With re-infection, it is necessary to identify the sources of infection.</p>
17	Intimate health	<p>The causative agents of diseases of the genitourinary system are many and often the cause of the disease is not one microorganism, but several. Infections can occur in acute or chronic form. Chronic infection over time can cause narrowing of the urethra, inflammation of the joints, that can lead to infertility and increase the risk of an ectopic pregnancy. Periodically, during a cold or emotional overload, exacerbations may occur. Infections are mainly transmitted through sexual intercourse, and some during birth from an infected mother to the child. Parasites such as ureaplasma can live in the body for years without causing any symptoms and not manifesting themselves. The main factor in protecting the body is the natural microflora, and modern DeVita technologies that adversely affect the parasitic flora, while maintaining the body's natural defenses. The program effectively fights against parasites such as ureaplasma, toxoplasma, mycoplasma.</p> <p><b>Application:</b> daily for a month.</p>
18	No schistosomes (blood-flukes) & trematods	<p>Trematodes are flatworms that can cause a variety of serious diseases in humans. They can parasitize in the blood (schistosomes), liver, intestines, lungs. Infection occurs when chewing blades of grass, especially near ponds or while bathing in them, when eating poorly washed greens or lettuce, and also not fully cooked cattle meat, fish. The program is aimed at eliminating most types of trematodes, including schistosomes, regardless of the stage of their development and localization in the body.</p> <p><b>Application:</b> daily for 14 days.</p>

19	No nematodes and ascarides	<p>With unwashed fruits, poorly prepared meat, and dirty hands, worm eggs that cannot be seen with the naked eye can get into our mouths. Worms, living in the host's body, feed on its tissues, "take food" from a person (that's why you want to eat a lot), and release toxic substances. A person is losing weight, hair is falling out, skin is peeling. The most common helminths are roundworms. They can cause damage to the wall of the small intestine, intestinal obstruction, inflammation of the peritoneum, bronchitis, focal pneumonia, hemoptysis, liver and pancreas diseases. Ascaris parasitism causes an allergic reaction, which can manifest itself in the form of bronchial asthma and skin rash. Under the action of the program, worms cease to feed, which leads to their death and the natural removal of parasites from the digestive tract. The program does not have a toxic effect on the human body. It is recommended that people with pets, small children and occupations related to working with animals regularly conduct preventive measures once every three months.</p> <p><b>Application:</b> the program is used daily. A total of at least 10 sessions. After 14 days, repeat the course.</p>
20	No parasites, General	<p>Parasites are infesting the human body through the skin, and can also enter the body through food, water, and dust. Many people who lead a healthy lifestyle experience health problems due to parasites. Improving the body (proper nutrition, exercise, hardening procedures) without ridding the body of parasites does not produce a pronounced positive effect. Parasites trigger many chronic diseases: cholecystitis, cholelithiasis, pancreatitis, colitis, diabetes mellitus, bronchial asthma, and dermatitis. Problems such as chronic fatigue, irritability and anxiety, hyperactivity in children, anemia, brittle nails and hair, problem skin, headaches, appetite disorders, decreased immunity can be signals of ongoing parasitosis. All parasitic diseases are terrible in that they perfectly mask as other diseases. Prevention of parasitic diseases is, first of all, observance of hygienic norms of social behavior; use of the "No parasites, General" program. These measures are aimed at maintaining your health and the health of others.</p> <p><b>Application:</b> the program is used daily. A total of at least 10 sessions. After 14 days, repeat the course.</p>
21	No parasites, Digestive tract	<p>The presence of parasites in the digestive tract and liver causes the development of diseases such as cholangitis, biliary dyskinesia, cholecystitis, hepatitis. Their presence in the intestine causes inflammation of the intestinal walls and a decrease in the absorption of nutrients. Toxic substances that secrete parasites cause pain in the joints and muscles, affect the central nervous system. Contribute to the onset of depressive states, nervousness, anxiety, and sleep disturbances. The program is aimed at cleansing the body from parasites of the gastrointestinal tract, liver and pancreas.</p> <p><b>Application:</b> the program is used daily for 7-10 days. After 14 days, repeat the course.</p>
22	Anti-Flu	<p>Influenza is a severe viral infection that affects people of any age. Up to 500 million people get the flu every year. In severe influenza, irreversible damage to the cardiovascular system, respiratory tract, and central nervous system occurs. The virus is transmitted very easily by airborne droplets (by sneezing, coughing, talking). The most common complication of influenza is pneumonia. The program is aimed at the destruction of various strains of influenza viruses, the removal of inflammatory manifestations.</p> <p><b>Application:</b> to prevent influenza during an epidemic, the program is applied once every 2 days. If necessary - daily up to 2-3 times a day, alternating it with the "Against flu, General" program. After that it is necessary to use the "No streptococcus" program. It can be used together with the programs "No staphylococcus", "No streptococcus", "No candida" in order to avoid bacterial and fungal burden, as well as the programs "Antivirus", "Anti-cold".</p>
23	Against flu - General	<p>Influenza is a disease that affects the upper and lower respiratory tract, is accompanied by severe intoxication and can lead to serious complications. Influenza always begins acutely and quickly (48 hours pass from the moment of infection to the first signs). The main symptom of intoxication is increased body temperature (40 °C and higher). This temperature can persist for more than 5 days, which causes: shortness of breath or respiratory failure, impaired consciousness, convulsions. In people with poor health, elderly, young children, this temperature leads to an infectious toxic shock (a condition when vital organs fail). Using the program "Against flu, General" you can alleviate the patient's condition, and with the integrated use of several programs against influenza viruses, you can quickly restore your health without complications.</p> <p><b>Application:</b> at the first sign of flu with increased body temperature, use the program 1-3 times in a row or alternating it with the "Anti-flu" program. Be sure to drink plenty of fluids and observe the rules of personal hygiene.</p>

24	Healthy gastrointestinal	<p>Intestinal infection - a disease that occurs when certain pathogens of viruses or bacteria enter the body. People of all ages are affected by this infection, most often children and those who take antibiotics. An intestinal infection is accompanied: nausea, vomiting, diarrhea, abdominal pain, fever. In case of viral infections, other symptoms can also be observed: runny nose and redness of the throat, as well as pain that occurs when swallowing, lack of appetite, loss of strength. The program is aimed at the destruction of the intestinal flu virus - rotavirus and pathogenic strains of E. coli, helps in the normalization of the condition and stops dehydration that is destructive for the body.</p> <p><b>Application:</b> Use on the first day 4-5 times in a row to achieve a stable result. Then within 3-5 days, 1-2 times daily, alternate with the program "Drainage, General" and "Antivirus". Be sure to drink as much water as possible.</p>
25	Help against poisoning	<p>When enterobacteria or enteroviruses are ingested to the body through food, water, through the mucous membranes of the respiratory tract, an infection develops, which is commonly known as poisoning. This infection can have a latent character or obvious signs of its presence: headache, muscle pain, in the abdominal region, fever, vomiting, impaired bowel movements. Subsequently, it can lead to serious organ damage. The program helps the body cope with the infection caused by enteroviruses and some types of salmonella, accelerates the recovery process, and reduces the likelihood of complications.</p> <p><b>Application:</b> Use on the first day 4-5 times in a row to achieve a sustainable result. Then, for 3-5 days use 1-2 times daily, alternate with the program "Drainage, General" and "Antivirus". Be sure to drink as much water as possible.</p>
26	Antivirus	<p>Viral infections are a large group of diseases caused by viruses. Viruses are a transitional form between animate and inanimate nature. Since viruses have genetic material, but do not have their own cell, they are only intracellular parasites. This is precisely the difficulty of their destruction. Viruses can be transmitted by airborne droplets, through various objects, kisses and sexually. A large role in the complete extermination of viruses in the body belongs to the immune system. The "Antivirus" program is aimed at combating various viruses, and in particular the cytomegalovirus and Epstein-Barr virus. It is also used to attenuate any strains of viruses and enhance the body's immunity.</p> <p><b>Application:</b> daily for a month.</p>
27	Drainage, General	<p>One of the main programs to cleanse the body of toxins. The lymphatic system is involved in the removal of toxic substances from tissues and organs. With various diseases, the amount of harmful substances increases. The lymphatic system does not have enough time to remove them and intoxication of the body may occur again. In addition, due to illness or injury, a violation of the natural outflow of contents from the wound cavity or organ is possible. Excessive accumulation of fluid containing pathogens can cause serious inflammation and deterioration of the human body. By using the "Drainage, General" program you protect your body from the excessive toxic load that it experiences with various diseases or injuries. The program is aimed at neutralizing and removing toxins from the body, while not loading the lymphatic system.</p> <p><b>Application:</b> once after using any of the DeVita AP device programs. Included in the «Cleansing» automatic program.</p>
28	Drainage Antitox	<p>The intercellular matrix is the main place of accumulation of toxins, which leads to metabolic disorders, the development of many serious chronic diseases, and premature aging of the body. Tissue drainage, is a way of detoxification of the tissue sector, helps to remove salts of heavy metals and radionuclides, improves metabolic processes and the general condition of the body.</p> <p><b>Application:</b> with severe intoxication daily and in conjunction with the "Drainage, General" program.</p>
29	Mexican detoxification	<p>The system of deep complete cleansing of the body, which was proposed by R. Rife. The program helps to eliminate metabolic products and toxins, improve the drainage function of the lymphatic system, increase the detoxification capabilities of the liver, activate the excretory function of the kidneys, and increase the level of immune defense.</p> <p><b>Application:</b> daily for a month with serious diseases of the liver, kidneys, lymphatic system. If there aren't any, it is recommended to use the automatic program "Complete Cleansing" on the DeVita Ritm device for gentle cleansing of the body</p>

30	No allergies	<p>Allergy is an increased sensitivity of the body to any substances - products, wool, dust, pollen, chemical, infectious agents. Currently, more than 85% of the entire population of our planet suffer from its various varieties. The cause of allergies is most often malnutrition and an unhealthy lifestyle. An allergy can be recognized by an unexpected runny nose, sneezing or lacrimation, redness and itching of the skin. The program removes the main manifestations of allergies, eases the general condition, and helps stabilize the immune system.</p> <p><b>Application:</b> daily if indicated and to alleviate the condition.</p>
31	Weight Control	<p>The program is aimed at combating adenovirus 36, a common cause of being overweight. In addition, when it is used, the ANS is regulated, and the processes of weight loss are launched. The program also helps in maintaining the body's energy resources.</p> <p><b>Application:</b> it is advisable to use according to the scheme of the complex "BODY FIGURE".</p>
32	Without Salmonella	<p>Salmonella is a representative of the enterobacteriaceae family. Resistant to the external environment: can be stored in house dust, in open water, in frozen meat and sausages, poultry eggs, dairy products from 3 to 6 months. They perish at a temperature of more than 70 degrees. Infection is accompanied by: high fever, abdominal pain, vomiting, diarrhea. However, the infection may be without symptoms.</p> <p><b>Application:</b> daily for 2 weeks.</p>
33	Without Yersinia	<p>One of the representatives of enterobacteria, which causes inflammation in the small intestine, but can also affect other organs: liver, kidneys, heart, joints. It can be found in water, soil, and food. It tolerates low temperatures, so it can be stored and multiplied in the refrigerator. In the human body, without causing signs of inflammation, it can't be located until the immune system is weakened. Then the exacerbation begins.</p> <p><b>Application:</b> daily for 2 weeks.</p>
34	Without Proteus	<p>Bacteria of the genus Proteus are part of the conditionally pathogenic intestinal flora. With active growth, they cause damage to the digestive tract, urinary tract, congenital malformations, and wound infections. The pathogenicity of this type of bacteria is observed in the case of uncontrolled intake of antibiotics. Also, when bacteria enter from the outside through food, water, unwashed hands.</p> <p><b>Application:</b> daily for 2 weeks.</p>
35	Without Klebsiella	<p>Conditionally pathogenic bacterium. With an increase in number occurs, it affects the organs of the gastrointestinal tract, respiratory tract, and urinary system. Causes severe intoxication with impaired functionality of these organs. Ways of infection are through food, water, when communicating with an infected person.</p> <p><b>Application:</b> daily for 2 weeks.</p>
36	Without Campylobacter	<p>Campylobacter bacteria are causative agents of intestinal infection. They enter into the mucosa and lymphatic tract of the intestine and other organs. Infection occurs through agricultural animals and birds, through the use of contaminated meat products, milk, water. Less often through damaged skin, with bites of infected animals. Also transplacental or during childbirth. At risk are rural residents, livestock and poultry farmers, as well as tourists visiting developing countries.</p> <p><b>Application:</b> daily for 2 weeks.</p>
37	Without Clostridia	<p>Clostridium difficile is part of the natural intestinal microflora, but due to their multiplication can cause damage to the walls of the colon. Clostridium botulinum is an inhabitant of soil and sea water, distributed everywhere and forms the most powerful bacterial toxin known - botulinum toxin.</p> <p><b>Application:</b> daily for 2 weeks.</p>
38	Without Shigella	<p>Bacteria that cause dysentery. Infection can occur by contact with an infected person or item, through water, food. Carriers can be flies and cockroaches.</p> <p><b>Application:</b> daily for 2 weeks.</p>
39	Without Rotavirus	<p>The intestinal flu virus multiplies in the cells of the gastrointestinal mucosa, which leads to severe diarrhea and dehydration of the body. The main route of transmission is food.</p> <p><b>Application:</b> daily for 2 weeks.</p>

40	Anti-virus Gastrointestinal	<p>Enterovirus and Koksaki virus belong to the intestinal group, are able to affect the central nervous system, muscles, heart, skin, respiratory and genitourinary systems. They are transmitted from person to person by airborne droplets and through contact with objects of everyday life. The wave of infection is higher during the warm season.</p> <p><b>Application:</b> daily for 2 weeks.</p>
41	Without Trematodes - General 1	<p>Trematodes, or flukes, are parasitic flatworms that infect the liver, pancreas, intestines, lungs, blood vessels, and brain.</p>
42	Without Trematodes - General 2	<p>Ways of infection: swimming in open waters of countries such as Asia, Africa, South America, eating poorly cooked or raw fish, shellfish, crustaceans, and meat.</p> <p><b>Application:</b> the programs are used together, daily for 2 weeks.</p>
43	Stop Amoeba	<p>Entamoeba histolytica is a representative of protozoa. It can cause amoebic dysentery, accompanied by damage to the colon, liver, lungs and other organs. It is ingested with food and water.</p> <p><b>Application:</b> daily for 2 weeks.</p>
44	Without Enterococcus	<p>Enterococcus is a part of the natural human microflora. When using antibiotics, when the barrier function of the skin and mucous membranes are impaired, when older people experience serious diseases, enterococcus can cause problems in the small intestine, urinary tract, vagina, heart, liver, skin.</p> <p><b>Application:</b> daily for 2 weeks.</p>
45	Without nematodes and ascariidae - New	<p>A new program to combat round parasitic worms. Affects the majority of nematode species that infect humans and animals.</p> <p><b>Application:</b> daily for 2 weeks.</p>
46	Stop Adenovirus - General	<p>Adenoviruses are a group of viruses that cause acute diseases, with damage to the respiratory system, eyes, and lymph nodes. Infection occurs by airborne droplets, often in the cold season.</p> <p><b>Application:</b> daily for 2 weeks.</p>
47	Stop Coronavirus	<p>The virus is spread in Europe, Asia, North America. It is transmitted by airborne droplets with damage to the respiratory tract, up to the destruction of the pulmonary alveoli.</p> <p><b>Application:</b> daily for 2 weeks.</p>
48	Stop Penicillium	<p>Genus of molds affecting the organs of the urinary and respiratory systems, lymph nodes, skin. They cause severe intoxication, and also are allergens.</p> <p><b>Application:</b> daily for 2 weeks.</p>
49	Without Rickettsia	<p>Intracellular parasite. It is able to affect any human organs, cause significant disorders of the central nervous system and circulatory disorders. Routes of transmission: through a tick bite, food, contact. Distributed in developing countries.</p> <p><b>Application:</b> daily for 2 weeks.</p>
50	Without Papilloma - General	<p>HPV affects the skin and mucous membranes with the formation of lesions. Some can cause malignant tumors, in the risk zone of women. Ways of transmission of this infection: contact with objects of everyday life, through sexual intercourse, during childbirth from the mother to the child.</p> <p><b>Application:</b> daily for 2 weeks.</p>
51	Without Kidney Papilloma	<p>With the defeat of HPV of the mucous surface of the urinary system, there are lesions left in the renal pelvis and bladder. With prolonged exposure to adverse factors, these lesions can not only reduce the functionality of organs, but also are capable of degeneration.</p> <p><b>Application:</b> daily for 2 weeks.</p>

52	Without Trichomonas	Trichomonas is one of the simplest pathogens of many sexual diseases. Infection occurs with any type of sexual contact, as well as through contact with everyday life objects: the rim of the toilet bowl, the pool, the use of common towels, washcloths, hygiene products, and fitting other people's clothes (underwear). Once in the bloodstream, they can disguise themselves as blood cells. The biggest threat, most often, comes from the infectious agents that it carries: chlamydia, candida, gonococcus, cytomegalovirus and herpes viruses. <b>Application:</b> daily for 2 weeks.
53	Without Leptospira	A bacterium that causes damage to capillaries, liver, kidneys, and muscles. Distributed in all regions. Infection occurs from domestic and wild animals, through damaged mucous membranes and skin, from natural water sources. <b>Application:</b> daily for 2 weeks.
54	Without Cryptococcus	Yeast fungi of the genus Cryptococcus, affect the central nervous system, lungs, skin, mucous membranes. The causative agent is present in soil contaminated with bird feces. The way of transmission is through airborne droplets. <b>Application:</b> daily for 2 weeks.
55	Without Toxoplasma	The prevalence of Toxoplasma infection in the world is very high, up to 90%. Being in the body often does not cause problems for its carrier. Infection occurs either in the uterus or is acquired. Toxoplasma affects the nervous system, eyes, heart, muscles, lungs, intestines, skin. <b>Application:</b> daily for 2 weeks.
56	Without Ureaplasma	Ureaplasma is a pathogenic form that causes diseases of the genitourinary system, joints, pathology of pregnancy and the fetus. The transmission of an infectious agent occurs through sexual contact, in the uterus and during childbirth. <b>Application:</b> daily for 2 weeks.
57	Without Echinococcus	Serious problems are caused by the larvae of this microorganism, which form cysts in the body. These cysts compress the surrounding tissue, which leads to their atrophy and their death. Larvae affect the liver, lungs, kidneys, bones, and brain. Infection occurs through contact with infected animals, when collecting berries and herbs, when drinking water from infected sources. Also, at risk are slaughterhouse workers, shepherds, tanners. <b>Application:</b> daily for 2 weeks.
58	Without Pseudomonas aeruginosa	Pseudomonas aeruginosa is a bacterium characterized by natural resistance to most antimicrobial agents. It affects the organs of the gastrointestinal tract, respiratory system, central nervous system, other organs and systems. Infection sources are people and animals, patients and carriers. Transmission ways: contact with everyday life objects, airborne droplets, food. <b>Application:</b> daily for 2 weeks.
59	Against Lyme Disease 1	Tick-borne borreliosis (Lyme disease) is the most common problem that occurs when a person gets bitten and infected by an ixodid tick. The central nervous system, heart, joints are affected. <b>Application:</b> the program is used together with the program "Against Lyme Disease 2", after bites every day for 2 weeks. For problems arising from bites in the past, see the "Lyme Disease" complex.
60	Against Lyme Disease 2	Tick-borne encephalitis virus is transmitted by bites of ixodid ticks, less often through the use of raw milk from cows and goats infected with the virus. Distribution: Siberia, the Far East, the Urals. The main organ affected is the central nervous system. <b>Application:</b> the program is used together with the program "Against Lyme Disease 1", after bites every day for 2 weeks. For problems arising from bites in the past, see the "Lyme Disease" complex.
61	Healthy Cell - 1	Various pathogenic factors contribute to the accumulation of toxins in the body that can cause cell degeneration. The programs "Healthy Cell - 1" and "Healthy Cell - 2" reduce the risk of oncological pathology, contribute to the correction of pre-tumor changes, reduce inflammatory processes and enhance immune defense. <b>Application:</b> programs are used together, see the "Healthy Cell" complex.
62	Healthy Cell - 2	



63	Healthy Gums	Gum disease (gingivitis) is caused by viruses, bacteria, and fungi. At risk may be people wearing dentures and braces. The program helps to reduce the pathogenic load, reduce the inflammatory process, and improve trophic tissue. <b>Application:</b> daily for 1-2 months.
64	Mucous Membrane Regeneration	The appearance of ulcers (aphthae) on the oral mucosa indicates the presence of an inflammatory process that can be caused by viral or bacterial agents, trauma or a lack of trace elements. The program removes the inflammatory process, helps restore the oral mucosa, and strengthens local immunity. <b>Application:</b> It is advised to use the program in conjunction with the "Healthy Mouth" program.
65	Healthy Mouth	Inflammation of the oral mucosa (stomatitis) occurs due to a number of factors: the presence of infectious agents, damage to the oral cavity, during dental interventions, and an unbalanced diet. It is localized on the tongue, the inner surface of the lips, cheeks, and palate. The program removes the inflammatory process, helps restore the oral mucosa, and strengthens local immunity. <b>Application:</b> daily for 2 weeks. It is advised to use the program in conjunction with the "Mucous Membrane Regeneration" program. Use in any dental procedures accompanied by mucosal damage.
66	Without Fistula	Suppuration of the gums is very often accompanied by the formation of a fistula. This can be fraught with serious consequences from tooth loss to deep tissue inflammation. The program helps prevent suppuration, reduces the inflammatory process. <b>Application:</b> daily for 2-4 weeks. Use after dental treatment or tooth extraction.
67	Stop Injury (teeth)	The program reduces pain, speeds up the healing process, and relieves from swelling. <b>Application:</b> daily for 1-2 weeks. Use immediately after dental procedures for the treatment, removal, implantation of teeth.
68	CNS Without Problems	The program helps to improve the functioning of the nervous system after exposure to various factors (infectious, traumatic, etc.). It helps to restore nerve connections between the central nervous system and peripheral, relieves swelling and pain. <b>Application:</b> for 1-6 months daily. Preferably together with the program "Healthy Nerves" (see the "Lyme Disease" complex).
69	Healthy Nerves	In case of violation of the integrity of the nerve fiber, functional disorders of the nervous system can occur. The program prevents the destruction of the myelin sheath of nerve fibers, thereby preserving the conductivity of the nerve fiber. <b>Application:</b> for 1-6 months daily. Preferably together with the program "CNS Without Problems" (see the complex "Lyme Disease").
70	Clean Sinuses	Inflammation of the mucous membrane of one or more sinuses - sinusitis. It can occur as a complication of infectious diseases of the nose and respiratory tract, as well as after facial trauma. The program helps to reduce the inflammatory process in the sinus area, better discharge of mucus, and reduce pain during acute processes. <b>Application:</b> daily for 2 weeks.
71	Clean Sinuses - General	The greatest problem for many people is chronic sinusitis, which occurs without obvious manifestations, but for a sufficiently long time. The causes may be infections of the upper respiratory tract, the presence of polyps, curvature of the nasal septum, and problems with teeth. The program helps to reduce the inflammatory process in the sinus area, reduce swelling, and prevent exacerbations. <b>Application:</b> daily for a month if there is a problem. And in order to prevent exacerbation during the cold season.
72	Ears Without Inflammation	The causes of ear inflammation can be: a pathological process in the area of the nasopharynx, curvature of the nasal septum, ear injury. The program helps to reduce swelling, reduce inflammation and pain. <b>Application:</b> daily for 2 weeks. Alternate with the DeVita Ritm program "Healthy Ears".

73	Throat Without Pain	<p>Tonsil inflammation of the pharyngeal ring is often observed after suffering tonsillitis or other infectious diseases of the pharynx. May occur in acute or chronic form. Transmission externally - by airborne droplets, or internally - from internal foci of chronic infection. The program helps reduce tonsil inflammation, swelling, pain, and toxic stress.</p> <p><b>Application:</b> daily for 2 weeks. Alternate with the DeVita Ritm program "Throat Without Problems".</p>
74	Against Cold - General	<p>Inflammation of the mucous membranes of the nasopharynx (nasopharyngitis) has an infectious nature, viral or bacterial. The highest frequency is observed in the offseason. The program reduces inflammation, swelling, prevents the development of exacerbations.</p> <p><b>Application:</b> daily for 2 weeks.</p>
75	Anti-inflammation - General	<p>Inflammation is a body reaction that occurs in response to damage or the action of a pathogenic agent. The program helps the body in the elimination of the source of inflammation, wherever it is located, in the reduction of the influence of pathogenic factors, and also in the restoration of homeostasis.</p> <p><b>Application:</b> daily for 1-4 weeks. It is advisable to combine with any other program in order to reduce the inflammatory process.</p>
76	Stop Migraine	<p>Migraine headaches are one of the most common types of headaches. They are characterized by acute throbbing pain localized in a certain part of the head. The program relieves spasms of blood vessels, reduces pain, helps normalize blood circulation in the head.</p> <p><b>Application:</b> as needed.</p>
77	Antiseptic	<p>Damage to the skin: scratches, cuts, burns, etc. are open gates for infectious agents. As a result, the wounds may not heal for a long time, become inflamed, or hurt. The program helps to disinfect the wound, reduce inflammation and quickly tighten it.</p> <p><b>Application:</b> daily for 1-4 weeks. It is advisable to combine with any other program in order to reduce the inflammatory process.</p>
78	Healthy Bile	<p>Improper and irregular nutrition (prolonged fasting, alcohol, spices, overly spicy fatty and overcooked meals, snacks, refusal to consume vegetable oil), stress, a small amount of fluids intake, infectious liver lesions are the causes of problems in the work of the bile. The program facilitates the outflow of bile, helps restore the work of the gallbladder, reduces pain.</p> <p><b>Application:</b> daily for 1-6 months, preferably in the evening. Otherwise, during the day.</p>

## LIST OF DEVITA RITM MANUAL PROGRAMS

№ p/	The name of the program	Program description
1	Active protection	<p>The human body exists in an ever-changing environment, to which it is necessary to adapt. It is adversely affected by changes in temperature and atmospheric pressure, electromagnetic radiation from household appliances and office equipment, and frequent stresses. The body is also affected by physical and mental stress, a constantly changing mode of operation. The program normalizes the functioning of the nervous and endocrine systems, increases working capacity, attentiveness, gives tone, relieves from fatigue and irritability. Protects from pathogenic radiation.</p> <p><b>Application:</b> The program is used daily in the morning or afternoon. It is part of the "Healthy Day" and "Complete Cleansing" automatic programs.</p> <p><b>Attention!</b> It is not recommended for use in the evening and at night.</p>
2	Blood circulation regulation	<p>A sedentary lifestyle, bad habits, poor nutrition, stress negatively affect the state of the circulatory system. The program restores local and general blood circulation, regulates oxygen uptake, improves blood supply to all organs.</p> <p><b>Application:</b> The program can be used at any time of the day if necessary.</p> <p><b>Attention!</b> It is not recommended to use the program for 2 months after myocardial infarction.</p>
3	Digestion	<p>Factors such as a violation of the diet, the transition to unusual foods have an increased load on the gastrointestinal tract. This causes indigestion, flatulence, pain and discomfort. The program has a beneficial effect on the work of the gastrointestinal tract, helps the body absorb food more easily, reduces the risk of eating disorders, reduces the feeling of heaviness and discomfort, and relieves heartburn and intestinal cramps.</p> <p><b>Application:</b> The program can be used at any time of the day, optimally - from 19:00 to 21:00. Reuse if necessary after 0,5 hours.</p>
4	Detoxification	<p>In the process of life, toxins accumulate in the body. They accumulate both in the intercellular substance and in the cells of organs. Toxins disrupt metabolic processes and create a favorable environment for the propagation of pathogens. The program activates the activity of excretory systems and helps to cope with the removal of toxins from the body, and also regulates the liver and gall bladder.</p> <p><b>Application:</b> The program can be used at any time of the day. Be sure to drink clean drinking water – 1-1.5 liters per day.</p>
5	Anti-stress	<p>Almost every day we experience stress: at work, on the road, at home. Usually the psyche neutralizes the effects of stress, but with a constant layering of negativity or with a sharp change in the rhythm of life, the body needs help. The program will help in a difficult situation, with psychological or physical stress. It normalizes the functioning of the nervous system, promotes quick falling asleep and good sleep. For severe nervous agitation, use in conjunction with the Healthy Sleep program.</p> <p><b>Application:</b> The program is used before bedtime. The device should be located under the pillow. If at the end of the program sleep does not occur, you must repeat the session.</p>
6	7 Chakra	<p>Humans are complex energy structures, the basis of which are seven chakras located along the spine. Each chakra is associated with certain organs; if the chakras are damaged, the functioning of their corresponding organs is disrupted. The chakras accumulate energy. Negative energy clogs the chakras and has a detrimental effect on the body as a whole. The program restores lost energy, distributes it correctly and thereby restores the functioning of organs and systems.</p> <p><b>Application:</b> The program is used in the daytime.</p>
7	Hormones normalization	<p>The endocrine system regulates the activity of organs with the help of hormones and ensures the constancy of the internal environment (homeostasis). The endocrine system is involved in the development of the body, its growth, reproductive function. It also regulates the mental state and emotional reactions. The program normalizes the function of the hypothalamus – pituitary – adrenal – gonadal system and regulates excretory processes in the body.</p> <p><b>Application:</b> The program can be used at any time of the day, daily for at least a month. Included in the automatic programs "Healthy Day" and "Complete Cleansing".</p> <p><b>Attention!</b> Not recommended for use by children under 13 years of age.</p>

8	Clean kidneys	<p>The functions performed by the kidneys are diverse. The main one is the regulation of water and electrolyte metabolism in the body, maintaining a constant volume and pressure of blood and other fluids. The kidneys also regulate the acid-base balance in the body. They excrete products of nitrogen metabolism - urea, creatinine, uric acid, etc. Nutritional disorders (especially excessive meat consumption), metabolism, vitamin A deficiency, hypothermia, and infection cause contamination of the kidneys, the formation of sand and stones. The program promotes the active removal of toxins from the body, regulates the functioning of the kidneys, and prevents the formation of kidney stones.</p> <p><b>Application:</b> The program can be used at any time of the day, optimally - from 17:00 to 19:00. In case of toxic or infectious lesions of the body - repeat after 0.5 hours. For the normal functioning of the kidneys, be sure to use clean drinking water - 1-1.5 liters per day.</p> <p><b>Caution:</b> During this program, increased urination may be observed.</p>
9	Clean lymph	<p>The lymphatic system is around the clock guarding the purity of the body. It plays an important role in the removal of metabolic products, as well as in maintaining immunity. Naturally, lymph is cleansed by passing through the lymph nodes, but sometimes this is not enough. In violation of the functions of the lymphatic system, up to 83 harmful substances accumulate in the intercellular space. This increases the load on all organs of excretion: liver, intestines, kidneys. The program helps to liquefy the lymph, remove toxins and excess water from the body, improves metabolism in cells, normalizes blood pressure, and boosts immunity.</p> <p><b>Application:</b> The program can be used at any time. During the operation of the program, it is recommended to carry out a light massage in the form of stroking (without effort).</p> <p><b>Caution:</b> For the most effective removal of toxins, use clean drinking water – 1-1.5 liters per day.</p>
10	Healthy eyesight	<p>A constant load on the eyes (reading, working on the computer), as well as a change in conditions and rhythm of life, cause eye fatigue and irritation, and reduced vision. The program reduces eye strain, relieves from fatigue and pain, and helps normalize vision.</p> <p><b>Application:</b> Use before bedtime daily for 7 days. After 7 days, use it every second day. Repeat sessions until the desired result is obtained. After the end of the program, do not strain your eyes. For best results, combine with exercises for the ciliary muscle.</p>
11	Female sexuality	<p>Sexual desire is an indicator of different conditions of a woman. If a woman ceases to be interested in sexual life, it can be assumed that the disease has arisen. In most cases, a decrease in sexual desire develops with a combination of several factors - somatic, mental and microsocial. The program will help in the absence or decrease in a woman's sexual desire, weaken the feeling of emotional stress, mental discomfort, and strengthen sexual desire.</p> <p><b>Application:</b> daily 1 time per day.</p>
12	Male sexuality	<p>The manifestation of male attraction depends on age, physical exertion, nervous and mental state, external stimuli, etc. Stress, problems at work, various complexes lead to psychological impotence in men. As a result, they become wary of the appearance and "quality" of an erection, compare, analyze. And the more a man prepares for this "event", the less chance he has of success. The program is aimed at psychological correction, regulation of the 6th chakra, relieves anxiety and strengthens confidence in the favorable completion of intimacy.</p> <p><b>Application:</b> daily 1 time per day.</p>
13	Recovery from injuries	<p>Injuries, dislocations, fractures lie in wait for both supporters of various sports and ordinary people. In combination with medical care, the program speeds up the healing process, reduces pain, relieves swelling, eliminates hemorrhages, and restores muscle elasticity.</p> <p><b>Application:</b> for injuries daily 2-3 times a day until the desired result is obtained. If necessary, combine with the "Joints without inflammation" and / or "Joints without pain" programs.</p>
14	Warm-up	<p>Often, gym visitors lack energy and motivation, especially after a working day. The Warm-up program energizes, increases the intensity of training, helps to overcome fatigue, increases metabolic processes in muscles and makes training more effective.</p> <p><b>Application:</b> before each workout.</p>

15	Workout	<p>Successful and active training often leads to fatigue and muscle pain, which do not make it possible to achieve greater results. The program increases muscle strength, relieves muscle fatigue, strengthens muscles, improves the supply of blood and oxygen to muscles, increases endurance, and allows you to achieve maximum results with minimal cost to the body.</p> <p><b>Application:</b> The program is used during training.</p> <p><b>Attention!</b> It is not a doping program.</p>
16	Healthy Heart	<p>The supply of organs with oxygen and nutrients depends on the work of the heart, which means the effective work of the whole organism. Malfunctions in the activity of the heart are caused by an improper lifestyle, stress, and physical overload. The program regulates complex blood circulation and blood pressure, provides normal blood supply to the heart muscle, and promotes self-regulation of the heart and its stable operation.</p> <p><b>Application:</b> Use at any time, once, twice, if necessary (optimally – from 11:00 to 13:00).</p> <p><b>Attention!</b> It is not recommended to use the program within 2 months after myocardial infarction.</p>
17	Cardiac rhythm normalization	<p>When the heart beats as it should, we do not notice its work. But if it begins to stray from its rhythm, it brings great discomfort. Arrhythmias can occur in healthy people for various reasons – as a reaction to tight clothing, insect bites, certain medications, stress, plentiful food, and problems with bowel movement. The program restores the heart rhythm, improves blood circulation to the heart, normalizes the nervous system.</p> <p><b>Application:</b> The program is used for heart rhythm disturbances 1 time per day.</p>
18	Anti-Cold	<p>The first signs of a cold are a stuffy nose and a sore throat. Often we do not pay attention to these symptoms which leads to complications, the most common of which are upper respiratory tract infections - sinusitis, bronchitis, and otitis media. The program activates the immune system, increases resistance, reduces the likelihood of complications, reduces swelling, promotes better discharge of mucus from the sinuses of the nose, and helps to recover in a short time.</p> <p><b>Application:</b> The program can be used at any time of the day, during the off-season or epidemic of colds, use daily. It alleviates the condition with allergies.</p>
19	Free breath	<p>Bronchospasm (cough) is a protective reaction of the body to an irritant, respiratory tract infection, and may also be a consequence of the effects of certain drugs on the body. Bronchospastic syndrome is characteristic of people with a tendency to allergies and those suffering from bronchial asthma. The program restores the functions of the bronchopulmonary system, normalizes lung function, regulates blood supply and oxygen absorption, and relieves from spasms.</p> <p><b>Application:</b> daily as needed at any time of the day. With cramping and severe coughing – up to 2-3 times in a row.</p>
20	Women's Health	<p>Bad ecology, bad habits, the sun adversely affect the female body, accelerate the aging process and cause problems in the female body. The program covers a wide range of women's issues and is useful at any age. Reduces premenstrual tension, pain, corrects menopausal disorders.</p> <p><b>Application:</b> Girls at puberty are recommended to use it 1-2 times a week. For women experiencing pain and discomfort during menstruation – daily during this period. For the effectiveness of the impact in the correction of problems in the sexual sphere and / or in violation of the organs of the female reproductive system – daily in combination with the program "Female sexuality".</p> <p><b>Attention!</b> It is not recommended to use the program during pregnancy or if a woman wants to become pregnant.</p>
21	Men's Health	<p>Man's health is not only a healthy potency, but also the proper functioning of the rest of the internal organs. It is necessary to take care of male health from youth. The program is designed to harmonize the male genitourinary system. It regulates the processes of urination, the work of the male sex glands, including the representative. The older the man, the more useful this program.</p> <p><b>Application:</b> The program is used once a week, optimally from 21:00 to 23:00. For problems with potency and / or prostate every day in combination with the program "Male sexuality", "Stop Prostatitis".</p>

22	Healthy sleep	<p>Healthy sound sleep is an important component of the life of any person. During sleep, the body restores the energy spent during the day. However, any psycho-emotional upheaval can serve as an impetus for the development of sleep disorders. Which, in turn, provokes the development of both physiological and mental disorders. Irritability, fatigue, inability to concentrate are the very first signs of a sleep disorder. The program helps to reduce the excitability of the nervous system, quickly fall asleep, normalize sleep phases.</p> <p><b>Application:</b> daily at night, it is possible to reuse the program up to 2-3 times.</p>
23	Good mood	<p>Depression, anxiety, fear - these concepts are interconnected. Anxiety and fear are a signs of depression. Very often, this condition leads to psychosomatic diseases, emotional exhaustion, self-dissatisfaction. Fears are powerful factors that prevent a person from discovering his abilities and achieving success. The "Good mood" program allows you to stop a depressive state, reduce anxiety, normalize sleep, and improve mood and mental activity.</p> <p><b>Application:</b> daily, 1 time, during the day..</p>
24	Joints without pain	<p>Joint or spinal pain can have several causes. Most often, this is a pinched or inflammation of the nerves, which occurs due to a sedentary lifestyle, bruises, falls, being overweight, spinal disc herniation, hypothermia, active physical exertion, and destruction of cartilage tissue or inflammation of the joint bag. The program helps to reduce pain, improve blood microcirculation, normalize muscle tone, ligaments, cartilage, activates metabolic processes in damaged tissues.</p> <p><b>Application:</b> in case of exacerbation, use the regimen a support to the main therapy and in combination with the «Joints without inflammation» program. Daily, at any time of the day, for 2 weeks, then every second day, then repeat the course until the desired results are obtained.</p>
25	Joints without inflammation	<p>Every third inhabitant of our planet experiences pain or a crunch in the joints and spine. The causes of this problem may be infections, injuries, allergic reactions, being overweight, professional sports, improper nutrition, etc. The program improves blood supply to the damaged area, relieves from inflammation, swelling and pain. It has established itself as a supportive and restorative means during the period of "active successful work", as well as in sprains and severe muscle fatigue.</p> <p><b>Application:</b> daily for a month, starting 2 weeks before active actions in the country or sporting events. Then 1-2 times a week until the end of the season. For people with periodically worsening problems of the back and joints, in order to maintain health 2 times a year daily for 1 month. If an exacerbation occurs, combine with the Joints Without Pain program.</p>
26	Intestine without problems	<p>How many of us have experienced symptoms such as bloating, nausea, bowel problems, decreased appetite, and abdominal pain? All these are indicators of intestinal abnormalities. These symptoms appear not at the initial stage, but when the disease has turned into a chronic form. To reduce pain, remove nausea and bloating, and normalize bowel movement, use the Intestine Without Problems program.</p> <p><b>Application:</b> daily as needed and depending on the condition 1-2 times a day.</p>
27	Digestion normalization	<p>Metabolism is a combination of chemical reactions that ensure the normal functioning of the body. The main factors of metabolic disorders are: malnutrition, the use of certain medications, stress, a sedentary lifestyle. Due to metabolic disorders, diseases such as diabetes, gout, obesity, allergies, cardiovascular disorders, growth retardation and development in children can occur. The program is aimed at the balanced functioning of the thyroid and pancreas, the proper assimilation of potassium, calcium, magnesium, sodium, and the alignment of acid-base balance.</p> <p><b>Application:</b> as needed 1-2 times a day.</p>
28	Beautiful skin	<p>Skin problems arise due to a malfunction of the internal organs - in this case, the skin takes on the function of cleaning. Acne, vascular defects, pigmentation, etc. appear on the face. This is faced by adolescents experiencing hormonal changes, and people suffering from disorders of the endocrine and digestive systems. The program normalizes the work of the skin glands, enhances the barrier function of the skin, improves the metabolic processes in the cells.</p> <p><b>Application:</b> 1-2 times a day until the desired result.</p>

29	Beautiful hair and nails	<p>Hair loss, dandruff, brittle nails - indicators of disorders of the body or the negative impact of the environment. The program helps improve blood microcirculation, enhance cellular metabolism and cell regeneration. This leads to the strengthening of hair follicles and nails, stimulating hair growth, increasing the protective functions of the skin.</p> <p><b>Application:</b> daily, until the desired result is obtained.</p>
30	Without allergens	<p>Allergens are substances to which the body is hypersensitive due to a disturbed immune response mechanism. They surround us everywhere: on the street, at home, at work. Common allergens are plant pollen, dust, animal hair, medicines, household chemicals, etc. The body can respond to these substances with symptoms such as lacrimation, swelling of the mucous membranes, skin rashes, etc. The program alleviates allergy symptoms, normalizes the functions of the immune system.</p> <p><b>Application:</b> daily as needed. People allergic to flowering plants are advised to start using 1 month before flowering.</p>
31	Immune protection	<p>The immune system is a very complex, self-regulating, self-learning and multi-level complex. Thanks to its well-coordinated work, biological protection is ensured. The program reduces the load on the liver, activates the work of the thymus, spleen and lymphatic system, and promotes rejuvenation and longevity.</p> <p><b>Application:</b> daily or every other day, at any time of the day. It is part of the "Health" automatic program.</p>
32	Healthy kidneys	<p>The program is used to remove excess fluid from the body, as well as salts that have the ability to retain it in the tissues. The program has a pronounced diuretic effect, and also helps to restore the physiological function of the kidneys and reduce swelling.</p> <p><b>Application:</b> daily as needed for a month.</p>
33	Vegetative regulation	<p>The autonomic nervous system controls many of the body's life support functions. Problems in its work can be expressed in violation of the heart rhythm, changes in blood pressure, excessive sweating, fatigue, weakness, brittle hair and nails, etc. The program is aimed at restoring the ANS functions, activating the brain, normalizing metabolism, heat transfer, appetite, cell regeneration.</p> <p><b>Application:</b> daily or every other day, at any time of the day. It is part of the "Health" automatic program.</p>
34	Cellular respiration	<p>Hypoxia - oxygen starvation, which may be caused by a lack of oxygen entering the body or the termination of its absorption by body tissues.</p> <p>The program helps to improve the blood supply to the head, heart and lungs, the active absorption of oxygen by the blood, relieve vasospasm, reduce acidification, and eliminate toxins.</p> <p><b>Application:</b> daily or every other day, at any time of the day. Included in the automatic program "Health".</p>
35	Healthy joints	<p>Pain in the joints can occur not only in older people, but also in young people. Pain is associated with inflammatory reactions, accompanied by swelling, redness of the skin, restriction of movement. The causes of this problem can be: infections, metabolic disorders, injuries, being overweight, large mechanical and functional loads. The program helps to improve the condition of damaged joint tissues, and to reduce pain.</p> <p><b>Application:</b> daily for a month. It is recommended to be used in conjunction with the "Joints without pain", "Joints without inflammation" programs. The most effective program for performing light exercises on a problem joint.</p>
36	Stop prostatitis	<p>The prostate gland is an organ that performs a number of necessary functions: reproductive, secretory, and barrier, controls the outflow of urine, and affects the psychoemotional state. The gland is very vulnerable to various types of pathological processes. Factors that disrupt the functioning of the prostate can be both external (poor ecology, ionizing radiation, stress, passive lifestyle, exposure to infectious agents), and internal (genetic predisposition, concomitant chronic and acute diseases, bad habits). The program is aimed at reducing inflammatory processes and normalizing the gland.</p> <p><b>Application:</b> daily for 1-3 months. It is recommended to combine with the "Men's health" program.</p>

37	Balanced nervous system	<p>Nervous overexcitation is a consequence of stress accumulated over a certain period of time. This can occur in the form of: insomnia, depression, muscle or headaches, irritability, mood instability, etc. The program removes nervous excitement, harmonizes the work of the nervous system.</p> <p><b>Application:</b> with great psychological stress every day at the end of the working day and / or before bedtime.</p>
38	Recovery after training	<p>During training, the body activates the metabolism, which leads to the accumulation of metabolites. The load on the autonomic system increases, microtraumas of muscles, their soreness, etc. can be observed. Not everyone can withstand such discomfort, and as a result, they stop further training. The program relieves muscle fatigue, promotes the rapid withdrawal of metabolites and the maintenance of the nervous system.</p> <p><b>Application:</b> after any kind of training.</p>
39	Healthy venous system	<p>We pay attention to the problem with veins when blue protruding "snakes" are formed on the legs (most often), delivering, initially, aesthetic discomfort. And further leading to other functional disorders. Most often, the cause of this problem is a person's lack of mobility, lack of training or regular walking, hours of sitting behind the wheel or at a computer monitor, and long-term wearing high-heeled shoes. The program helps to relieve pain and swelling, strengthen the venous wall and increase its tone, reduce inflammatory reactions, improve blood circulation in the capillaries.</p> <p><b>Application:</b> if necessary (with severity, leg fatigue), or daily (in the presence of varicose veins). It is advisable to use the program at night or before bedtime. In the second case, put a folded blanket or pillow under your feet so that your legs are in a raised state. After the session, do not get up for 30 minutes.</p>
40	Without headache	<p>A headache can cause a lot of problems. The causes of a headache are different: blood pressure drops, intoxication, autonomic dysfunctions, osteochondrosis, etc. The program is aimed at eliminating pain, relieving vascular spasms, improving blood supply and nutrition of brain cells.</p> <p><b>Application:</b> as needed at any time of the day. If at the end of the program the pain persists, repeat the session after 20-30 minutes.</p>
41	Normal blood sugar	<p>Elevated blood sugar can be observed in humans under the influence of a number of factors: hard physical or intense mental work; carbohydrate-rich foods; severe stress, panic attacks; taking medications, functional circulatory disorders. Monitor your blood sugar regularly. To prevent the development of a more serious disease, use the "Normal blood sugar" program, which will help normalize the pancreas.</p> <p><b>Application:</b> daily for a month, then as needed. It can be combined with the programs "Digestion normalization" and "Healthy Pancreas".</p>
42	Functional memory +	<p>The human brain contains a large amount of knowledge. With age, as well as with great mental stress, a decrease in the process of memorization can be observed and memory is not always capable of issuing the necessary information. The program activates brain activity and the process of memorization, and, in combination with a healthy lifestyle, sets up the psycho-emotional component of a person to be positive.</p> <p><b>Application:</b> The program is recommended to be applied daily in the morning and afternoon, preferably in the process of studying the necessary material.</p>
43	Nervous System - General	<p>After tick bites, very often there may be signs of a malfunction of the nervous system, expressed as headache, fatigue, insomnia, depression, impaired movement of the facial muscles, malfunctions in the functioning of the optic nerve, etc. The program maintains the nervous system in a harmonious state, brings all its departments are coordinated to work together, establishes the functional and psychological state of a person, and regulates mental processes and mood.</p> <p><b>Application:</b> within 1-6 months. See the "Lyme Disease" complex.</p>
44	Healthy Ears	<p>Inflammation and pain in the ears can occur as a result of injuries, damage to various infectious agents: viruses, bacteria, fungi. The problem can occur at any age, more often observed in children and the elderly. The program is aimed at restoring the mucous membrane of different parts of the hearing organ, reducing discomfort in ear diseases, and improving hearing.</p> <p><b>Application:</b> daily for 1-2 weeks. Alternation with the DeVita AP program "Ears without inflammation" is desirable.</p>



45	Throat Without Problems	<p>Sore throat is a problem that occurs in both adults and children in different periods of the year. The reasons may be: decreased immunity, infectious diseases of the upper respiratory tract, carious teeth, gum disease. The program helps maintain local immunity, restore the nasopharyngeal mucosa, reduce pain, swelling, and hoarseness.</p> <p><b>Application:</b> daily for 1-2 weeks. Alternation with the DeVita AP program "Throat without pain" is recommended.</p>
46	Strong teeth	<p>Gum disease is a very big problem for a person, causing not only discomfort, but also capable of leading to tooth loss. The reasons are often: decreased immunity, impaired blood circulation, bad habits. The program helps to normalize blood microcirculation, strengthen gums, and repair damaged mucous membranes, thereby improving the condition of the oral cavity and teeth.</p> <p><b>Application:</b> daily for 1 month, see the complex «Healthy teeth».</p>
47	Healthy Pancreas	<p>The pancreas performs two main functions in the body: externally and intracretory. An imbalance in the functioning of the gland can lead to serious problems in terms of the functioning of the digestive system and metabolic disorders. The program helps to normalize the functions of the gland in various disorders.</p> <p><b>Application:</b> daily for 1-6 months. It can be combined with the programs "Digestion normalization" and "Normal blood sugar".</p>
48	Healthy Thyroid	<p>The thyroid gland is part of the human endocrine system, which plays a very large role in the body. If its functioning is disturbed, an imbalance in the work of the cardiovascular, respiratory, urinary, reproductive systems, as well as metabolic and fluid disorders may appear. The program provides support to the thyroid gland under various stresses, normalizing its functioning and metabolic processes in the body.</p> <p><b>Application:</b> daily for 1-6 months. It can be combined with the program "Digestion normalization".</p>



## ELIXIR OF HEALTHY LIFE

# DELIXIR

## HOW FORMULAS WORK



Fill with energy and vitality. A balanced complex of vitamins and minerals, packed in an easily digestible formula, gives vigor and energy, relieves from fatigue and stress.



Maintain the acid-base balance of the internal environment of the organism. Modern nutrition, modern lifestyle and stress lead to acidification of the organism, and this as a result leads to early aging and the development of diseases



Gently, but effectively cleanse the body of toxins. Natural sorbents and plant antioxidants remove accumulated toxins, allowing organs to work more efficiently.



Maintain the internal and external youth intact. A special formula with collagen and hyaluronic acid makes the skin smooth and the muscles and ligaments elastic and mobile.



They support a healthy metabolism, accelerating the conversion of nutrients into energy, thereby contributing to weight loss. Protect the excretory system from toxins and remove excess fluid.

## MULTI-ENERGY



Energizes with vitamins and minerals.



Contains Coenzyme Q10, Selenium and Lutein.



Relieves eye strain and fatigue.



For those who want to lead an active lifestyle.

ORANGE FLAVOUR

## DETOX



Cleanses from toxins and restores the body.



Complexly strengthens and revitalizes.



Contains aloe vera, noni juice, MSM, goji berries and ginger.







For those who live in the rhythm of a modern city.

KIWI FLAVOUR

## COLLAGEN +







-  Supports youthful skin, reduces age-related manifestations.
-  Contains Collagen, Hyaluronic Acid, Acai Berries.
-  Maintains the elasticity of ligaments, joints and tendons.
-  For those who want to keep their internal and external youth intact.

BLACK CURRANT FLAVOUR

## PH BALANCE







-  Harmonizes the acid-base balance of the organism.
-  Contains spirulina, turmeric, marine minerals.
-  It has an antioxidant effect.
-  For those who care about their health.

LEMON FLAVOUR

## ULTRA SLIM+







-  Normalizes energy metabolism and the process of splitting fats.
-  Restores the functioning of the excretory system, helps in excess fluid removal.
-  Reduces the absorption of fats and glucose in the intestines, contributes to weight loss.
-  Improves metabolism, ensures proper processing of food, creating the perfect figure

CRANBERRY FLAVOUR

## TOP FORM



-  Promotes general body enhancement and improves performance efficiency
-  Maintains metabolism, skin and muscles tone, new healthy cells growth
-  Provides quicker recovery after cold and intensive exercises
-  Protects DNA cells, lipids and protein free radicals

TANGERINE FLAVOUR

## EXPLANATIONS FOR THE USE OF PROGRAMS DEVITA RITM

DeVita Ritm (manual programs)				
Name of the program	In which cases the program can be used	Effect	First results	Lasting effect
Active protection	With physical and psychological stress, decreased performance, attentiveness and tone, with chronic fatigue and irritability and daily work with household and computer equipment. To remove the bioenergetic load when working with people.	Immunostabilizing Gives energy Anti-stress	Within an hour	1 month
Blood circulation regulation	Violation of blood microcirculation (blood supply to organs), lack of oxygen in the blood, headaches, impaired blood pressure and/or vascular wall tone.	Relieves from spasms Supports tone Expands blood vessels Normalizes blood pressure	Within an hour	1 month
Digestion	Overeating, heaviness in the stomach, heartburn, bloating.	Relieves from spasms Anti-inflammatory Normalizes the acid balance	Within an hour	1 month
Detoxification	High physical activity, malnutrition, food poisoning, pain in the right hypochondrium, alcohol consumption.	Antitoxic Protects the liver Improves the outflow of bile Relieves from spasms	Within an hour	1 month
Anti-stress	Stress, high physical and psychological stress (study, work, sport)	Anti-stress Calming effect	Within an hour	1 month
7 chakras	Chronic fatigue, headaches, decreased immunity, high physical and psychological stress.	Gives tonus Harmonizing effect Anesthetic	Within an hour	1 month
Hormones normalization	In case of puberty disorders in adolescents from 12 years old, after surgical interventions, with hormonal failures, with reproductive system dysfunction in men and women.	Normalizes the hormonal background	Within 1 month	6 months
Clean kidneys	Inflammatory processes in the kidneys, bladder, swelling.	Improves urine output Decongestant Anti-inflammatory	Within 1 day	1-2 months
Clean lymph	Impaired lymphatic system, swelling, intoxication, impaired cellular metabolism, impaired blood pressure, decreased immunity, skin diseases.	Lymphatic drainage effect Decongestant Antitoxic Immunostabilizing	Within 1 day	1 month
Healthy eyesight	Decreased muscle tone of the eyes, visual impairment, visual fatigue. For prevention of visual impairment during prolonged work with computer equipment or texts on paper.	Correction effect Gives tonus Strengthening	Within 1 day	1-3 months
Female sexuality	Decreased libido, increased emotional excitability, a tendency to nervous excitement.	Antidepressant Activating	2-3 days	1 - 2 months
Male sexuality	Increased stress, excitability of the nervous system, decreased potency, fear of "failure".	Antidepressant Activating	2-3 days	1 - 2 months

Recovery from injuries	Any injuries of bones, cartilage and soft tissues, ligamentous apparatus: fractures, bruises, dislocations, cuts. After surgery.	Anesthetic Decongestant Anti-inflammatory Regenerating	Within an hour	1 month
Warm -up	Used before physical activity.	Gives tonus Gives energy	Within an hour	1 month
Workout	Used during exercise.	Gives tonus Gives energy Regulatory	Within an hour	1 month
Healthy heart	Disorder of the functional activity of the heart, impaired heart rate, impaired blood pressure, high stress loads.	Relieves from spasms Normalizes blood pressure Gives tonus Anesthetic	Within an hour	1-2 months
Cardiac rhythm normalization	Heart rhythm disturbances: tachycardia, arrhythmia, including respiratory and atrial fibrillation.	Relieves from spasms Normalizes the cardiac rhythm	Within an hour	1-2 months
Anti-cold	Diseases of the upper respiratory tract. Allergy, accompanied by swelling of the nasopharynx, lacrimation.	Anti-inflammatory Decongestant Immunostabilizing	Within an hour	Within 1 month
Free breath	Diseases of the lower respiratory system. Allergy. Spasms of the bronchopulmonary system of various origins.	Anti-inflammatory Relieves from spasms Expands blood vessels	Within an hour	1 month
Women's health	Dysfunction and inflammatory processes of the reproductive system.	Anti-inflammatory Expands blood vessels Hormone-normalizing Anti-aging	2-3 days	6-12 months
Men's health	Reproductive system dysfunctions.	Anti-inflammatory Expands blood vessels Hormone-normalizing Anti-aging	2-3 days	6-12 months
Healthy sleep	Insomnia, sleep phase disturbance, increased nervous irritability.	Antidepressant Calming Normalizes sleep	Within 1 day	1 month
Good mood	Depression, irritability, psycho-emotional stress, fatigue.	Antidepressant Regulating mood	Within 1 day	1 month
Joints without pain	Inflammation and pain of the musculoskeletal system (joints, muscles, spine), high physical exertion, including sports.	Anesthetic Regulatory	Within 1 day	1 month
Joints without inflammation	Inflammation of the musculoskeletal system, swelling, trauma, with high physical exertion, including sports.	Anti-inflammatory Decongestant Regulatory	Within 1 day	1 month
Intestine without problems	Pain, nausea, bloating, dysbiosis, decreased tone of the intestines, gall bladder, violation of the gastric mucosa and duodenum, liver function.	Relieves from spasms Anesthetic Improves the outflow of bile Anti-inflammatory Normalizes metabolism Decongestant	Within 1 day	6-12 months
Digestion normalization	Thyroid, pancreas and gastrointestinal dysfunction, swelling, impaired absorption of K, Ca, Mg, Na, impaired acid-base balance.	Regulatory Harmonizing	5-7 days	6-12 months

Beautiful skin	Decreased barrier function of the skin and cellular metabolism. Acne, dryness or excessive oily skin, relieves from psoriasis and vitiligo.	Regulates the immune and barrier functions of the skin Normalizes cell nutrition	5-7 days	1 month
Beautiful hair and nails	Hair loss and fragility, brittle nails. Dry and oily seborrhea.	Regulating the function of the skin glands Normalizes cell nutrition	5-7 days	1 month
Without allergens	Various types of allergies (respiratory - odors, dust, pets, hay fever, skin, food), swelling of the upper respiratory tract.	Antitoxic Antiallergic Decongestant Anti-inflammatory Immunostabilizing	Within 1 day	1-2 months
Immune protection	Decrease in immune defense, increased physical and mental stress.	Immunostabilizing Anti-inflammatory	Within 1 month	1 month
Healthy kidneys	Violation of metabolic processes associated with the withdrawal of certain acids and salts, impaired functional activity of the kidneys (formation of kidney stones, deposition of salts in the joints, high blood pressure).	Improves urine output Decongestant	Within 1 day	1-2 months
Cellular respiration	Headaches, problems with the vessels of the head, lethargy, drowsiness, lack of oxygen in the blood, a sedentary lifestyle. For residents of large cities.	Relieves from spasms Expands blood vessels Normalizing Fills with oxygen	Within 1 day	1-2 months
Vegetative regulation	When doing hard work for prolonged periods of time, psychological stress, imbalance of the ANS.	Regulatory	Within 1 day	1-2 months
Healthy venous system	Decreased tone of veins, stagnation in the veins of the lower extremities.	Anti-inflammatory Anesthetic	Within an hour	1-2 months
Without headache	Spasms of cerebral vessels, intense mental work, instability of blood pressure, impaired blood microcirculation in the head.	Anesthetic Relieves from spasms	Within an hour	1-2 months
Normal blood sugar	Pancreatic dysfunction associated with an increase in blood sugar.	Regulatory Normalizing	Within 1 day	6-12 months
Healthy joints	Inflammatory changes in the shell of joints and periarticular tissues, hypothermia. Physical overload of joints of a professional, domestic or sports nature. Joint injuries.	Anti-inflammatory Anesthetic	Within an hour	3-6 months
Stop prostatitis	Dysfunction of the prostate gland, decreased functional activity of the testicles, impaired urination. With irregular sexual relations,	Anti-inflammatory Regulatory	2-3 days	6-12 months
Recovery after training	Post-training period. To relieve general discomfort, headaches, lethargy, drowsiness, irritability, reduce the concentration of lactic acid, normalize respiratory processes.	Regulatory Restores energy Relieves from spasms Antitoxic	Within an hour	1 month
Functional memory +	Decrease in the speed of memorization, concentration of attention, the duration of information storage.	Stimulating	5-7 days	1 month



Balanced nervous system	Stress, irritability, emotional instability, high physical and psychological stress (study, work, sport)	Anti-stress Calming	Within an hour	1 month
Nervous System - General	Functional problems of the nervous system, a violation of the integrity of the tissue and the connection between the departments of the central nervous system, leading to a violation of psychological stability.	Anti-inflammatory Regenerating Regulatory Anti-stress	Within 1 month	3-4 months
Healthy Ears	Pain and discomfort in the ear area arising from contact with infectious agents, water, pressure drop. Also in case of hearing loss.	Regulatory Relieves from spasms Anti-inflammatory Anesthetic	Within 1 day	1 month
Throat Without Problems	Pain and swelling in the throat, nasopharynx, feeling of stuffiness, hoarseness or loss of voice arising from the inflammatory process.	Anti-inflammatory Anesthetic Decongestant	1-2 days	Within 1 month
Strong Teeth	Bleeding, friability of the gums, violation of the integrity of the mucous membranes (microtrauma, dental procedures).	Anti-inflammatory Regenerating Normalizes microcirculation	5-7 days	Within 1 month
Healthy Pancreas	Violation of the externally secretory (digestive) or within the secretory (endocrine) functions of the gland.	Normalizing Regulating Anti-inflammatory	1 month	6 - 12 months
Healthy Thyroid	Violation of metabolic processes associated with the functional imbalance of the gland.	Normalizing Regulating Anti-inflammatory	1 month	6 - 12 months

DeVita Ritm (automatic programs)				
Name of the program	In which cases the automatic program can be used	Effect	First results	Lasting effect
Healthy Day - automatic	Improves blood supply to organs and tissues, activates the work of excretory systems, improves metabolism, the correct distribution of energy balance, relieves nervous tension, strengthens the protective functions of the body.	Immunostabilizing Gives energy Anti-stress Antitoxic Normalizes the hormonal background	5-7 days	1 month
Health - automatic	Corrects the state of the immune system for colds and viral diseases, improves cellular metabolism, which contributes to faster cleansing of toxins. It also restores the work of systems and organs that are affected during the illness.	Anti-inflammatory Immunostabilizing Relieves from spasms Decongestant	5-7 days	1 month
Complete cleansing - automatic	Activates the work of excretory systems, removes toxins, increases adaptive reserves and the defense mechanism of the organism. The daily biological activity of organs was taken into account.	Gives energy Antitoxic Decongestant Normalizes the hormonal background	5-7 days	1 month

High tonus - automatic	Relieves the toxic load from the liver, kidneys and lymphatic system, activates their work, improves blood circulation. Starts natural rejuvenation processes.	Anti-inflammatory Improves urine output Decongestant Immunostabilizing	Within 1 month	6-12 months
Healthy metabolism - automatic	Harmonizes the work of the organism, regulates the activity of the autonomic nervous system, stabilizes the appetite, activates the digestive and excretory organs. It contributes to the normalization of metabolic processes and psycho-emotional state, weight loss and improves the quality of life.	Anti-inflammatory Normalizes metabolism Decongestant Regulates appetite	Within 1 month	6-12 months
No poisoning - automatic	Reduces the harmful effects of toxins on the organism (food, drug, alcohol). Neutralize and removes them from the intestines and liver, reducing the burden on these organs.	Antitoxic Protects the liver Improves the outflow of bile Relieves from spasms	Within 1 day	5-7 days
Normal blood pressure - automatic	Relieves from spasm of blood vessels, normalizes blood pressure. Restores the work of heart and kidneys, normalizes cellular metabolism, reduces hypoxia and removes toxins. Improves blood supply to organs, reduces swelling.	Normalizes blood pressure Relieves from spasms Improves urine output	Within 1 day	1 month
Healthy sleep - automatic	Relieves from psychological stress, relaxes the central nervous system, improves the quality and duration of sleep.	Anti-stress Calming	Within 1-2 hours	1 month
Healthy joints - automatic	Relieves from inflammation and pain in the joints and spine. Improves their blood supply and nutrition. Restores damaged joint tissues, including articular cartilage and synovial membrane, as well as surrounding tissues: ligaments, tendons and muscles. Effective for any problems with joints: arthritis, arthrosis, stiffness	Anti-inflammatory Anesthetic Decongestant Regulatory Supports the function of joints	Within 1 day	3-6 months
Joints without pain - automatic	Relaxes tense muscles, improves the condition of the ligamentous apparatus, activates blood supply, and reduces pain and inflammation in the joints and spine. Effective with active load on the musculoskeletal system, injuries.	Anti-inflammatory Decongestant Anesthetic	Within 2-3 hours	1 month
Healthy skin and hair - automatic	Improves blood circulation to the skin, corrects the functioning of the sebaceous and sweat glands, which enhances skin immunity and reduces the manifestation of bacterial and fungal infections. It also improves blood circulation in the capillaries of the hair follicles and the bed of nail plates, contributing to enhanced nutrition of hair and nails.	Regulates the immune and barrier functions of the skin, the function of the glands Normalizes cell nutrition	5-7 days	1 month

Anti cold - automatic	This program is used during colds and allergic reactions (hay fever). Relieves from inflammation and reduces the respiratory manifestation of colds and allergies: swelling, nasal congestion, cough.	Anti-inflammatory Decongestant Immunostabilizing Relieves from spasms	1-2 hours	1 month
Men's Health - automatic	Contributes to the normalization of the prostate gland, the overall improvement of the genitourinary system in men, increases activity and sexual stamina. It also prevents problems associated with the impaired function of the organs of the genitourinary system.	Anti-inflammatory Expands blood vessels Hormone-normalizing Anti-aging	2-3 days	6-12 months
Without gout - automatic	Normalizes the metabolism in the organism, regulates the balance of K, Ca, Na, P. Relieves from inflammation in the joints, reduces swelling, and removes salts and toxins. Improves the functional state of the kidneys.	Anti-inflammatory Anesthetic Decongestant Improves urine output Regulates salt metabolism	1-2 days	3-6 months

## EXPLANATIONS FOR THE USE OF PROGRAMS DEVITA AP

DeVita AP (manual programs)				
Name of the program	In which cases the program can be used	Effect	First results	Lasting effect
No staphylococcus	Purulent wound infections, infection, purulent inflammatory diseases of various organs (upper respiratory tract, throat, skin, musculoskeletal system, external genital organs, urinary system, meninges, foodborne toxicosis). Prevention of fetal infection during pregnancy.	Antibacterial Anti-inflammatory	2-3 hours	10-14 days
No streptococcus	Diseases of the ear, throat, nose, inflammatory processes of the musculoskeletal system, heart, kidneys, skin. Prevention of fetal infection during pregnancy.	Antibacterial Anti-inflammatory	2-3 hours	10-14 days
No chlamydia	Infections of the genitourinary system, joints, eyes, upper and lower respiratory tracts, as well as the cardiovascular system. Prevention of occupational diseases among people working with birds. As well as prevention of infection through water during bathing and transmission of the infection to the fetus during pregnancy.	Antibacterial Anti-inflammatory	1-2 days	10-14 days
No helicobacter	Inflammation of the stomach, duodenum, damage to the mucous membranes of these organs.	Antibacterial Anti-inflammatory	5-7 days	3-5 weeks
No bacteria	Problems caused by Klebsiella: inflammation of the respiratory system, urinary tract, meninges, eyes, intestinal infections, blood poisoning. Problems caused by listeria: inflammation of the digestive tract, meninges, joints, respiratory organs. In order to prevent occupational diseases of people involved in animal husbandry, poultry farming, primary processing of meat and skins. Problems caused by pneumococcus: inflammation of the respiratory tract, ears, joints, meninges, blood poisoning.	Antibacterial Anti-inflammatory	Within 1 day	14-21 days
No hepatitis A	Prevention of the virus infection when visiting countries of Asia and Africa, reduces the risk of complications.	Antiviral Anti-inflammatory	1-2 weeks	6 months
Without toxins	Mandatory program for any toxic load (viral, bacterial, protozoal, mycotic, parasitic)	Antitoxic Energy saving	1-2 hours	1 month
No herpes zoster	Shingles - skin rashes, accompanied by severe pain along the sensory nerve.	Antiviral Anti-inflammatory Anesthetic	Within 1 day	1 month
No herpes simplex	Genital herpes and herpes (cold sore) on the lips.	Antiviral Anti-inflammatory	Within 1 day	1 month

No herpes	As an addition to any of the antiherpetic programs to reduce the time of recovery and reduce inflammatory and pain manifestations.	Antiviral Anti-inflammatory Antipruritic	Within 1 day	1 month
Healthy lymph	Disorders in the lymphatic system manifested by the following symptoms: swollen fingers and a tendency to swelling, chronic fatigue, constipation, discolored skin on the legs, chest pain and swelling, dry irritated skin, inability to concentrate, cold hands and feet, cellulite.	Anti-inflammatory Decongestant Antitoxic	Within 1 day	1 month
Anti-Cold	Diseases of the nose, throat, bronchi. In conjunction with the programs: "No Staphylococcus", "No Streptococcus", "No Bacteria".	Antiviral Anti-inflammatory Decongestant	2-3 hours	5-7 days
No mould & fungi	Lesions of the skin, subcutaneous fat, mucous membranes, visual and hearing apparatus, bones, joints, as well as internal organs (bronchopulmonary system, liver, spleen, heart, blood vessels, gastrointestinal tract, kidneys, central nervous system), allergic reactions in the form of cough, nasal congestion, skin rashes.  Effective in prevention of fungal infections after visiting the sauna, pool, beach, with prolonged exposure to damp rooms.	Antifungal Anti-inflammatory	1-2 days	2 months
No candida	Diseases of the mucous membranes of the internal organs, skin, accompanied by rashes, dysbiosis, burning sensation in the genitals.	Antifungal Anti-inflammatory	2-3 days	1 month
No fungi	Damage to the skin, hair, nails, mucous membranes of the internal organs by fungi of the genus: trichophyton, microsporum, mucor. Effective in prevention of fungal infections after visiting the sauna, pool, beach, getting in contact with an infected person.	Antifungal Anti-inflammatory Antipruritic	2-3 days	3-6 months
No protozoans	Inflammatory processes of the digestive tract caused by a number of protozoa: giardia, amoeba, etc. Prevention of infection through water when swimming.	Anti-protozoan Anti-inflammatory	2-3 days	1 month
Intimate health	Diseases of the mucous membrane of the respiratory tract, genitourinary tract caused by pathogens such as: ureaplasma, mycoplasma, toxoplasma. Prevention of transmission of infection to the fetus, sexual partner.	Antibacterial Anti-protozoan Anti-inflammatory	2-3 days	1 month
No schistosomes (blood-flukes) & trematods	Damage to the liver, pancreas, lungs, circulatory and genitourinary systems, toxic-allergic reactions caused by a number of flukes: feline, hepatic, pulmonary, lanceolate, as well as schistosomes.	Antiparasitic	1-5 days	1 - 2 months

No nematodes and ascarides	Infection with roundworms, such as ascarides, pinworms, etc. Causing weight loss, peeling of the skin, rash, hair loss, bowel obstruction, inflammation of the respiratory tract, liver and pancreas. The program fights parasites at different stages of their life cycle.	Antiparasitic	2-3 days	1-2 months
No parasites, General	Suppression of the vital activity of various types of parasites. Used in combination with other helminth extermination programs.	Antiparasitic	2-3 days	1 - 2 months
No parasites, Digestive tract	Parasites of the liver and pancreas, duodenum ulcer. Used in combination with other helminth extermination programs.	Antiparasitic	2-3 days	1 -2 months
Anti-Flu	The program aims to destroy various strains of influenza viruses (A, B, C).	Antiviral	2-3 hours	5-10 days
Against flu - General	Mutated forms of flu, fever, pain.	Antiviral Antipyretic Anesthetic Anti-inflammatory	1-2 hours	5-10 days
Healthy gastrointestinal	Normalizes the condition and stops dehydration in case of dysbacteriosis caused by rotavirus infection, E. coli.	Antibacterial Antiviral Anti-inflammatory Protects cells Normalizes metabolism	2-3 hours	5-7 days
Help against poisoning	Normalizes the condition with Enterovirus infection and Salmonella infection.	Antibacterial Antiviral Anti-inflammatory Protects cells Normalizes metabolism Antitoxic	2-3 hours	10 days
Antivirus	Reduces the total viral load, cytomegalovirus infection, proceeding in the form of a persistent cold with headaches and joint pain, inflammation of the salivary glands.	Antiviral Anti-inflammatory	2-3 hours	1 month
Drainage, General	Infections of the upper respiratory tract, lymphoid tissue, liver, spleen caused by the Epstein-Barr virus.	Antiseptic Antitoxic Lymphatic drainage	Within 1 day	1 month
Drainage Antitox	With chemical toxic loads of the tissue sector with salts of heavy metals, radionuclides, food toxins.	Antitoxic Lymphatic drainage	1-2 days	14-21 days
Mexican detoxification	After surgical interventions, with chronic and drug intoxication.	Lymphatic drainage Decongestant Antitoxic Stabilizing Regulates excretory functions	6-7 days	1 month
No allergies	Various types of allergies (respiratory - odors, dust, pets, hay fever, skin, food), stabilizes the immune system.	Antitoxic Antiallergic Decongestant Anti-inflammatory Immunostabilizing	1-2 days	4-6 months

Weight Control	Presence of excess weight, impaired appetite, destroys one of the causes of obesity - Ad-36 adenovirus.	Antiviral Antitoxic	2-4 weeks	4-6 months
Against Lyme Disease 1	Tick bites. Those programs are used in pairs. If the bite happened earlier, you must take the full application scheme of the complex "Lyme Disease".	Antibacterial	1-2 days	14 – 28 days
Against Lyme Disease 2		Antiviral		
Against Cold - General	Diseases of the nose, throat, bronchi. Use together with the programs: "No staphylococcus", "No streptococcus", "No bacteria".	Decongestant Anti-inflammatory Anesthetic	1-2 days	1 month
Anti-inflammation - General	With any inflammatory processes in the body.	Decongestant Anti-inflammatory Anesthetic	1-2 days	1 – 2 months
Antiseptic	After surgical interventions, with any damage to the integrity of the skin and mucous membranes, with inflammatory processes.	Antiseptic Decongestant Anti-inflammatory Anesthetic Accelerates regeneration	2-3 hours	7 – 21 days
Anti-virus Gastrointestinal	Prevention of Coxsackie virus disease when visiting countries of Asia and Africa, reducing the risk of complications. The following symptoms can indicate the presence of the virus: dehydration, sore throat, the presence of a vesicle rash on the skin.	Antiviral Protects cells Normalizes metabolism	1-2 days	14 – 21 days
Clean Sinuses	Acute inflammation of the nasal mucosa and sinuses (sinusitis, sinusitis, etc.)	Decongestant Anti-inflammatory Anesthetic	1-2 days	14 - 21 days
Clean Sinuses - General	Together with antiviral and antibacterial programs.	Decongestant Anti-inflammatory Anesthetic Stabilizing	1-2 days	14 days
CNS Without Problems	For chronic problems with the nasal mucosa and sinuses. As one of the means of preventing an acute process in the upper respiratory tract for colds.	Decongestant Anesthetic Protects cells Revitalizing	2-3 days	6-12 months
Ears Without Inflammation	Pain and discomfort in the ear area arising from contact with infectious agents, water, pressure drop.	Decongestant Anti-inflammatory Anesthetic	1-2 days	14 - 21 days
Healthy Bile	Violation of the timely and uniform separation of bile into the intestine, which may be a consequence of irregular nutrition, consumption of harmful products, excessive consumption of fats, medicines, and the presence of parasites.	Anti-inflammatory Anesthetic Cholagogic	2-3 hours	1-4 months
Healthy Cell - 1	Prevents the development of cell mutations, reduces the negative effects of the environment and toxic load. Particular attention should be paid to those who live in adverse ecological areas, experience great emotional stress and the influence of negative production factors. Programs should be used in the "Healthy Cell" complex.	Antiviral Anti-inflammatory Protects cells	3-4 days	6-12 months
Healthy Cell - 2				

Healthy Gums	Presence of dentures, braces, crowns, as well as infections of the oral cavity that cause gum disease.	Anti-inflammatory Anesthetic Antitoxic Improves nutrition	2-3 hours	1 month
Healthy Mouth	Inflammation of the oral mucosa (sores, mucosal damage) that occurs when there are infectious agents, damage to the oral cavity, during dental interventions, an unbalanced diet.	Antiviral Anti-inflammatory Anesthetic Accelerates regeneration	1-2 days	1 month
Stop Injury (teeth)	The best effect will be acquired if used in combination with the program "Mucous Membrane Regeneration".	Anesthetic Decongestant Accelerates regeneration	2-3 hours	7-14 days
Mucous Membrane Regeneration	Inflammation of the oral mucosa (sores, mucosal damage) that occurs when there are infectious agents, damage to the oral cavity, during dental interventions, an unbalanced diet.	Anti-inflammatory Accelerates regeneration Immunostabilizing	1-2 days	1 month
Without Fistula	The best effect will be acquired if used in combination with the program "Healthy Mouth".	Anti-inflammatory Decongestant Antitoxic	1-2 days	7-14 days
Stop Adenovirus - General	After dental interventions (removal, treatment), in the presence of swelling of the gums, soreness.	Antiviral Anti-inflammatory	1-2 days	14-21 days
Healthy Nerves	Inflammation of the mucous membranes of the nose, throat, eyes, intestines.	Anti-inflammatory Protects cells	2-3 days	6-12 months
Stop Amoeba	Normalizes the condition with dysentery, gastroenteritis, intestinal ulcers, caused by the presence of protozoa. It is recommended to use to prevent infection when visiting tropical and subtropical countries of Asia, Africa, South America.	Anti-protozoan Anti-inflammatory Normalizes metabolism	2-3 days	1 month
Stop Coronavirus	Inflammation of the mucous membranes of the nose, throat, lungs.	Antiviral Anti-inflammatory	1-2 days	14-21 days
Without Klebsiella	Inflammation of the respiratory system, urinary tract, meninges, eyes, intestinal infections, blood poisoning. It is recommended as an addition for viral load, malnutrition, also for people with reduced immunity, people suffering from alcoholism, diabetes, chronic lung diseases, liver diseases.	Antibacterial Anti-inflammatory Antitoxic	1-2 days	14-28 days
Stop Migraine	Headaches of various origins.	Anesthetic Antispasmodic Normalizes blood circulation	Within 1 day	1 month
Stop Penicillium	Reduces toxic infections when using products affected by the fungus, in case of chronic diseases of the respiratory tract, kidneys, allergic reaction to penicillium.	Antifungal Anti-inflammatory	2-3 days	1 month



Throat Without Pain	Pain and swelling in the throat, nasopharynx, feeling of stuffiness, hoarseness or loss of voice arising from the inflammatory process.	Anti-inflammatory Decongestant Anesthetic Antitoxic	Within 1 day	14-21 days
Without Campylobacter	As an additional method to prevent infection and the development of diarrhea, for people working in livestock and poultry farms, as well as tourists visiting developing countries, people who use contaminated meat products, milk, water.	Antibacterial Anti-inflammatory	1-2 days	14-28 days
Without Clostridia	In case of damage to the integrity of the skin with contamination (cuts, bites, other injuries). When using expired or doubtful foods, canned foods, etc. Alternate with the "Without toxins" program throughout the day. Make sure to contact/see a doctor!	Antibacterial Anti-inflammatory	1-2 days	14-28 days
Without Cryptococcus	As one of the methods for preventing the disease, for people working in poultry farms engaged in breeding birds in an apartment or on their own farm. Particular attention is given to people with reduced immunity and chronic lung diseases located in places where birds gather.	Antifungal Anti-inflammatory	2-3 days	1 month
Without Enterococcus	Frequent and painful urination, pain during bowel movements, decreased potency, increased fatigue, depression, gastritis, dysbiosis, food poisoning.	Antibacterial Anti-inflammatory Antitoxic	1-2 days	14-28 days
Without Echinococcus	Prevention of helminthic invasion and damage to internal organs to people raising animals, livestock, collecting berries, fruits, vegetables, hunting. Especially living in South America and Europe, North Africa, New Zealand, Australia, in Russia in areas: Volga region, Western Siberia, Krasnodar Territory, in the Far East.	Antiparasitic	2-3 days	1-2 months
Without Leptospira	Jaundice, hemorrhagic fever, occupational diseases of pastoralists and agricultural workers.	Antibacterial Anti-inflammatory	1-2 days	14-28 days
Without nematodes and ascarididae - New	Infection with roundworms, such as ascariides, pinworms, etc. Causing weight loss, peeling of the skin, rash, hair loss, bowel obstruction, inflammation of the respiratory tract, liver and pancreas. An enhanced, expanded program fights parasites at different stages of their life cycle.	Antiparasitic Antitoxic	2-3 days	1-2 months
Without Kidney Papilloma	HPV damage to the mucous surface of the urinary system (renal pelvis and bladder).	Antiviral Protects cells	2-3 days	1 month
Without Papilloma - General	Damage to the skin and mucous membranes caused by the Human Papilloma Virus. It is recommended to use with the "Clean Skin" complex.	Antiviral Protects cells	2-3 дня	1 month

Without Proteus	Gastritis, diarrhea, abdominal pain, vomiting, exacerbation of dermatitis, allergies, problems in the genitourinary system, long non-healing wounds.	Antibacterial Anti-inflammatory Antitoxic	1-2 days	14-28 days
Without Pseudomonas aeruginosa	Local purulent-inflammatory processes in the soft tissues, organs of hearing, lungs, joints, etc. Also for people who are in hospital for a long time.	Antibacterial Anti-inflammatory Antitoxic	1-2 days	14-28 days
Without Rickettsia	Tick bites, especially in Asia and Australia, typhus, tick-borne fevers, etc.	Antibacterial Anti-inflammatory Protects cells	1-2 days	14-28 days
Without Rotavirus	Normalizes the condition and stops dehydration in case of dysbacteriosis caused by rotavirus infection.	Antiviral Protects cells	1-2 days	14-21 days
Without Salmonella	Normalizes the condition during toxic infection, which appeared as a result of eating foods with pathogenic microorganisms.	Antibacterial Anti-inflammatory Antitoxic	1-2 days	14-28 days
Without Shigella	Normalizes the condition and stops dehydration with diarrhea caused by the presence of this bacterium.	Antibacterial Anti-inflammatory Protects cells	1-2 days	14-28 days
Without Toxoplasma	Parasitic invasion caused by intracellular protozoa and accompanied by the development of the inflammatory process in the lymph nodes, liver, brain, lungs, heart, muscles, etc.	Anti-protozoan Anti-inflammatory	2-3 days	1 month
Without Trematodes - General 1	Damage to the liver, pancreas, lungs, circulatory and genitourinary systems, toxic-allergic reactions caused by a number of flukes: feline, hepatic, pulmonary, lanceolate, etc. Programs are used together with each other.	Antiparasitic Protects cells Antitoxic	2-3 days	1-2 months
Without Trematodes - General 2				
Without Trichomonas	Inflammatory processes of the genitourinary sphere, pain during intercourse, itching and burning, discharge and unpleasant odor.	Antiprotozoal Anti-inflammatory	2-3 days	1 month
Without Ureaplasma	Inflammatory processes of the genitourinary sphere, impaired reproductive function, pain in the lower abdomen, rapid and painful urination.	Anti-protozoan Anti-inflammatory	2-3 days	1 month
Without Yersinia	Toxic-allergic reaction that occurs when hygiene measures are not followed and accompanied by skin rashes, a violation of the digestive tract, inflammation of the joints.	Antibacterial Anti-inflammatory	1-2 days	14-28 days
Stop Coronavirus - New	Inflammation of the mucous membranes of the nose, throat, lungs, weakness, decreased sense of smell and hearing. Improved and advanced program. It is advisable to use according to the methodology of the "Stop Coronavirus" complex.	Antiviral Protects cells Anti-inflammatory	1-2 days	1 month

Bacterial pneumonia	Inflammation of the respiratory system caused by Streptococcus pneumoniae, Chlamydomphila pneumoniae, Mycoplasma pneumoniae, Klebsiella pneumoniae, etc., accompanied by fever, weakness, headache, cough, shortness of breath, pain in the chest, muscles and joints.	Antibacterial Anti-inflammatory	1-2 days	1 month
Without cough	In the presence of both dry and wet cough with various infections. To get the best effect, it is recommended to use in conjunction with the "Anti-cold" and "Free Breath" programs on the DeVita Ritm device.	Anti-inflammatory Decongestant Relieves from spasms	Within 1 day	1 month
Without pneumonia	In the presence of respiratory tract infection, cough, general weakness, fatigue, drowsiness, shortness of breath, chest pain. It is recommended to use with other anti-inflammatory programs aimed at normalizing the function of the respiratory system.	Anti-inflammatory Decongestant	2-3 days	1 month

DeVita AP (automatic programs)				
Name of the program	In which cases the program can be used	Effect	First results	Lasting effect
Mode 1 - automatic	Viral and protozoan loads. Optimal disease prevention when changing climatic zones, visiting Asian countries.	Antiviral Anti-inflammatory Anti-protozoan Antitoxic	1-3 days	2-4 weeks
Mode 2 - automatic	Toxic stress caused by streptococcus and staphylococcus. Effective both for colds and for the prevention of coccal infections.	Antibacterial Anti-inflammatory Antitoxic	2-3 hours	10-14 days
Mode 3 - automatic	Bacterial burden in the digestive tract and the genitourinary system.	Antiviral Anti-inflammatory Antibacterial Antitoxic	1-3 days	2-4 weeks
Mode 4 - automatic	Presence of mycotoxins in the organism formed by mold and yeast. Reduces the risk of accumulation of carcinogenic factors.	Antifungal Anti-inflammatory Antitoxic	2-3 days	1-2 months
Cleansing - automatic	Violation of the lymphatic system, swelling, decreased immunity.	Antiseptic Antitoxic Lymphatic drainage Decongestant	Within 1 day	1 month
No parasites - automatic	Parasitic burden in the liver, pancreas, intestines, as well as toxic load caused by the vital activity of these microorganisms.	Antiparasitic Antitoxic	2-3 days	1 - 2 months
No flu - automatic	Inflammatory processes caused by various strains of the influenza virus.	Antiviral Antipyretic Anesthetic Anti-inflammatory Antitoxic	2-3 hours	5-10 days

Stop herpes S - automatic	In case of vesicles forms on the mucous membranes of the lips, genitals, eyes, mouth, skin, which are formed due to activation of the Herpes Simplex virus.	Antiviral Anti-inflammatory Antipruritic Antitoxic	Within 1 day	1 month
Stop herpes Z - automatic	Shingles - skin rashes, accompanied by severe pain along the sensory nerve.	Antiviral Anti-inflammatory Antipruritic Anesthetic Antitoxic	Within 1 day	1 month
Without cold - automatic	Inflammation of the nose, throat, bronchi.	Antiviral Anti-inflammatory Antitoxic	Within 1 day	5-7 days
Support immune system - automatic	Effective for preventing colds infections during hypothermia and during mass illnesses.	Antiviral Antitoxic	2-3 hours	5-7 days
Without helminths - automatic	Effective for the prevention of helminthiasis in children and adults.	Antiparasitic Antitoxic	2-3 days	1 - 2 months
Weight Reduction - automatic	Excess weight, resulting from viral complications, Puffiness.	Antiviral Antitoxic	1-2 days	4-6 months
Reduce Points - automatic	Violation of the lymphatic system, fluid retention, decreased immunity due to increased intoxication.	Antitoxic Lymphatic drainage Decongestant	1-2 days	4-6 months
Cell Protection 1	For removal of the general viral, herpetic and bacterial toxic loads. Reduces the activity of viruses and the negative effect of streptococcal toxins.	Antiviral Antibacterial Antitoxic	3-4 days	4-6 months
Cell Protection 2	Wrong lifestyle, great stress and physical exertion, living in adverse ecological areas. It is recommended to use the program within the application scheme of the "Healthy Cell" complex.	Antiviral Anti-protozoan Anti-inflammatory Protects cells	Within 1 week	6-12 months
Cell Protection 3	The program is effective against staphylococcal and chlamydial infections, and also reduces the negative impact of protozoan lesions of the genitourinary tract.	Antiviral Antibacterial Anti-protozoan Anti-inflammatory	3-4 days	4-6 months
Cell Protection 4	Fungal load and in order to prevent mutations and negative cell degeneration. It is recommended to use the program within the application scheme of the "Healthy Cell" complex.	Antifungal Anti-inflammatory Antitoxic Protecting cells	Within 1 week	6-12 months
Lyme 1 - automatic	Removal of the general viral, herpetic and bacterial toxic load. Reduces the activity of viruses and the negative effect of streptococcal toxins.	Antiviral Antibacterial Antitoxic	3-4 days	4-6 months
Lyme 2 - automatic	Tick bites, the detection of tick-borne borreliosis or Lyme disease. It is recommended to use the program within the application scheme of the "Lyme disease" complex.	Antiviral Antibacterial Anti-protozoan Anti-inflammatory	3-4 days	6-12 months

Lyme 3 - automatic	The program is effective against staphylococcal and chlamydial infections, and also reduces the negative impact of protozoan lesions of the genitourinary tract.	Antiviral Antibacterial Anti-protozoan Anti-inflammatory	3-4 days	4-6 months
Lyme 4 - automatic	Tick bites, the detection of tick-borne borreliosis or Lyme disease. It is recommended to use the program within the application scheme of the "Lyme disease" complex.	Antibacterial Antifungal Anti-inflammatory Antitoxic	3-4 days	6-12 months
Restoration 1 - automatic	Inflammation and damage to the mucous cavity, with dental interventions, an unbalanced diet.	Anesthetic Decongestant Anti-inflammatory Accelerates regeneration	1-2 days	1 month
Restoration 2 - automatic	After surgical interventions, with any damage to the integrity of the skin and mucous membranes, with inflammatory processes.	Anesthetic Decongestant Anti-inflammatory Accelerates regeneration	Within 1 day	2-4 weeks
Healthy Gums - Automatic	Infections and injuries of the oral cavity. It is recommended to use according to the application scheme of the "Healthy Teeth" complex	Decongestant Anti-inflammatory Accelerates regeneration	1-2 days	4-6 months
Nerve Regeneration - automatic	Negative disorder of the nervous system, a violation of the integrity of the nerve fiber caused by tick bites or other factors.	Anti-inflammatory Decongestant Anesthetic Protects cells	Within 1 week	6-12 months
Healthy Breathing - automatic	Diseases of the respiratory tract: nose, throat, bronchi, lungs.	Antiviral Antibacterial Antitoxic Anti-inflammatory	1-2 days	1 month

## THE BIOLOGICAL ROLE OF DELIXIR

Organ system	DeLixir	Biological Role
Cardiovascular system	Multi-Energy	Regulates the functions of blood coagulation, hematopoiesis (for women during menstruation), myocardium, relieves from swelling, prevents thinning and fragility of blood vessels. Normalizes blood pressure, improves blood composition.
	pH Balance	Strengthens blood vessels, prevents the development of varicose veins, hemorrhoids, activates the process of hematopoiesis, regulates blood pressure.
	Collagen+	Regulates the process of hematopoiesis (blood formation), supports the work of the heart muscle, lowers cholesterol, and improves blood circulation.
	Detox	Prevents blood thickening, normalizes blood pressure, cleanses the walls of blood vessels and strengthens them.
Gastrointestinal tract	Multi-Energy	Relieves from inflammation of the mucous membrane, improves the process of digestion and assimilation of nutrients, regulates the secretory function.
	pH Balance	Normalizes the acidity of gastric juice, improves digestion, prevents the development of intestinal diseases, inhibits fermentation, and stabilizes the flora.
	Collagen+	Improves the functioning of the gastrointestinal tract, normalizes appetite.
	Detox	It facilitates the digestion of proteins, helps with disorders of the intestines (atony, constipation) and indigestion. Removes toxins, suppresses pathogenic flora.
Nervous system	Multi-Energy	Improves brain function, relieves from nervous tension. Favorably affects the processes of inhibition and excitation in the cerebral cortex. Improves mood.
	pH Balance	It contributes to the formation of stress resistance, helps with irritability, sleep disturbances.
	Collagen+	It improves the state of depression, relieves from fatigue, and has a positive effect on brain activity.
	Detox	Relieves from nervous tension, improves mood.
Endocrine system, metabolism	Multi-Energy	It activates tissue respiration reactions and prevents disturbances in carbohydrate metabolism. In particular, pyruvic acid, which reduces tissue acidosis and hypoxia. Protects cells from destruction. Stimulates cell metabolism. Increases vitality. Facilitates the start of a diet for weight loss. Lowers cholesterol. Improves protein and lipid metabolism. Regulates the hormonal background.
	pH Balance	Stimulates carbohydrate and protein metabolisms, improves the absorption of calcium, iron, activates the metabolism, promotes weight loss. Lowers cholesterol, stabilizes sugar levels.
	Collagen+	Regulates protein and fat metabolism, slows down aging.
	Detox	It accelerates protein biosynthesis, cell metabolism, carbohydrate metabolism, inhibits the aging process, protects against free radicals.
Reproductive system	Multi-Energy	It normalizes hormonal levels (infertility), enhances potency, increases libido.
	pH Balance	Soothes menstrual ailments.
	Collagen+	Normalizes hormonal balance.
	Detox	Enhances sexual activity.
Urinary system	Multi-Energy	Removes water, eliminates swelling.
	pH Balance	Mild diuretic effect, prevention of urolithiasis.
	Collagen+	Reduces swelling.
	Detox	It removes excess water, eliminates swelling. Soothes kidney and bladder problems.

Respiratory system	Multi-Energy	Soothes asthma attacks.
	pH Balance	Enhances local immunity.
	Collagen+	Reduces inflammation, restores vascular elasticity.
	Detox	Strengthens the mucous membranes.
Skin / mucous membranes	Multi-Energy	It speeds up the regeneration of the skin, participates in the formation of collagen in the skin and connective tissue, improves the quality of hair and skin, strengthens the mucous membrane, eliminates bleeding gums, nosebleeds. Slows down the aging process.
	pH Balance	It activates the renewal of connective tissue, the growth of hair and nails, reduces bleeding gums, protects the mucous membranes.
	Collagen+	Accelerates skin regeneration, supports skin immunity (dermatitis).
	Detox	Revitalizes the color of hair and skin, strengthens the mucous membrane, reduces bleeding gums.
Musculoskeletal system	Multi-Energy	It improves the absorption of calcium and phosphorus, prevents the deformation of the pelvic bones in adolescent girls, a decrease in bone strength in adults, and muscle atrophy. Activates their development.
	pH Balance	Participates in bone and teeth strengthening, in renewing the connective tissue, and activates muscle work.
	Collagen+	Maintains the structure and function of bones, joints and muscles.
	Detox	Relieves from muscle spasms, strengthens the connective tissue, reduces joint pain.
Immune system	Multi-Energy	Restorative effect.
	pH Balance	Prevents colds, acidosis with infections.
	Collagen+	It has an anti-inflammatory, general strengthening effect.
	Detox	Protects against colds, restores the barrier function of the mucosa, inhibits inflammatory processes.
Organs of vision	Multi-Energy	It is necessary for the formation of visual pigment, the formation of visual acuity, and the reduction of dryness of the cornea and mucous membranes.
	pH Balance	Improves vision.
	Collagen+	Strengthens the muscles of the eye, inhibits the process of vision reduction.
	Detox	Increases visual acuity, reduces night blindness. Relieves from swelling and inflammation of the eyelids.
Energy in the body / Rejuvenation	Multi-Energy	Gives strength and energy. Protects cells from destruction (antioxidants).
	pH Balance	Reduces fatigue, promotes rejuvenation (cell renewal).
	Collagen+	Relieves from fatigue, improves performance.
	Detox	Reduces toxic load, slows down aging (antioxidants).

## RECOMMENDED DEVITA PROGRAMS FOR USE WITH VARIOUS FUNCTIONAL PROBLEMS

Organ system	Disorders of the organs of the system	Programs on DeVita Ritm	Programs on DeVita AP
Cardiovascular system	Cardiac / coronary insufficiency (coronary heart disease, angina pectoris, post-infarction state, toxic damage to the heart muscle, functional overload of the heart, etc.)	Healthy Day - automatic Health - automatic Complete cleansing - automatic Healthy sleep Healthy heart	Antivirus Mode 2 – automatic No chlamydia No schistosomes (blood-flukes) & trematods Cleansing - automatic Intimate health Anti-virus Gastrointestinal Restoration 2 - automatic
	Heart rhythm disorders (arrhythmia, tachycardia)	Healthy Day - automatic Health - automatic Complete cleansing - automatic Healthy heart Cardiac rhythm normalization Healthy sleep - automatic	Antivirus Mode 2 – automatic No chlamydia No schistosomes (blood-flukes) & trematods Cleansing - automatic Anti-virus Gastrointestinal Intimate health
	Hypertension, hypotension, varicose veins, atherosclerosis and other diseases of the vascular system	Healthy Day - automatic Health - automatic Complete cleansing - automatic Normal blood pressure - automatic Healthy venous system Without headache Digestion normalization Healthy sleep - automatic	Antivirus Mode 2 – automatic Mode 3 – automatic No chlamydia, No schistosomes (blood-flukes) & trematods Cleansing - automatic Anti-virus Gastrointestinal Without Yersinia Stop Migraine
Digestive system	Gastritis, esophagitis, heartburn, ulcer	Healthy Day - automatic Health – automatic Healthy sleep - automatic	No helicobacter Mode 4 – automatic No parasites - automatic Cleansing - automatic Restoration 2 - automatic
	Duodenitis, colitis, enteritis, ulcer, irritable bowel syndrome, dysbiosis, atony.	Healthy Day - automatic Health - automatic Complete cleansing - automatic	Mode 1 - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic No parasites – automatic Without nematodes and ascarididae - New Cleansing - automatic Without Salmonella Without Proteus Without Klebsiella Stop Amoeba Without Clostridia
	Hepatitis, biliary dyskinesia, cholecystitis, pancreatitis, diabetes, cholelithiasis.	Healthy Day - automatic Health - automatic Complete cleansing - automatic Digestion normalization Normal blood sugar Healthy Pancreas	Mode 1 - automatic Mode 3 - automatic No parasites – automatic Without nematodes and ascarididae - New No bacteria Cleansing - automatic Healthy Bile Without Trematodes – General 1.2 Without Salmonella Without Leptospira Restoration 2 - automatic



Digestive system	Acute diseases (poisoning, allergies)	Healthy Day - automatic Digestion normalization Without allergens No poisoning - automatic	Healthy Gastrointestinal Help against poisoning Drainage Antitox No allergies Cleansing - automatic Stop Penicillium Mexican detoxification Without Clostridia Without Shigella Without Proteus Without Rotavirus
Nervous system	Neurosis	Healthy Day - automatic Health – automatic Good mood Healthy sleep - automatic Functional memory +	Mode 1 - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic Cleansing - automatic Stop Migraine No parasites – automatic No bacteria No fungi
	Autonomic dysfunction	Healthy Day - automatic Health – automatic Good mood Healthy sleep - automatic Healthy heart Healthy kidneys Without headache	Mode 1 - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic Cleansing - automatic No parasites – automatic No bacteria No fungi No chlamydia Anti-virus Gastrointestinal Stop Migraine
	Stroke	Healthy Day - automatic Health – automatic Recovery from injuries Healthy sleep - automatic Healthy heart Without headache	Mode 1 - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic Cleansing - automatic No parasites – automatic No bacteria, No fungi No chlamydia Anti-virus Gastrointestinal Stop Migraine
	Radiculitis, neuritis, shingles	Healthy Day - automatic Health – automatic Recovery from injuries Healthy sleep – automatic Joints without pain – automatic Nervous System - General	Mode 1 - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic Cleansing - automatic No parasites – automatic No bacteria, No fungi Stop herpes Z - automatic Anti-virus Gastrointestinal Restoration 2 - automatic
Endocrine system	Hypothalamic-pituitary disorders	Healthy Day - automatic Health - automatic Complete cleansing - automatic Women's health, Men's health	Mode 1 - automatic Mode 4 - automatic Cleansing - automatic Drainage Antitox Anti-virus Gastrointestinal Without Cryptococcus

Endocrine system	Pancreatic disease	Healthy Day - automatic Health - automatic Complete cleansing - automatic Digestion normalization Normal blood sugar Healthy Pancreas	Mode 1 - automatic No parasites – automatic Cleansing - automatic Without Trematodes - General 1,2 Without Salmonella Restoration 2 - automatic
	Thyroid disease	Healthy Day - automatic Health - automatic Complete cleansing - automatic Women's health, Men's health Digestion normalization Healthy Thyroid	Mode 1 - automatic Mode 4 - automatic Cleansing - automatic Drainage Antitox Restoration 2 - automatic
	Adrenal disease	Healthy Day - automatic Health - automatic Complete cleansing - automatic Women's health, Men's health	Mode 1 - automatic Mode 2 - automatic Cleansing - automatic Drainage Antitox Restoration 2 - automatic.
	Genital Diseases	Healthy Day - automatic Health - automatic Complete cleansing - automatic Women's health, Men's health	Mode 1 - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic Cleansing - automatic No chlamydia Stop herpes S - automatic Anti-virus Gastrointestinal, Without Papilloma - General Without Trichomonas Without Enterococcus Restoration 2 - automatic
Reproductive system	Female Genital Disorders (cycle disorders, painful menstruation, menopause, neoplasms)	Healthy Day - automatic Health – automatic Women's health Female sexuality Complete cleansing - automatic	Mode 1 - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic Cleansing - automatic No chlamydia Stop herpes S - automatic Without Proteus Anti-virus Gastrointestinal Without Papilloma - General Without Trichomonas Without Enterococcus Restoration 2 – automatic.
	Male Genital Disorders (prostatitis, adenoma, problems with potency)	Healthy Day - automatic Health – automatic Men's Health - automatic Male sexuality Complete cleansing - automatic	Mode 1 - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic Cleansing - automatic No chlamydia Stop herpes S - automatic Without Proteus Anti-virus Gastrointestinal Without Papilloma - General Without Trichomonas Without Enterococcus Restoration 2 - automatic

Urinary system	Urethritis, Cystitis, Urolithiasis, gout, pyelonephritis	Healthy Day - automatic Health – automatic Complete cleansing – automatic Healthy joints - automatic Without gout - automatic	Mode 1 - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic Cleansing - automatic No chlamydia Stop herpes S - automatic Without Proteus Anti-virus Gastrointestinal Without Enterococcus Stop Penicillium Without Kidney Papilloma Without Trichomonas Restoration 2 - automatic
Respiratory system	Problems of the upper respiratory tract (acute respiratory infections, rhinitis, tracheitis, laryngitis, pharyngitis, sinusitis)	Healthy Day - automatic Health – automatic Without allergens Anti-cold	Without cold - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic Cleansing - automatic No allergies No bacteria Clean Sinuses Against Cold - General, Restoration 2 - automatic Anti-virus Gastrointestinal Healthy Breathing - automatic
	Lower respiratory problems (bronchitis, bronchial asthma, pneumonia)	Healthy Day - automatic Health – automatic Without allergens Anti-cold - automatic	Without cold - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic Cleansing - automatic No allergies No bacteria No parasites - automatic No nematodes and ascarides Healthy Breathing - automatic Without pneumonia Without cough Restoration 2 - automatic Without Pseudomonas aeruginosa
Skin / mucous membranes	Contact dermatitis, seborrheic dermatitis	Healthy Day - automatic Health – automatic Complete cleansing – automatic Healthy sleep – automatic Healthy skin and hair - automatic	Mode 1 - automatic Mode 2 - automatic Mode 4 - automatic No fungi No allergies Drainage Antitox Cleansing - automatic Restoration 2 - automatic
	Neurodermatitis, eczema	Healthy Day - automatic Health – automatic Complete cleansing – automatic Healthy sleep – automatic Without allergens Women's health, Men's health Beautiful skin	Mode 2 - automatic Mode 3 - automatic No allergies Drainage Antitox Cleansing - automatic Without Proteus Without Pseudomonas aeruginosa Restoration 2 - automatic

Skin / mucous membranes	Allergy	Healthy Day - automatic Health - automatic Complete cleansing - automatic Without allergens Anti-cold	Mode 4 - automatic No parasites - automatic Without nematodes and ascarididae - New No allergies Drainage Antitox Cleansing - automatic
Musculoskeletal system	Osteochondrosis, arthritis, arthrosis, gout, protrusion, trauma, bursitis, myositis	Healthy Day - automatic Health – automatic Complete cleansing – automatic Joints without pain - automatic Healthy joints – automatic Without gout - automatic Recovery from injuries	Mode 1 - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic Cleansing - automatic No chlamydia Drainage Antitox Against Lyme Disease 1,2 Without Pseudomonas aeruginosa Without Enterococcus
Immune system	Allergy, problems of the thymus, spleen, lymph nodes	Healthy Day - automatic Health – automatic Complete cleansing – automatic Without allergens Clean lymph Without headache	Without cold - automatic Mode 2 - automatic Mode 4 - automatic No parasites - automatic Without nematodes and ascarididae - New Drainage Antitox Cleansing - automatic Mexican detoxification No allergies Against Cold - General Throat Without Pain Restoration 2 - automatic Stop Penicillium
Sensory organs	Eye problems (optic atrophy, blepharitis, keratitis, conjunctivitis, cataracts, glaucoma, myopia, hyperopia, astigmatism)	Healthy Day - automatic Health – automatic Complete cleansing – automatic Healthy eyesight Without allergens Nervous System - General	Support immune system - automatic Mode 2 - automatic Mode 3 - automatic Cleansing - automatic Stop herpes S - automatic Restoration 2 - automatic No parasites - automatic Without nematodes and ascarididae - New No chlamydia No fungi No allergies Drainage Antitox
	Hearing problems (otitis media, hearing loss)	Healthy Day - automatic Health – automatic Complete cleansing – automatic Healthy Ears Nervous System - General Recovery from injuries Anti-cold	Without cold - automatic Mode 2 - automatic Mode 4 - automatic, Intimate health No allergies Anti-virus Gastrointestinal Ears Without Inflammation Clean Sinuses Against Cold - General Restoration 2 - automatic

Sensory organs	Problems of the sense of smell (acute respiratory disease, sinusitis, rhinitis, trauma)	Healthy Day - automatic Health – automatic Complete cleansing – automatic Nervous System - General Anti-cold Without allergens Recovery from injuries	Without cold - automatic Mode 2 - automatic Intimate health No allergies No bacteria Anti-virus Gastrointestinal Clean Sinuses - General Clean Sinuses Against Cold - General Stop Coronavirus - New Restoration 2 - automatic
Dental problems	Damage and inflammation of the mucosa, bleeding gums, pain after dental interventions	Healthy Day - automatic Health – automatic Complete cleansing – automatic Strong Teeth Recovery from injuries	Mode 1 - automatic Mode 2 - automatic Mode 4 - automatic No fungi Healthy Gums - Automatic Restoration 1-2 - automatic Stop Injury (teeth) Without Fistula Drainage Antitox

**Note:** the necessary programs for the BASE series devices can be recorded via the DelInfo software.





The combined use of all three DeVita devices is the key to a healthy and active life. Only by using all devices will your organism be completely cleansed of toxins, your organs will begin to function harmoniously, and you will feel cheerfulness and a surge of strength.

Also, the necessary smart cellular nutrition will establish vital processes.

The combined use of the technologies and smart cell nutrition supplements made by DEHolding is the key to health and longevity.

[WWW.DEHOLDING.INFO](http://WWW.DEHOLDING.INFO)

20200525

