Deta
Elis
Holding
Gesundheit
Innovationen
Mensch

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A wellness device for
the harmonization and
correction of natural
electromagnetic radiation of
natural biosystems

DeVita Ritm+
Model Mini

Intended use

The wellness device "DeVita Ritm+ Model MINI" is a generator of faint electromagnetic radiation, similar to the radiation of natural organic systems.

It is intended to strengthen our own protective forces and to correct natural electromagnetic radiation in natural biosystems. The principle of corrective action is based on the theories of use of the frequencies of R. Rife, R.Faull, F. Morell.

Caution! The device can be used by children only under the supervision of an adult. Children shall not play with the device. People with physical or mental disabilities may use the device only under the supervision of a competent person.

Set

| | |
|-------------------------------------------------|---------|
| Wellness device "DeVita Ritm+ Model MINI" | 1 piece |
| Passport / Instructions | 1 piece |
| Packaging | 1 piece |
| USB charging cable / data transfer..... | 1 piece |
| Neck strap..... | 1 piece |

Technical Data

| | |
|---------------------------------------------------------|---------------|
| 1. Maximum number of programs | 50 |
| 2. Maximum number of Complexes (folders)..... | 50 |
| 3. Maximum number of frequencies per program | 80 |
| 4. Frequency range, Hz | 0.1-100 |
| 5. Continuous operating time, no less than, hours | 20 |
| 6. Keyboard Type | Touch sensors |
| 7. Screen type, pixels | TFT 128x160 |
| 8. Updates via the Internet: | |
| - programs - Complexes..... | Yes |

9. Power supply:

| | |
|--------------------------------------------------|-----------|
| - Li-Po 3.7V battery, built-in device, pcs | 1 |
| 10. Overall dimensions, not more than, mm | 103x45x16 |
| 11. Device weight, not more than, kg | 0.1 |
| 12. Average life span, not less, years | 5 |

Safety Information

The wellness device is:

- Safe to use
- Non-toxic
- Not flammable
- Not explosive
- Not a source of dangerous radiation

When charging the battery, the device is powered by an exceptionally safe low-voltage power adapter.

Environmental Protection

The built-in battery contains pollutants that pollute the environment. After the device has ended its use and before transferring it for disposal, the battery must be removed from the service center. The battery should be disposed of in special areas.

Instructions for care, storage and transport

Device without packaging must be stored in a room temperature from 10 to 35 ° C and a relative humidity of 80%.


The device must be protected from direct humidity and impact of active chemicals.

Do not disassemble the device. Do not drop the device. Do not let other objects fall on the device.

We recommend keeping the packaging container within the period of the warranty.

Device description

1. ON/OFF button
2. Display
3. Charge/Data transfer port
4. Navigation keys
5. Selection key (OK)
6. Return key
7. Red LED indicator light - execution of program

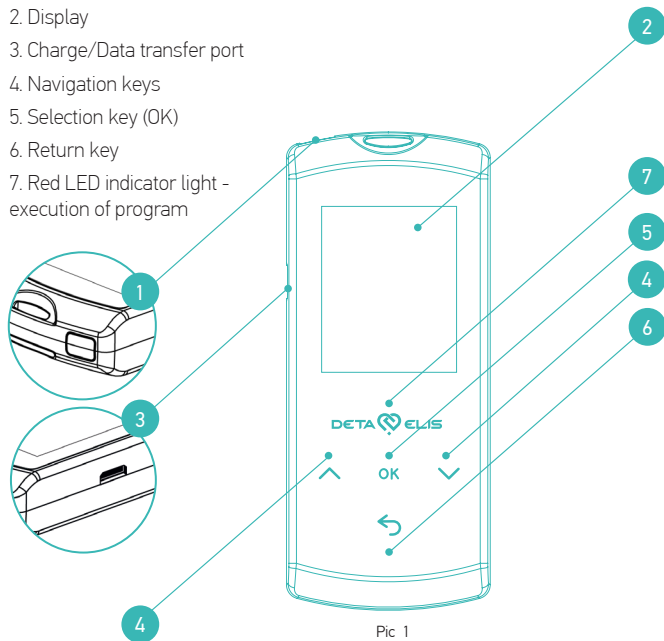
Pic 1

Working with the device


To turn on the device, press the ON/OFF button 1 of the device. **1** (pic.1).

The display illumination time is 30 seconds since the last key was pressed. The display can be turned off earlier than this time by pressing the ON/OFF button 1 **1** again. In the event of the display turning off before you enter the necessary commands, it must be turned on again by pressing ON/OFF button 1. **1**

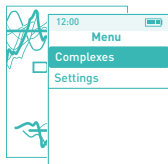After pressing the ON/OFF button 1 of the device (pic. 1), the home page will appear on the display, and after 2 seconds, the device will automatically enter the main menu mode.

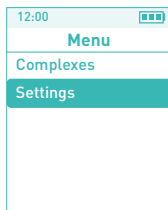
Working with the main menu

Submenu "Settings"

-  Use the navigation keys to select "Settings"
-  Press the selection key

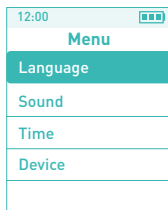
The submenu "Settings" appears on the display

Language selection



-  Use the navigation keys to select "Language"
-  Press the selection key

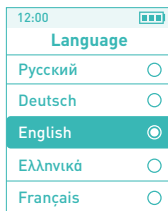The submenu "Language" appears on the display



-  Use the navigation keys to select the desired language
-  Press the selection key
-  To return to the "Settings" submenu press the return key

Attention! The default language is English. If there is no need for changing the language you can skip this step.

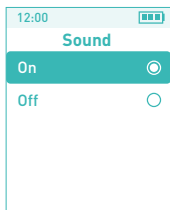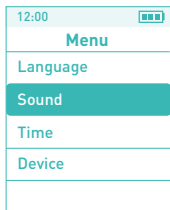If you select "Off", the sound of pressing any key is turned off.

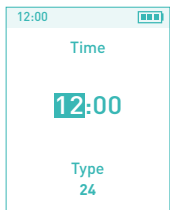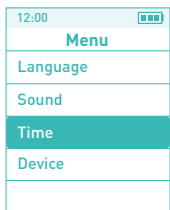Sound settings

 Use the navigation keys to select "Sound"

 Press the selection key

The submenu "Sound" appears on the display

 Use the navigation keys to select "On" or "Off"

 Press the selection key

 To return to the "Settings" submenu press the return key

Time settings

 Use the navigation keys to select "Time"


 Press the selection key

The submenu "Time" appears on the display

 Use the navigation keys to set hours, minutes and type

 Use the selection key to select "Hours", "Minutes" and "Type"

 To return to the "Settings" submenu press the return key

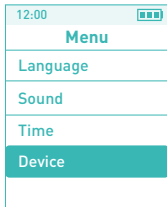About the device


Use the navigation keys to select "Device"

Press the selection key

The display will show the model name of the device, the manufacturer serial number of your device and the version of the program.

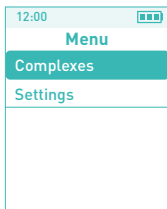Submenu "Complexes"


Use the navigation key to select "Complexes"

Press the selection keys

The submenu "Complexes" appears on the display.

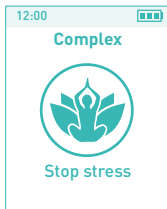Select a complex


Use the navigation keys to select the desired Complex (eg Stop Stress)

Press the selection key

The submenu of the selected Complex appears on the display

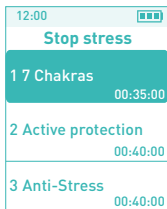Select a program


Use the navigation keys to select the desired program


Press the selection key

The program will start, which will be signaled by a flashing red LED indicator light pic.1 (7).

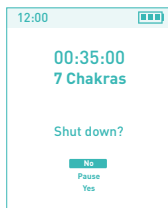
The submenu for running the current program appears on the display.

The timer at the top of the display shows the remaining time.

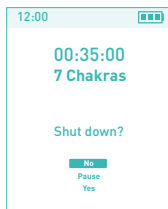Information about the program in progress

If you need to get information about the program in progress when the screen is turned-off, press the ON/OFF button 1 **1** of the device (pic. 1) and the submenu for running the current program appears on the display.

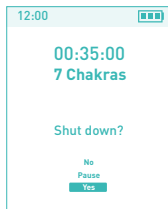Interrupting the program

If you need to shut down the program


 If you need to stop the execution of the program use the navigation keys to select "Yes"


 Press the selection key

If you need to pause the program

 Use the navigation keys to select "Pause"

 Press the selection key

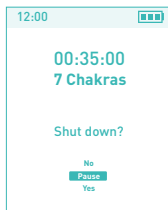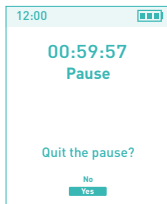The "Pause" submenu appears on the display
If you need to exit the pause

Use the navigation keys to select
"Yes"

Press the selection key

If you do not select the "Yes" button under the "Quit the pause?" option, the device will automatically turn off after 1 hour.

Attention! When the device is paused, the red LED indicator light pic. 1 (7) will flash at a slower rate!

Using the device in program mode


After starting the program, place the device with its back facing your body, for example, in your pocket.

Charge the device

Battery charge indicator

The battery indicator is highlighted in white. If three or two internal lines of the indicator are visible - the battery is charged.  

The indicator of the battery charge is highlighted in red, one internal line of the indicator is visible - the charge of the battery is running out. 

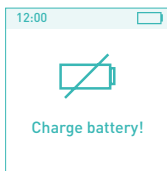The battery charge indicator is highlighted in red, none of the internal lines of the indicator are visible, the display shows "Charge battery!" - the battery is completely discharged.

Charging the battery

In the case of low charge or a complete discharge of the battery, connect the micro-USB connector of the cable to the Charge/Data transfer port pic 1 (3) of the device, and the USB connector to the USB port of a working computer, or an external charger with output 5V-1A. The battery charge indicator will start signaling the charging process.

Programming the device

The device is programmed using the "DeInfo Mini" program, recorded within a memory of the device. There is also the "User's manual of the DeInfo Mini program". To open the manual, connect the device with the charging/ data transfer cable to the USB input of a computer. Select the "RITM_MINI" disk and open the folder "Manual Deinfo Mini". For further instructions, follow the steps according to the Deinfo Mini manual.

Manufacturer's warranty

The manufacturer guarantees full compliance with the technical specifications of the device "DeVita Ritm+ Model Mini", provided that the user strictly observes the instructions for operation, transport and storage.

The warranty period of the device is 24 months from the date of sale through the retail network.

In the absence of the sale date and the stamp of the trade organization in the passport, the warranty period is calculated from the date of production of the device by the manufacturer.

During the warranty period, the owner is entitled to free repair. Warranty repairs are carried out at service centers authorized by the manufacturer. Delivery of the defective device to the repair site is carried out at the expense of the buyer.

If there is a violation of the external covers and / or in the absence of the certificate warranty, no claims are accepted and no warranty repair is performed.

The following faults are not covered by the warranty:

- defects due to improper operation;
- defects caused by natural disasters;
- in case of violation of the integrity of seals;
- in case of presence of external defects (cracks, fractures etc.)

The buyer has the right to replace the defective device with a new one during the warranty period in the following cases:

- in case of repairing the device for the fourth time
- if the device is unrepairable

Attention! The manufacturer reserves the right to change the design and software, without notice, without degrading existing characteristics.

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EC Declaration of Conformity Nr. W13

The company **DETA-ELIS EUROPA GMBH** with a seat at Kaiserstr. 74, 63065 Offenbach, Germany

Herewith to DECLARE

That the equipment described below complies with all applicable at the time of construction essential requirements of the **Directive 2011/65/EU (RoHS)** and the **Directive 2014/30/EU (EMC)**.

This declaration relates exclusively to the equipment in the state in which it was placed on the market, and excludes alterations or manipulations, components which are added and/or operations carried out subsequently by the final user without our official consent.

Equipment: DeVita Ritm
Type / Model: Mini
Standards Applied: EN 55014-1:2006+A1:2009+A2:2011
EN 55014-2:1997+A1:2001+A2:2008
EN 61000-3-2:2006+A1:2009+A2:2009
EN 61000-3-3:2008

Date: 23 January 2018

Authorized Representative: The General Manager Ms. Konopleva Tatyana

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CONCLUSION: IN REFER TO SGS TEST REPORT NO.: H0FA0005
NO DEVIATION FROM THE STANDARDS HAVE
BEEN IDENTIFIED. THEREFORE THE TESTED
PRODUCT IS IN COMPLIANCE WITH THE
APPROPRIATED STANDARDS.

This Verification of Compliance has been granted to the applicant based on the results of tests, performed by Laboratory of SGS-CSTC Standards Technical Services Co., Ltd. on sample of the above-mentioned product in accordance with the provisions of the relevant specific standards. The CE marking as shown below can be affixed, under the responsibility of the manufacturer, after completion of an EC Declaration of Conformity and compliance with all relevant EC Directives.

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SGS Germany GmbH, Hofmannstr. 50, D-81379 Munich is testing facility for
ELECTRO-MECH. COMPONENTS TESTING ENVIRONMENTAL SIMULATION ELECTROMAGNETIC COMPATIBILITY
MEDICAL DEVICE TESTING BATTERY TESTING PRODUCT SAFETY TELECOM CONFORMANCE TESTS

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Member of the SGS Group (Société Générale de Surveillance)

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DeTa-Elis Europa GmbH
Deta-Elis Europa GmbH